



**US  
DENTAL  
DEPOT**

**US DENTAL DEPOT**

6555 NW 9<sup>th</sup> Avenue suite 212  
Fort Lauderdale 33309  
Florida USA

USA - 1 1866 625 5558 – 1 954 874 6325

Colombia 57 1 489 9051

Peru 51 1 707 3444

Venezuela 58 212 729 1639

Argentina 54 11 5258 7148 Cba 54 351 569 1409

[www.USDENTALDEPOT.com](http://www.USDENTALDEPOT.com)



## **Cordless Gutta Percha Obturation System**

Item # 102222

### **System includes:**

- \*Obturation Gun
- \*Oburation Pen
- \*Gun Needle (3pcs)
- \*Charging Base
- \*Adapter
- \*Heat insulator (2pcs)
- \*Stand for gun
- \*Stand for pen
- \*Pen tips (4 pcs)
- \*Battery (3pcs)
- \*Accessory tool
- \*Brush (2pcs)
- \*English manual

### **Features:**

- \*Perfect cordless 3- dimensional obturation system
- \*Fantastic combination of filling method of vertical compaction and backfill
- \*Easy and quick obturation with accurate and predictable filling
- \*Excellent hermetic sealing with 3-dimensional filling  
(Total obturation of the entire canal including lateral canals)
- \*Excellent apical control & obturation of anatomically complicated canals
- \*Condensed gutta percha makes fill obturation



# US DENTAL DEPOT

## US DENTAL DEPOT

6555 NW 9<sup>th</sup> Avenue suite 212  
Fort Lauderdale 33309  
Florida USA

USA - 1 1866 625 5558 – 1 954 874 6325

Colombia 57 1 489 9051

Peru 51 1 707 3444

Venezuela 58 212 729 1639

Areentina 54 11 5258 7148 Cba 54 351 569 1409

[www.USDENTALDEPOT.com](http://www.USDENTALDEPOT.com)

- \*Digitizing temperature display contributes to convenient setting of temperature
- \*Easy to learn and operate
- \*Compact design comply with standard of human engineering

### Charging system

- \*Operating Voltage DC 3.6V
- \*Adaptor Input: AC 100-240V Output: DC 4.5V
- \*Charging time 3 hours
- \*Continuous heat-up time for gun after fully charging 60minutes
- \*Continuous heat-up time for pen after fully charging 40 minutes

### Obturation Gun

- \*Heater DC3.6 V 8W
- \*Temperature sensor Thermocouple
- \*Heating time 2 mintues
- \*Temperature range 50-200 degrees celsius

### Obturation Pen

- \*Temperature sensor K type
- \*Heating time 2-10 seconds
- \*Temperature range 150-350 degrees Celsius

### Battery (Li-on)

- \*Capacity DC3 3.6V 2.2AH
- \*Dimensions 25 x 75 ( 3 pcs total )
- \*Operating condition temperature 0-40 degrees celsius



# US DENTAL DEPOT

## US DENTAL DEPOT

6555 NW 9<sup>th</sup> Avenue suite 212  
Fort Lauderdale 33309  
Florida USA

USA - 1 1866 625 5558 – 1 954 874 6325

Colombia 57 1 489 9051

Peru 51 1 707 3444

Venezuela 58 212 729 1639

Argentina 54 11 5258 7148 Cba 54 351 569 1409

[www.USDENTALDEPOT.com](http://www.USDENTALDEPOT.com)




Charging Base



Accessory Tools




Battery

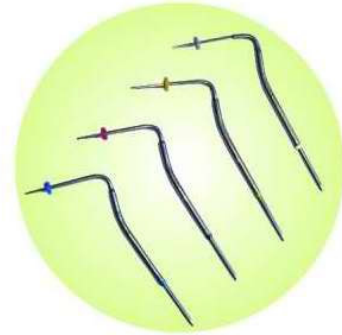

Heat Tips



Heat Insulator

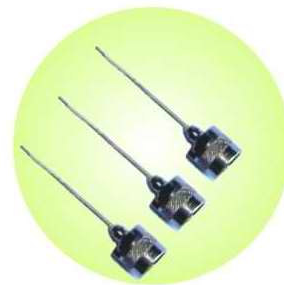

Injection Needle