

CAD/CAM

MILLING DISCS & BLOCKS



ZIRCONIA DISCS & BLOCKS

Dental zirconia discs are pre-sintered discs made of zirconium dioxide (ZrO₂) that are used in the fabrication of dental restorations such as crowns, bridges, and implant abutments. These discs are milled into the desired restoration shape using Computer-Aided Design/Computer-Aided Manufacturing (CAD/CAM) systems and then sintered in a furnace to achieve their final hardness and translucency.

Dental zirconia discs and their applications:



1. Types of Zirconia: There are various types of zirconia used in dentistry based on their aesthetic and strength characteristics. The most common types include:

- High-translucency zirconia (HT): Used for anterior restorations due to their superior aesthetics.
- High-strength zirconia: Typically used for posterior restorations because of their durability.
- Multilayer zirconia: These discs have varying layers of translucency and color, mimicking the natural gradation of teeth.

2. Milling Process: The zirconia disc is placed into a CAD/CAM milling machine. The design of the dental restoration is digitally sent to the machine, which then precisely mills the restoration out of the disc.

3. Sintering: After milling, the zirconia restoration is "green" or in a softer state. It's then placed into a sintering furnace where it's heated to very high temperatures. This process densifies the zirconia and gives it its final strength and translucency.

4. Characteristics:

- Strength: Zirconia is one of the strongest ceramic materials used in dentistry, making it ideal for restorations that require durability.
- Biocompatibility: Zirconia is biocompatible, meaning it's well-tolerated by the human body and doesn't cause allergic reactions.
- Aesthetics: Modern zirconia materials offer a balance of strength and aesthetics, making them suitable for both anterior and posterior restorations.

5. Applications: Dental zirconia is used for various applications like:

- Single crowns
- Bridges
- Implant abutments
- Inlays and onlays

6. Advantages:

- Durability: Zirconia restorations have a long lifespan and are resistant to wear and fractures.
- Aesthetics: With advancements in dental ceramics, zirconia can be made to closely match the appearance of natural teeth.
- Metal-free: For patients who want metal-free restorations or have allergies to certain metals, zirconia is an excellent alternative.

ZIRCONIA DISCS



ZIRCONIA DISCS | DIGITECH



MONOLAYER ZIRCONIA DISCS

98.5mm MONOLAYER ZIRCONIA DISCS

98.5mm MONOLAYER - DIGITECH

* Our Zirconia blocks, Standard Color
Zirconia Blocks Chemical compositions:

- * $ZrO_2 + HfO_2 + Y_2O_3$ 99.6 %
- * Y_2O_3 5.2 %
- * Al_2O_3 0.5 %
- * Bending Strength 1200 Mpa
- * Vickers Hardness 12.8 Gpa
- * COLD ISOSTATIC PRESSING with hydrostatic pressure.

Firing Chart:

- * Sinterization temperature
- HS and ST 1530 Celsius 2 hours
- AT 1450 Celsius 2 hours

CE & FDA CERTIFIED

#USDD98.5



98.5mm SUPER TRANSLUCENT Zirconia Discs

98.5mm PRE-SHADED ST SUPER TRANSLUCENT Zirconia Discs

St PRE 98.5mm 30Mm

466410	DIGITECH - ST PRE A1 Dental Zirconia discs Super Translucent 98.5mm (98mm) x...	US\$ 97.60	+Info
--------	---	------------	-------

ST PRE 98.5mm 20mm

437067	DIGITECH - ST PRE A1 Dental Zirconia discs Super Translucent 98.5mm (98mm) x...	US\$ 67.20	+Info
437055	DIGITECH - ST PRE A2 Dental Zirconia discs Super Translucent 98.5mm (98mm) x...	US\$ 67.20	+Info
437054	DIGITECH - ST PRE BL1 Dental Zirconia discs Super Translucent 98.5mm (98mm) ...	US\$ 54.40	+Info

ST PRE 98.5mm 14mm

124264	DIGITECH - ST PRE A1 Dental Zirconia discs Super Translucent 98.5mm (98mm) x...	US\$ 54.40	+Info
124265	DIGITECH - ST PRE A2 Dental Zirconia discs Super Translucent 98.5mm (98mm) x...	US\$ 54.40	+Info
417394	DIGITECH - ST PRE A3 Dental Zirconia discs Super Translucent 98.5mm (98mm) x...	US\$ 54.40	+Info
124268	DIGITECH - ST PRE A3.5 Dental Zirconia discs Super Translucent 98.5mm (98mm)...	US\$ 50.05	+Info
437053	DIGITECH - ST PRE BL1 Dental Zirconia discs Super Translucent 98.5mm (98mm) ...	US\$ 54.40	+Info

ST PRE 98.5mm 18mm

124266	DIGITECH - ST PRE A2 Dental Zirconia discs Super Translucent 98.5mm (98mm) x...	US\$ 65.60	+Info
124269	DIGITECH - ST PRE A3.5 Dental Zirconia discs Super Translucent 98.5mm (98mm)...	US\$ 65.60	+Info

ST PRE 98.5mm 16mm

417080	DIGITECH - ST PRE A1 Dental Zirconia discs Super Translucent 98.5mm (98mm) x...	US\$ 62.40	+Info
416997	DIGITECH - ST PRE A2 Dental Zirconia discs Super Translucent 98.5mm (98mm) x ...	US\$ 60.80	+Info

ST PRE 98.5mm 25mm

437056	DIGITECH - ST PRE A2 Dental Zirconia discs Super Translucent 98.5mm (98mm) x...	US\$ 81.60	+Info
--------	---	------------	-----------------------

98.5mm WHITE ST SUPER TRANSLUCENT Zirconia Discs

ST MONO WHITE

ST MONO WHITE

110670	Zirconia Bck ST 98.5mm x 16	US\$ 67.20	+Info
112912	Zirconia Bck ST 98.5mm x 18	US\$ 60.80	+Info
112945	Zirconia Bck ST 98.5mm x 14	US\$ 51.20	+Info

113551	Zirconia Blck ST 98.5mm x 20	US\$ 64.00	+Info
113552	Zirconia Blck ST 98.5mm x 22	US\$ 70.40	+Info
114937	Zirconia Blck ST 98.5mm x 25	US\$ 79.68	+Info
116013	Zirconia Blck ST 98.5mm x 12	US\$ 43.20	+Info
116041	Zirconia Blck ST 98.5mm x 10	US\$ 35.20	+Info

Super Translucent

Super Translucent

421161	Zirconia Blck ST 98.5mm x 20	US\$ 69.88	+Info
421162	Zirconia Blck ST 98.5mm x 22	US\$ 82.59	+Info
421163	Zirconia Blck ST 98.5mm x 25	US\$ 70.40	+Info

98.5mm AT WHITE ANTERIOR TRANSLUCENT Zirconia Discs

Code	Name	Price	+
123227	DIGITECH - MONO AT Dental Zirconia discs Anterior Translucency AT 98.5mm (98m...	US\$ 97.60	+Info
123228	DIGITECH - MONO AT Dental Zirconia discs Anterior Translucency AT 98.5mm (98m...	US\$ 91.52	+Info
123229	DIGITECH - MONO AT Dental Zirconia discs Anterior Translucency AT 98.5mm (98m...	US\$ 73.97	+Info

98.5mm HS HIGH STRENGTH Zirconia Discs

Code	Name	Price	+
101972	DIGITECH - HS MONO Dental Zirconia discs 98.5mm (98mm) x 10mm - 1 p/box - #...	US\$ 30.40	+Info
111119	DIGITECH - HS MONO Dental Zirconia discs 98.5mm (98mm) x 12mm - 1 p/box - #...	US\$ 36.80	+Info
101973	DIGITECH - HS MONO Dental Zirconia discs 98.5mm (98mm) x 14mm - 1 p/box - #...	US\$ 45.90	+Info
101974	DIGITECH - HS MONO Dental Zirconia discs 98.5mm (98mm) x 16mm - 1 p/box - #...	US\$ 49.60	+Info
101976	DIGITECH - HS MONO Dental Zirconia discs 98.5mm (98mm) x 20mm - 1 p/box - #...	US\$ 56.00	+Info
113481	DIGITECH - HS MONO Dental Zirconia discs 98.5mm (98mm) x 25mm - 1 p/box - #...	US\$ 65.60	+Info

93mm ZIRCONIA DISCS MONOLAYER

ZIRCONIA DISCS MONOLAYER for AMANN GIRRBACH - DIGITECH

* Our Zirconia discs , Standard Color

Zirconia Blocks Chemical compositions:

- * ZrO₂+HfO₂+Y₂O₃ 99.6 %
- * Y₂O₃ 5.2 %
- * Al₂O₃ 0.5 %
- * Bending Strength 1200 Mpa
- * Vickers Hardness 12.8 Gpa
- * COLD ISOSTATIC PRESSING with hydrostatic pressure.



Firing Chart:

- * 44% Translucency
- * 1200 Mpa
- * Sinterization temperature
- HS and ST 1530 Celsius 2 hours
- AT 1450 Celcius 2 hours

#UsddAmann

AMANN GIRRBACH is a trademark of Amann Girrbach AG

SUPER TRANSLUCENT 1200Mpa

ST MONO WHITE 93mm

Code	Name	Price	+
121052	DIGITECH - ST Dental Zirconia discs Super Translucent AG 93mm x 71mm x 12mm - ...	US\$ 28.80	+Info
115346	DIGITECH - ST Dental Zirconia discs Super Translucent AG 93mm x 71mm x 14mm - ...	US\$ 33.41	+Info

PRE-SHADED SUPER TRANSLUCENT

ST 93mm 18mm

421058	DIGITECH - ST PRE A3.5 Zirconia discs Super Translucent AG 93mm x 71mm x 18mm ...	US\$ 44.93	+Info
--------	---	------------	-------

ST 93mm 20mm

124515	DIGITECH - ST PRE A2 Dental Zirconia discs Super Translucent AG 93mm x 71mm x ...	US\$ 48.38	+Info
124518	DIGITECH - ST PRE A3.5 Zirconia discs Super Translucent AG 93mm x 71mm x 20mm ...	US\$ 48.38	+Info

95mm ZIRCONIA DISCS MONOLAYER

95mm ZIRCONIA Discs CAD/CAM for ZIRKINZAHN

Characterisics:

- * First grade Zirconia powder
- * Shrinkage adjusted for 20% -
- * Compensation 1.25

Physical Qualities:

- * Bending Strength:
- HS High Strength 1400 Mpa
- ST Super Translucent 1200 Mpa
- AT Anterior Translucency 600 Mpa (TT Total Translucency)
- * Vickers Hardness: 12.8 Gpa



Colors:

- * Our Zirconia blocks, are Standard LL0.5 .
- * For darker shades ECCO liquid Coloring is available in the following: *LL1
- *LL2 *LL3 *LL4 *LL5 as the darkest shade.
- * Coloring for Full contour is also available in 16 Vita shades

Firing Chart:

- * Starting temperature: Room Temperature
- * Temperature Heating rate 4 Celcius per minute
- * Sinterization temperature
- HS and ST 1530 Celsius 2 hours
- AT 1450 Celcius 2 hours
- * Temperature Coling Rate 4 Celcius per minute
- * Total Cycle time: 16 hours .Then open door, don?t forget to decontaminate your sinter oven once a week, or at least once a month. This advice is for every sinter oven in the market, regardless of brand name, or origin.

#Usdd95

#ZirkonZahn

ZIRKONZAHN is a trademark of ZIRKONZAHN GMBH

HS HIGH STRENGTH

Code	Name	Price	+
112451	DIGITECH - HS MONO Dental Zirconia discs ZZ 95mm x 14mm - (HS MONO) - 1 p/box...	US\$ 44.80	+Info

SUPER TRANSLUCENT

ST WHITE 95mm

Code	Name	Price	+
116085	DIGITECH - ST Dental Zirconia discs Super Translucent - ZZ 95mm x 12mm - 1 p/b...	US\$ 43.20	+Info
112450	DIGITECH - ST Dental Zirconia discs Super Translucent - ZZ 95mm x 14mm - 1 p/b...	US\$ 51.20	+Info
116084	DIGITECH - ST Dental Zirconia discs Super Translucent - ZZ 95mm x 16mm - 1 p/b...	US\$ 57.60	+Info
417138	DIGITECH - ST Dental Zirconia discs Super Translucent - ZZ 95mm x 18mm - 1 p/b...	US\$ 62.40	+Info
114145	DIGITECH - ST Dental Zirconia discs Super Translucent - ZZ 95mm x 20mm - 1 p/b...	US\$ 64.00	+Info
114403	DIGITECH - ST Dental Zirconia discs Super Translucent - ZZ 95mm x 25mm - 1 p/b...	US\$ 76.80	+Info
417137	DIGITECH - ST Dental Zirconia discs Super Translucent - ZZ 95mm x 30mm - 1 p/b...	US\$ 89.60	+Info

PRE-SHADED ST SUPER TRANLUCENT

PRE-SHADED ST SUPER TRANLUCENT Zirconia Discs

ST 95mm 20mm

419954	DIGITECH - ST PRE A2 Dental Zirconia discs Super Translucent - ZZ 95mm X 20mm ...	US\$ 67.20	+Info
--------	--	------------	-----------------------

ST 95mm 22mm

419955	DIGITECH - ST PRE A2 Dental Zirconia discs Super Translucent - ZZ 95mm X 22mm ...	US\$ 75.20	+Info
--------	--	------------	-----------------------

ST 95mm 25mm

419953	DIGITECH - ST PRE A1 Dental Zirconia discs Super Translucent - ZZ 95mm X 25mm ...	US\$ 81.60	+Info
417392	DIGITECH - ST PRE A2 Dental Zirconia discs Super Translucent - ZZ 95mm X 25mm ...	US\$ 81.60	+Info

ST 95mm 30mm

419956	DIGITECH - ST PRE A1 Dental Zirconia discs Super Translucent - ZZ 95mm X 30mm ...	US\$ 97.60	+Info
--------	--	------------	-----------------------

ST A2 95mm 14mm

124527	DIGITECH - ST PRE A1 Dental Zirconia discs Super Translucent - ZZ 95mm x 14mm ...	US\$ 54.40	+Info
420168	DIGITECH - ST PRE A2 Dental Zirconia discs Super Translucent - ZZ 95mm x 14mm ...	US\$ 53.31	+Info
124529	DIGITECH - ST PRE A3 Dental Zirconia discs Super Translucent - ZZ 95mm x 14mm ...	US\$ 54.40	+Info
124530	DIGITECH - ST PRE A3.5 Dental Zirconia discs Super Translucent - ZZ 95mm x 14mm...	US\$ 54.40	+Info

ST A3 95mm 18mm

124528	DIGITECH - ST PRE A2 Dental Zirconia discs Super Translucent - ZZ 95mm X 18mm ...	US\$ 64.00	+Info
124531	DIGITECH - ST PRE A3.5 Dental Zirconia discs Super Translucent - ZZ 95mm x 18m...	US\$ 65.60	+Info

AT ANTERIOR TRANSLUCENT

12mm

12mm

Code	Name	Price	+
123221	ZirkonZahn AT 95x12 w/marks	US\$ 62.40	+Info

14mm

14mm

Code	Name	Price	+
123222	ZirkonZahn AT 95x14 w/marks	US\$ 94.67	+Info

16mm

16mm

Code	Name	Price	+
123223	ZirkonZahn AT 95x16 w/marks	US\$ 83.20	+Info

18mm

18mm

Code	Name	Price	+
123224	ZirkonZahn AT 95x18 w/marks	US\$ 103.98	+Info

20mm

20mm

Code	Name	Price	+
123225	ZirkonZahn AT 95x20 w/marks	US\$ 111.74	+Info

98mm ZIRCONIA DISCS MONOL AYER

ZIRCONIA BLOCKS 98mm - CAD/CAM - DIGITECH

* Our Zirconia blocks, Standard Color

Zirconia Blocks Chemical compositions:

* ZrO₂+HfO₂+Y₂O₃ 99.6 %

* Y₂O₃ 5.2 %

* Al₂O₃ 0.5 %

* Bending Strength 1200 Mpa

* Vickers Hardness 12.6 Gpa

* COLD ISOSTATIC PRESSING with hydrostatic pressure.

Firing Chart:

* Sinterization temperature

HS and ST 1530 Celsius 2 hours

AT 1450 Celsius 2 hours



HIGH STRENGTH 1200Mpa

Code	Name	Price	+
113482	DIGITECH - Dental Zirconia discs High Strength HS MONO 98mm (is NOT 98.5mm) ...	US\$ 54.40	+Info
113534	DIGITECH - HS MONO Dental Zirconia discs High Strength HS MONO 98mm (is NOT 9...	US\$ 61.15	+Info
113533	DIGITECH - HS MONO Dental Zirconia discs High Strength 98mm (is NOT 98.5mm) ...	US\$ 65.60	+Info

SUPER TRANSLUCENT 1200Mpa

Code	Name	Price	+
117369	DIGITECH -ST MONO Dental Zirconia discs ST MONO Super TranslucentST MONO 98mm (...)	US\$ 36.80	+Info
116061	DIGITECH -ST MONO Dental Zirconia discs Super TransST MONO 98mm (is NOT 98.5m...	US\$ 43.92	+Info
116062	DIGITECH -ST MONO Dental Zirconia discs Super TransST MONO 98mm (is NOT 98.5m...	US\$ 52.80	+Info

116063	DIGITECH -ST MONO Dental Zirconia discs Super TransST MONO 98mm (is NOT 98.5m...	US\$ 57.50	+Info
121051	DIGITECH -ST MONO Dental Zirconia discs Super TranslucentST MONO 98mm (is NOT...	US\$ 62.40	+Info

MULTILAYER ZIRCONIA DISCS

4D MULTILAYER Zirconia Discs

MULTILAYER 4D TRANSLUCENT ZIRCONIA - DIGITECH

4D GRADUAL TRANSLUCENCY :

- * 39 to 50% Translucency
- * 900 Mpa
- * Suitable for full contour crown and bridge
- * Multi layer gradient translucency with 3 popular shades

* Same sintering process as regular ST blocks



4D 98.5mm Multilayer Zirconia Discs

4D AT ML AESTHETIC 98mm Zirconia Discs

14mm

124138	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 14mm A2 Multi-Layer Anterior...	US\$ 80.00	+Info
124139	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 14mm C2 Multi-Layer Anterior...	US\$ 77.60	+Info
124141	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 14mm A1 Multi-Layer Anterior...	US\$ 80.00	+Info
124142	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 14mm B1 Multi-Layer Anterior...	US\$ 80.00	+Info
416769	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 14mm A3 Multi-Layer Anterior...	US\$ 80.00	+Info
416770	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 14mm C3 Multi-Layer Anterior...	US\$ 80.00	+Info
416771	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 14mm B2 Multi-Layer Anterior...	US\$ 78.40	+Info
416772	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 14mm A3.5 Multi-Layer Anteri...	US\$ 80.00	+Info
417364	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 14mm C1 Multi-Layer Anterior...	US\$ 80.00	+Info
417365	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 14mm D2 Multi-Layer Anterior...	US\$ 80.00	+Info
434838	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 14mm BL1 Multi-Layer Translu...	US\$ 80.00	+Info

437064	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 14mm BL2 Multi-Layer Translu...	US\$ 80.00	+Info
474498	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 14mm D3 Multi-Layer Anterior...	US\$ 80.00	+Info
474923	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 14mm A4 Multi-Layer Anterior...	US\$ 80.00	+Info

16mm

474937	DIGITECH - ML 4D AESTHETIC AT Dental Zirconia discs 98.5mm x 16mm A3 Multi-Laye...	US\$ 88.00	+Info
474938	DIGITECH - ML 4D AESTHETIC AT Dental Zirconia discs 98.5mm x 16mm B2 Multi-Laye...	US\$ 88.00	+Info
474939	DIGITECH - ML 4D AESTHETIC AT Dental Zirconia discs 98.5mm x 16mm A2 Multi-Laye...	US\$ 88.00	+Info

18mm

437066	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 18mm A1 Multi-Layer Anterior...	US\$ 94.40	+Info
476662	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 18mm BL1 Multi-Layer Anterior...	US\$ 94.40	+Info
476663	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 18mm BL2 Multi-Layer Anterior...	US\$ 94.40	+Info
417346	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 18mm A2 Multi-Layer Anterior...	US\$ 94.40	+Info
124137	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 18mm A3 Multi-Layer Anterior...	US\$ 94.40	+Info
124140	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 18mm A4 Multi-Layer Anterior...	US\$ 94.40	+Info
417335	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 18mm B1 Multi-Layer Anterior...	US\$ 94.40	+Info
417334	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 18mm C1 Multi-Layer Anterior...	US\$ 94.40	+Info
474916	DIGITECH - ML 4D AT Dental Zirconia discs 98.5mm x 18mm C2 Multi-Layer Anterior...	US\$ 94.40	+Info

20mm

4D ML Anterior Translucent

4D ML Anterior Translucent

417336	Zirc Disc ML 98.5x20 4D A1	US\$ 99.00	+Info
417337	Zirc Disc ML 98.5x20 4D A2	US\$ 99.00	+Info
417338	Zirc Disc ML 98.5x20 4D B1	US\$ 97.45	+Info
437065	Zirc Disc ML 98.5x20 4D BL1	US\$ 99.00	+Info
437070	Zirc Disc ML 98.5x20 4D BL2	US\$ 99.00	+Info
470502	Zirc Disc ML 98.5x20 4D A3	US\$ 99.00	+Info

4D AT ML PRO 98mm Zirconia Discs

16mm PRO

474917	DIGITECH - ML 4D-PRO AT Dental Zirconia discs 98.5mm x 16mm A1 Multi-Layer Ant...	US\$ 88.00	+Info
474918	DIGITECH - ML 4D-PRO AT Dental Zirconia discs 98.5mm x 16mm A2 Multi-Layer Ant...	US\$ 88.00	+Info
474919	DIGITECH - ML 4D-PRO AT Dental Zirconia discs 98.5mm x 16mm B1 Multi-Layer Ant...	US\$ 88.00	+Info
421056	DIGITECH - ML 4D-PRO AT Dental Zirconia discs 98.5mm x 16mm BL1 BLEACH Multi-L...	US\$ 88.00	+Info
421057	DIGITECH - ML 4D-PRO AT Dental Zirconia discs 98.5mm x 16mm BL2 BLEACH Multi-...	US\$ 92.55	+Info

4D ML FUNCTIONAL 98mm Zirconia discs

14mm

437057	ZircDisc ML98.5x14 4D FUNCTIONALBL1	US\$ 80.00	+Info
437058	Zirc Disc ML 98.5x14 4D FUNCTIONAL BL2	US\$ 80.00	+Info

16MM

474928	DIGITECH - ML 4D FUNCTIONAL 98x16 BL1 Dental Zirconia discs Multi-Layer 98.5mm ...	US\$ 92.80	+Info
474929	DIGITECH - ML 4D FUNCTIONAL 98x16 BL2 Dental Zirconia discs Multi-Layer 98.5mm ...	US\$ 92.80	+Info
474930	DIGITECH - ML 4D FUNCTIONAL 98x16 A1 Dental Zirconia discs Multi-Layer 98.5mm (...)	US\$ 92.80	+Info
474931	DIGITECH - ML 4D FUNCTIONAL 98x16 A2 Dental Zirconia discs Multi-Layer 98.5mm (...)	US\$ 92.80	+Info
474932	DIGITECH - ML 4D FUNCTIONAL 98x16 A3 Dental Zirconia discs Multi-Layer 98.5mm (...)	US\$ 92.80	+Info
475035	DIGITECH - ML 4D FUNCTIONAL 98x16 B1 Dental Zirconia discs Multi-Layer 98.5mm (...)	US\$ 92.80	+Info

20mm

474354	Zirc Disc ML 98.5x20 4D FUNC BL3	US\$ 99.00	+Info
474329	Zirc Disc ML 98.5x20 4D FUNC A1	US\$ 97.45	+Info
474330	Zirc Disc ML 98.5x20 4D FUNC A2	US\$ 99.00	+Info
474331	Zirc Disc ML 98.5x20 4D FUNC A3	US\$ 99.00	+Info
474361	Zirc Disc ML 98.5x20 4D FUNC A4	US\$ 99.00	+Info
475034	Zirc Disc ML 98.5x20 4D FUNC B1	US\$ 99.00	+Info
474332	Zirc Disc ML 98.5x20 4D FUNC B2	US\$ 99.00	+Info
437059	Zirc Disc ML 98.5x20 4D FU	US\$ 99.00	+Info
437061	Zirc Disc ML 98.5x20 4D FUNCTIONAL BL2	US\$ 99.00	+Info

25MM

474924	Disc ML 98.5x25 4D A3	US\$ 115.84	+Info
474925	Disc ML 98.5x25 4D A1	US\$ 115.84	+Info
476414	Disc ML 98.5x30 4D A1	US\$ 121.27	+Info
476415	Disc ML 98.5x30 4D BL1	US\$ 121.27	+Info

4D MULTILAYER 95mm Zirconia Discs

4D AT ANTERIOR TRANSLUCENCY 95mm Zirconia Discs

14mm

416634	Z Disc 4D ML AT 95x14 4D A1	US\$ 80.00	+Info
416635	Zirc Disc 4D AT ML 95x14 A2	US\$ 80.00	+Info
416636	Zirc Disc 4D AT ML 95x14 B1	US\$ 76.80	+Info
417344	Zirc Disc 4D AT ML 95x14 A3	US\$ 77.60	+Info
417345	Zirc Disc 4D AT ML 95x14 B2	US\$ 80.00	+Info

18mm

417342	Zirc Disc 4D AT ML 95x18 4D A1	US\$ 94.40	+Info
416632	Zirc Disc 4D AT ML 95x18 4D A2	US\$ 70.80	+Info
417341	Zirc Disc 4D AT ML 95x18 4D A3	US\$ 94.40	+Info
416633	Zirc Disc ML 95x18 4D A3.5	US\$ 94.40	+Info

4D FUNCTIONAL 95mm Zirconia Discs

Code	Name	Price	+Info
474356	Zirc Disc ML 95x20 4D FUNCT BL1	US\$ 89.10	+Info
474357	Zirc Disc ML 95x20 4D FUNCT BL2	US\$ 89.10	+Info
474358	Zirc Disc ML 95x20 4D FUNCT A1	US\$ 99.00	+Info
474359	Zirc Disc ML 95x20 4D FUNCT B1	US\$ 93.06	+Info

4D MULTILAYER AMANN Zirconia Discs

4D AT ANTERIOR TRANSLUCENCY 95mm Zirconia Discs

14mm

416628	DIGITECH - ML 4D AT Dental Zirconia discs 93mm x 71mm x 14mm A1 Multi-Layer AM...	US\$ 58.06	+Info
--------	---	------------	-----------------------

16mm

417302	DIGITECH - ML 4D AT Dental Zirconia discs 93mm x 71mm x 16mm A2 Multi-Layer AM...	US\$ 39.60	+Info
420833	DIGITECH - ML 4D AT Dental Zirconia discs 93mm x 71mm x 16mm B2 Multi-Layer AM...	US\$ 39.60	+Info

420832	DIGITECH - ML 4D AT Dental Zirconia discs 93mm x 71mm x 16mm C2 Multi-Layer AM...	US\$ 39.60	+Info
420835	DIGITECH - ML 4D AT Dental Zirconia discs 93mm x 71mm x 16mm D2 Multi-Layer AM...	US\$ 39.60	+Info

18mm

417343	DIGITECH - 4D ML AT Dental Zirconia discs 93mm x 71mm x 18mm A2 Multi-Layer AM...	US\$ 67.29	+Info
--------	---	------------	-----------------------

20mm

420836	DIGITECH - ML 4D AT Dental Zirconia discs 93mm x 71mm x 20mm A3 Multi-Layer AM...	US\$ 71.28	+Info
420838	DIGITECH - ML 4D AT Dental Zirconia discs 93mm x 71mm x 20mm B2 Multi-Layer AM...	US\$ 71.28	+Info
420837	DIGITECH - ML 4D AT Dental Zirconia discs 93mm x 71mm x 20mm C2 Multi-Layer AM...	US\$ 71.28	+Info
420839	DIGITECH - ML 4D AT Dental Zirconia discs 93mm x 71mm x 20mm D2 Multi-Layer AM...	US\$ 71.28	+Info

AT ANTERIOR TRANSLUCENCY MULTILAYER Zirconia Discs

MULTILAYER AT TRANSLUCENT ZIRCONIA - DIGITECH



AT - Anterior Translucent - Multilayer (Also Called TT)

- * 48% Translucency
- * 600 Mpa
- * Suitable for full contour crown and bridage
- * Multi layer gradient translucency with 3 popular shades
- * Same sintering process as regular ST blocks

AT MULTILAYER 93mm Zirconia Discs

14mm

Code	Name	Price	+
461129	DIGITECH - ML AT Dental Zirconia discs 93mm x 14mm A2 Multi-Layer Anterior Tran...	US\$ 84.48	+Info

16mm

4D ML Ant. Translucent

4D ML Ant. Translucent

416768	Zirc Disc AG 4D AT ML 93x16 A3	US\$ 39.60	+Info
--------	--------------------------------	------------	-----------------------

ML Anterior Translucent

ML Anterior Translucent

416629	Zirc Disc AG AT ML 93x16 A2	US\$ 91.44	+Info
416630	Zirc Disc AG AT ML 93x16 -B2	US\$ 97.28	+Info
416631	Zirc Disc AG AT ML 93x16 D2	US\$ 91.44	+Info

18mm

Code	Name	Price	+
416626	DIGITECH - ML AT Dental Zirconia discs 93mm x 18mm A3 Multi-Layer Anterior Tran...	US\$ 99.84	+Info

AT MULTILAYER 98mm Zirconia Discs

ML Anterior Translucent

Code	Name	Price	+
420164	Zirc Disc ML 98.5x20 B3	US\$ 99.79	+Info
420165	Zirc Disc ML 98.5x20 C1	US\$ 97.60	+Info

ST MULTILAYER SUPERTRANSLUCENT Zirconia Discs

93mm ST MULTILAYER SUPER TRANSLUCENT Zirconia Discs

MULTILAYER ST SUPER TRANSLUCENT AMANN 93mm - DIGITECH

Features:

- * 46% Translucency
- * 1200 Mpa
- * Suitable for full contour crown and bridge
- * Multi layer gradient translucency with 3 popular shades
- * Same sintering process as regular ST blocks

Available in color: A1 - A2 - A3 - A3.5
- A4 - B2 - B3 - C2 - C3



ML ST 93mm 14mm

Multi-Layer Super Translucent

Multi-Layer Super Translucent

124519	Zirc Disc ML A2 93x75x14	US\$ 44.93	+Info
--------	--------------------------	------------	-----------------------

MultiLayer Super Translucent

MultiLayer Super Translucent

416720	Zirc Disc ML A3 93x75x14	US\$ 44.93	+Info
--------	--------------------------	------------	-------

ML ST 93mm 16mm

ML Super Translucent

ML Super Translucent

416781	Zirc Disc ML B2 93x75x16	US\$ 49.54	+Info
--------	--------------------------	------------	-------

MultiLayer Super Translucent

MultiLayer Super Translucent

416722	Z discs ML C3 93x75x16	US\$ 49.54	+Info
--------	------------------------	------------	-------

ML ST 93mm 18mm

Code	Name	Price	+
124521	DIGITECH - ML ST A2 Dental Zirconia discs 93mmx71mmx18mm - for Ceramill Motion ...	US\$ 54.14	+Info
124525	DIGITECH - ML ST A3 Dental Zirconia discs 93mmx71mmx18mm - for Ceramill Motion ...	US\$ 51.43	+Info

ML ST 93mm 20mm

Code	Name	Price	+
416747	DIGITECH - ML ST A3 Dental Zirconia discs 93mmx71mmx20mm - for Ceramill Motion ...	US\$ 55.81	+Info

95mm ST MULTILAYER SUPER TRANSLUCENT Zirconia Discs MULTILAYER ST SUPER TRANSLUCENT ZIRKONZAHN 95mm - DIGITECH

Features:

- * 46% Translucency
- * 1200 Mpa
- * Suitable for full contour crown and bridge
- * Multi layer gradient translucency with 3 popular shades
- * Same sintering process as regular ST blocks

Available in color: A1 - A2 - A3 - A3.5
- A4 - B2 - B3 - C2 - C3



ML ST 95mm 18mm

Code	Name	Price	+
124532	DIGITECH - ML ST Dental Zirconia discs ML A2 super Translucent - ZZ 95mmX18mm -...	US\$ 75.20	+Info
124533	DIGITECH - ML ST Dental Zirconia discs ML A3.5 Super Translucent - ZZ 95mmX18mm...	US\$ 66.18	+Info
124672	DIGITECH - ML ST Dental Zirconia discs ML C2 Super Translucent - ZZ 95mmX18mm -...	US\$ 70.69	+Info
124673	DIGITECH - ML ST Dental Zirconia discs ML D2 Super Translucent - ZZ 95mmX18mm -...	US\$ 74.45	+Info
470860	DIGITECH - ML ST Dental Zirconia discs ML A2 super Translucent - ZZ 95mmX16mm -...	US\$ 73.60	+Info

98.5mm ST MULTILAYER SUPER TRANSLUCENT Zirconia Discs



ML ST 98.5mm 14mm

Code	Name	Price	+
124128	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 14mm C2 Multi-Layer Su...	US\$ 64.00	+Info
124129	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 14mm A1 Multi-Layer Su...	US\$ 64.00	+Info
124130	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 14mm B1 Multi-Layer Su...	US\$ 64.00	+Info
124131	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 14mm A2 Multi-Layer Su...	US\$ 64.00	+Info
124135	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 14mm A3.5 Multi-Layer ...	US\$ 64.00	+Info
416713	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 14mm A3 Multi-Layer Su...	US\$ 64.00	+Info
416715	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 14mm B2 Multi-Layer Su...	US\$ 64.00	+Info

ML ST 98.5mm 16mm

Code	Name	Price	+
416714	DIGITECH - ML ST Dental Zirconia Block 98.5mm (98mm) x 16mm A4 MultiLayer Supe...	US\$ 73.60	+Info
416716	DIGITECH - ML ST Dental Zirconia Disc 98.5mm (98mm) x 16mm A3.5 MultiLayer Sup...	US\$ 73.60	+Info
416717	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 16mm C3 Multi-Layer Su...	US\$ 73.60	+Info
437062	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 16mm A2 Multi-Layer Su...	US\$ 73.60	+Info
470494	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 16mm A1 Multi-Layer Su...	US\$ 73.60	+Info
470496	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 16mm B1 Multi-Layer Su...	US\$ 73.60	+Info
470497	DIGITECH - ML ST Dental Zirconia Disc 98.5mm (98mm) x 16mm A3 MultiLayer Super...	US\$ 73.60	+Info

ML ST 98.5mm 18mm

Code	Name	Price	+
124132	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 18mm A1 Multi-Layer Su...	US\$ 75.20	+Info
124133	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 18mm A2 Multi-Layer Su...	US\$ 75.20	+Info
124134	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 18mm A3.5 Multi-Layer ...	US\$ 75.20	+Info
124136	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 18mm B1 Multi-Layer Su...	US\$ 75.20	+Info
124522	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 18mm C2 Multi-Layer Su...	US\$ 75.20	+Info
400598	DIGITECH - ML ST Dental Zirconia Discs 98.5mm (98mm) x 18mm D2 Multi-Layer Sup...	US\$ 80.00	+Info
451809	DIGITECH - ML ST Dental Zirconia Discs 98.5mm (98mm) x 18mm A3 Multi-Layer Sup...	US\$ 80.00	+Info

ML ST 98.5mm 20mm

Code	Name	Price	+
124523	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 20mm A3.5 Multi-Layer ...	US\$ 80.00	+Info
124524	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 20mm A2 Multi-Layer Su...	US\$ 80.00	+Info
124666	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 20mm C2 Multi-Layer Su...	US\$ 80.00	+Info
437060	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 20mm A1 Multi-Layer Su...	US\$ 80.00	+Info
470500	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 20mm A3 Multi-Layer Su...	US\$ 80.00	+Info
470501	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 20mm B1 Multi-Layer Su...	US\$ 80.00	+Info

ML ST 98.5mm 25mm

Code	Name	Price	+
124667	DIGITECH - ML ST Dental Zirconia discs 98.5mm (98mm) x 25mm A2 Multi-Layer Su...	US\$ 97.60	+Info

ZIRCONIA DISCS | ACONIA



SHTML Super High Translucent Multi Layer 98.5mm

SHTML

14mm

14mm

Code	Name	Price	+
476131	Aconia SHTML 98*14MM C2	US\$ 84.70	+Info
476140	Aconia SHTML 98*14MM OM3	US\$ 84.70	+Info
476139	Aconia SHTML 98*14MM OM2	US\$ 84.70	+Info
476138	Aconia SHTML 98*14MM OM1	US\$ 84.70	+Info
476137	Aconia SHTML 98*14MM D4	US\$ 84.70	+Info
476136	Aconia SHTML 98*14MM D3	US\$ 84.70	+Info
476135	Aconia SHTML 98*14MM D2	US\$ 84.70	+Info
476134	Aconia SHTML 98*14MM D1	US\$ 84.70	+Info
476133	Aconia SHTML 98*14MM C4	US\$ 84.70	+Info
476132	Aconia SHTML 98*14MM C3	US\$ 84.70	+Info
476121	Aconia SHTML 98*14MM A1	US\$ 84.70	+Info
476130	Aconia SHTML 98*14MM C1	US\$ 84.70	+Info
476129	Aconia SHTML 98*14MM B4	US\$ 84.70	+Info
476128	Aconia SHTML 98*14MM B3	US\$ 84.70	+Info
476127	Aconia SHTML 98*14MM B2	US\$ 84.70	+Info
476126	Aconia SHTML 98*14MM B1	US\$ 84.70	+Info
476125	Aconia SHTML 98*14MM A4	US\$ 84.70	+Info
476124	Aconia SHTML 98*14MM A3.5	US\$ 84.70	+Info
476123	Aconia SHTML 98*14MM A3	US\$ 84.70	+Info
476122	Aconia SHTML 98*14MM A2	US\$ 84.70	+Info

16mm

16mm

Code	Name	Price	+
476151	Aconia SHTML 98*16MM C2	US\$ 91.00	+Info
476160	Aconia SHTML 98*16MM OM3	US\$ 91.00	+Info
476159	Aconia SHTML 98*16MM OM2	US\$ 91.00	+Info
476158	Aconia SHTML 98*16MM OM1	US\$ 91.00	+Info
476157	Aconia SHTML 98*16MM D4	US\$ 91.00	+Info

476156	Aconia SHTML 98*16MM D3	US\$ 91.00	+Info
476155	Aconia SHTML 98*16MM D2	US\$ 91.00	+Info
476154	Aconia SHTML 98*16MM D1	US\$ 91.00	+Info
476153	Aconia SHTML 98*16MM C4	US\$ 91.00	+Info
476152	Aconia SHTML 98*16MM C3	US\$ 91.00	+Info
476141	Aconia SHTML 98*16MM A1	US\$ 91.00	+Info
476150	Aconia SHTML 98*16MM C1	US\$ 91.00	+Info
476149	Aconia SHTML 98*16MM B4	US\$ 91.00	+Info
476148	Aconia SHTML 98*16MM B3	US\$ 91.00	+Info
476147	Aconia SHTML 98*16MM B2	US\$ 91.00	+Info
476146	Aconia SHTML 98*16MM B1	US\$ 91.00	+Info
476145	Aconia SHTML 98*16MM A4	US\$ 91.00	+Info
476144	Aconia SHTML 98*16MM A3.5	US\$ 91.00	+Info
476143	Aconia SHTML 98*16MM A3	US\$ 91.00	+Info
476142	Aconia SHTML 98*16MM A2	US\$ 91.00	+Info

18mm

18mm

Code	Name	Price	+
476171	Aconia SHTML 98*18MM C2	US\$ 102.20	+Info
476180	Aconia SHTML 98*18MM OM3	US\$ 102.20	+Info
476179	Aconia SHTML 98*18MM OM2	US\$ 102.20	+Info
476178	Aconia SHTML 98*18MM OM1	US\$ 102.20	+Info
476177	Aconia SHTML 98*18MM D4	US\$ 102.20	+Info
476176	Aconia SHTML 98*18MM D3	US\$ 102.20	+Info
476175	Aconia SHTML 98*18MM D2	US\$ 102.20	+Info
476174	Aconia SHTML 98*18MM D1	US\$ 102.20	+Info
476173	Aconia SHTML 98*18MM C4	US\$ 102.20	+Info
476172	Aconia SHTML 98*18MM C3	US\$ 102.20	+Info
476161	Aconia SHTML 98*18MM A1	US\$ 102.20	+Info
476170	Aconia SHTML 98*18MM C1	US\$ 102.20	+Info
476169	Aconia SHTML 98*18MM B4	US\$ 102.20	+Info
476168	Aconia SHTML 98*18MM B3	US\$ 102.20	+Info
476167	Aconia SHTML 98*18MM B2	US\$ 102.20	+Info
476166	Aconia SHTML 98*18MM B1	US\$ 102.20	+Info
476165	Aconia SHTML 98*18MM A4	US\$ 102.20	+Info
476164	Aconia SHTML 98*18MM A3.5	US\$ 102.20	+Info
476163	Aconia SHTML 98*18MM A3	US\$ 102.20	+Info
476162	Aconia SHTML 98*18MM A2	US\$ 102.20	+Info

20mm

20mm

Code	Name	Price	+
476191	Aconia SHTML 98*20MM C2	US\$ 113.40	+Info
476200	Aconia SHTML 98*20MM OM3	US\$ 113.40	+Info
476199	Aconia SHTML 98*20MM OM2	US\$ 113.40	+Info
476198	Aconia SHTML 98*20MM OM1	US\$ 113.40	+Info
476197	Aconia SHTML 98*20MM D4	US\$ 113.40	+Info
476196	Aconia SHTML 98*20MM D3	US\$ 113.40	+Info

476195	Aconia SHTML 98*20MM D2	US\$ 113.40	+Info
476194	Aconia SHTML 98*20MM D1	US\$ 113.40	+Info
476193	Aconia SHTML 98*20MM C4	US\$ 113.40	+Info
476192	Aconia SHTML 98*20MM C3	US\$ 113.40	+Info
476181	Aconia SHTML 98*20MM A1	US\$ 113.40	+Info
476190	Aconia SHTML 98*20MM C1	US\$ 113.40	+Info
476189	Aconia SHTML 98*20MM B4	US\$ 113.40	+Info
476188	Aconia SHTML 98*20MM B3	US\$ 113.40	+Info
476187	Aconia SHTML 98*20MM B2	US\$ 113.40	+Info
476186	Aconia SHTML 98*20MM B1	US\$ 113.40	+Info
476185	Aconia SHTML 98*20MM A4	US\$ 113.40	+Info
476184	Aconia SHTML 98*20MM A3.5	US\$ 113.40	+Info
476183	Aconia SHTML 98*20MM A3	US\$ 113.40	+Info
476182	Aconia SHTML 98*20MM A2	US\$ 113.40	+Info

22mm

22mm

Code	Name	Price	+
476211	Aconia SHTML 98*22MM C2	US\$ 123.20	+Info
476220	Aconia SHTML 98*22MM OM3	US\$ 123.20	+Info
476219	Aconia SHTML 98*22MM OM2	US\$ 123.20	+Info
476218	Aconia SHTML 98*22MM OM1	US\$ 123.20	+Info
476217	Aconia SHTML 98*22MM D4	US\$ 123.20	+Info
476216	Aconia SHTML 98*22MM D3	US\$ 123.20	+Info
476215	Aconia SHTML 98*22MM D2	US\$ 123.20	+Info
476214	Aconia SHTML 98*22MM D2	US\$ 123.20	+Info
476213	Aconia SHTML 98*22MM C4	US\$ 123.20	+Info
476212	Aconia SHTML 98*22MM C3	US\$ 123.20	+Info
476201	Aconia SHTML 98*22MM A1	US\$ 123.20	+Info
476210	Aconia SHTML 98*22MM C1	US\$ 123.20	+Info
476209	Aconia SHTML 98*22MM B4	US\$ 123.20	+Info
476208	Aconia SHTML 98*22MM B3	US\$ 123.20	+Info
476207	Aconia SHTML 98*22MM B2	US\$ 123.20	+Info
476206	Aconia SHTML 98*22MM B1	US\$ 123.20	+Info
476205	Aconia SHTML 98*22MM A4	US\$ 123.20	+Info
476204	Aconia SHTML 98*22MM A3.5	US\$ 123.20	+Info
476203	Aconia SHTML 98*22MM A3	US\$ 123.20	+Info
476202	Aconia SHTML 98*22MM A2	US\$ 123.20	+Info

25mm

25mm

Code	Name	Price	+
476231	Aconia SHTML 98*25MM C3	US\$ 133.00	+Info
476240	Aconia SHTML 98*25MM OM3	US\$ 133.00	+Info
476239	Aconia SHTML 98*25MM OM2	US\$ 133.00	+Info
476238	Aconia SHTML 98*25MM OM1	US\$ 133.00	+Info
476237	Aconia SHTML 98*25MM D4	US\$ 133.00	+Info
476236	Aconia SHTML 98*25MM D3	US\$ 133.00	+Info
476235	Aconia SHTML 98*25MM D2	US\$ 133.00	+Info

476234	Aconia SHTML 98*25MM D2	US\$ 133.00	+Info
476233	Aconia SHTML 98*25MM C4	US\$ 133.00	+Info
476232	Aconia SHTML 98*25MM C3	US\$ 133.00	+Info
476221	Aconia SHTML 98*25MM A1	US\$ 133.00	+Info
476230	Aconia SHTML 98*25MM C1	US\$ 133.00	+Info
476229	Aconia SHTML 98*25MM B4	US\$ 133.00	+Info
476228	Aconia SHTML 98*25MM B3	US\$ 133.00	+Info
476227	Aconia SHTML 98*25MM B2	US\$ 133.00	+Info
476226	Aconia SHTML 98*25MM B1	US\$ 133.00	+Info
476225	Aconia SHTML 98*25MM A4	US\$ 133.00	+Info
476224	Aconia SHTML 98*25MM A3.5	US\$ 133.00	+Info
476223	Aconia SHTML 98*25MM A3	US\$ 133.00	+Info
476222	Aconia SHTML 98*25MM A2	US\$ 133.00	+Info

ST-PRE Super Translucent Mono Pre-Shaded 98.5mm

ST-PRE

12mm

12mm

Code	Name	Price	+
476249	Aconia ST 98*12MM B4	US\$ 49.00	+Info
476257	Aconia ST 98*12MM D4	US\$ 49.00	+Info
476256	Aconia ST 98*12MM D3	US\$ 49.00	+Info
476255	Aconia ST 98*12MM D2	US\$ 49.00	+Info
476254	Aconia ST 98*12MM D1	US\$ 49.00	+Info
476253	Aconia ST 98*12MM C4	US\$ 49.00	+Info
476252	Aconia ST 98*12MM C3	US\$ 49.00	+Info
476251	Aconia ST 98*12MM C2	US\$ 49.00	+Info
476250	Aconia ST 98*12MM C1	US\$ 49.00	+Info
476241	Aconia ST 98*12MM A1	US\$ 49.00	+Info
476248	Aconia ST 98*12MM B3	US\$ 49.00	+Info
476247	Aconia ST 98*12MM B2	US\$ 49.00	+Info
476246	Aconia ST 98*12MM B1	US\$ 49.00	+Info
476245	Aconia ST 98*12MM A4	US\$ 49.00	+Info
476244	Aconia ST 98*12MM A3.5	US\$ 49.00	+Info
476243	Aconia ST 98*12MM A3	US\$ 49.00	+Info
476242	Aconia ST 98*12MM A2	US\$ 49.00	+Info

14mm

14mm

Code	Name	Price	+
476266	Aconia ST 98*14MM B4	US\$ 61.60	+Info
476274	Aconia ST 98*14MM D4	US\$ 77.00	+Info
476273	Aconia ST 98*14MM D3	US\$ 61.60	+Info
476272	Aconia ST 98*14MM D2	US\$ 61.60	+Info
476271	Aconia ST 98*14MM D1	US\$ 61.60	+Info
476270	Aconia ST 98*14MM C4	US\$ 61.60	+Info
476269	Aconia ST 98*14MM C3	US\$ 61.60	+Info

476268	Aconia ST 98*14MM C2	US\$ 61.60	+Info
476267	Aconia ST 98*14MM C1	US\$ 61.60	+Info
476258	Aconia ST 98*14MM A1	US\$ 61.60	+Info
476265	Aconia ST 98*14MM B3	US\$ 61.60	+Info
476264	Aconia ST 98*14MM B2	US\$ 61.60	+Info
476263	Aconia ST 98*14MM B1	US\$ 61.60	+Info
476262	Aconia ST 98*14MM A4	US\$ 61.60	+Info
476261	Aconia ST 98*14MM A3.5	US\$ 61.60	+Info
476260	Aconia ST 98*14MM A3	US\$ 61.60	+Info
476259	Aconia ST 98*14MM A2	US\$ 61.60	+Info

16mm

16mm

Code	Name	Price	+
476283	Aconia ST 98*16MM B4	US\$ 65.80	+Info
476291	Aconia ST 98*16MM D4	US\$ 65.80	+Info
476290	Aconia ST 98*16MM D3	US\$ 65.80	+Info
476289	Aconia ST 98*16MM D2	US\$ 65.80	+Info
476288	Aconia ST 98*16MM D1	US\$ 65.80	+Info
476287	Aconia ST 98*16MM C4	US\$ 65.80	+Info
476286	Aconia ST 98*16MM C3	US\$ 65.80	+Info
476285	Aconia ST 98*16MM C2	US\$ 65.80	+Info
476284	Aconia ST 98*16MM C1	US\$ 65.80	+Info
476275	Aconia ST 98*16MM A1	US\$ 65.80	+Info
476282	Aconia ST 98*16MM B3	US\$ 65.80	+Info
476281	Aconia ST 98*16MM B2	US\$ 65.80	+Info
476280	Aconia ST 98*16MM B1	US\$ 65.80	+Info
476279	Aconia ST 98*16MM A4	US\$ 65.80	+Info
476278	Aconia ST 98*16MM A3.5	US\$ 65.80	+Info
476277	Aconia ST 98*16MM A3	US\$ 65.80	+Info
476276	Aconia ST 98*16MM A2	US\$ 65.80	+Info

18mm

18mm

Code	Name	Price	+
476300	Aconia ST 98*18MM B4	US\$ 78.40	+Info
476308	Aconia ST 98*18MM D4	US\$ 78.40	+Info
476307	Aconia ST 98*18MM D3	US\$ 78.40	+Info
476306	Aconia ST 98*18MM D4	US\$ 78.40	+Info
476305	Aconia ST 98*18MM D3	US\$ 78.40	+Info
476304	Aconia ST 98*18MM C4	US\$ 78.40	+Info
476303	Aconia ST 98*18MM C3	US\$ 78.40	+Info
476302	Aconia ST 98*18MM C2	US\$ 78.40	+Info
476301	Aconia ST 98*18MM C1	US\$ 78.40	+Info
476292	Aconia ST 98*18MM A1	US\$ 78.40	+Info
476299	Aconia ST 98*18MM B3	US\$ 78.40	+Info
476298	Aconia ST 98*18MM B2	US\$ 78.40	+Info
476297	Aconia ST 98*18MM B1	US\$ 78.40	+Info
476296	Aconia ST 98*18MM A4	US\$ 78.40	+Info

476295	Aconia ST 98*18MM A3.5	US\$ 78.40	+Info
476294	Aconia ST 98*18MM A3	US\$ 78.40	+Info
476293	Aconia ST 98*18MM A2	US\$ 78.40	+Info

20mm

20mm

Code	Name	Price	+
476317	Aconia ST 98*20MM B4	US\$ 88.20	+Info
476325	Aconia ST 98*20MM D4	US\$ 88.20	+Info
476324	Aconia ST 98*20MM D3	US\$ 88.20	+Info
476323	Aconia ST 98*20MM D2	US\$ 88.20	+Info
476322	Aconia ST 98*20MM D1	US\$ 88.20	+Info
476321	Aconia ST 98*20MM C4	US\$ 88.20	+Info
476320	Aconia ST 98*20MM C3	US\$ 88.20	+Info
476319	Aconia ST 98*20MM C2	US\$ 88.20	+Info
476318	Aconia ST 98*20MM C1	US\$ 88.20	+Info
476309	Aconia ST 98*20MM A1	US\$ 88.20	+Info
476316	Aconia ST 98*20MM B3	US\$ 88.20	+Info
476315	Aconia ST 98*20MM B2	US\$ 88.20	+Info
476314	Aconia ST 98*20MM B1	US\$ 88.20	+Info
476313	Aconia ST 98*20MM A4	US\$ 88.20	+Info
476312	Aconia ST 98*20MM A3	US\$ 88.20	+Info
476311	Aconia ST 98*20MM A3.5	US\$ 88.20	+Info
476310	Aconia ST 98*20MM A2	US\$ 88.20	+Info

22mm

22mm

Code	Name	Price	+
476334	Aconia ST 98*22MM B4	US\$ 98.00	+Info
476342	Aconia ST 98*22MM D4	US\$ 98.00	+Info
476341	Aconia ST 98*22MM D3	US\$ 98.00	+Info
476340	Aconia ST 98*22MM D2	US\$ 98.00	+Info
476339	Aconia ST 98*22MM D1	US\$ 98.00	+Info
476338	Aconia ST 98*22MM C4	US\$ 98.00	+Info
476337	Aconia ST 98*22MM C3	US\$ 98.00	+Info
476336	Aconia ST 98*22MM C2	US\$ 98.00	+Info
476335	Aconia ST 98*22MM C1	US\$ 98.00	+Info
476326	Aconia ST 98*22MM A1	US\$ 98.00	+Info
476333	Aconia ST 98*22MM B3	US\$ 98.00	+Info
476332	Aconia ST 98*22MM B2	US\$ 98.00	+Info
476331	Aconia ST 98*22MM B1	US\$ 98.00	+Info
476330	Aconia ST 98*22MM A4	US\$ 98.00	+Info
476329	Aconia ST 98*22MM A3.5	US\$ 98.00	+Info
476328	Aconia ST 98*22MM A3	US\$ 98.00	+Info
476327	Aconia ST 98*22MM A2	US\$ 98.00	+Info

25mm

25mm

Code	Name	Price	+
476351	Aconia ST 98*25MM B4	US\$ 106.40	+Info
476359	Aconia ST 98*25MM D4	US\$ 106.40	+Info
476358	Aconia ST 98*25MM D3	US\$ 106.40	+Info
476357	Aconia ST 98*25MM D2	US\$ 106.40	+Info
476356	Aconia ST 98*25MM D1	US\$ 106.40	+Info
476355	Aconia ST 98*25MM C4	US\$ 106.40	+Info
476354	Aconia ST 98*25MM C3	US\$ 106.40	+Info
476353	Aconia ST 98*25MM C2	US\$ 106.40	+Info
476352	Aconia ST 98*25MM C1	US\$ 106.40	+Info
476343	Aconia ST 98*25MM A1	US\$ 106.40	+Info
476350	Aconia ST 98*25MM B3	US\$ 106.40	+Info
476349	Aconia ST 98*25MM B2	US\$ 106.40	+Info
476348	Aconia ST 98*25MM B1	US\$ 106.40	+Info
476347	Aconia ST 98*25MM A4	US\$ 106.40	+Info
476346	Aconia ST 98*25MM A3.5	US\$ 106.40	+Info
476345	Aconia ST 98*25MM A3	US\$ 106.40	+Info
476344	Aconia ST 98*25MM A2	US\$ 106.40	+Info

TTML Total Translucent Multi Layer 98.5mm

TTML

14mm

14mm

Code	Name	Price	+
476360	Aconia TTML 98*14MM A1	US\$ 107.80	+Info
476361	Aconia TTML 98*14MM A2	US\$ 107.80	+Info
476362	Aconia TTML 98*14MM A3	US\$ 107.80	+Info
476363	Aconia TTML 98*14MM A4	US\$ 107.80	+Info
476364	Aconia TTML 98*14MM B1	US\$ 107.80	+Info
476365	Aconia TTML 98*14MM B2	US\$ 107.80	+Info
476366	Aconia TTML 98*14MM B3	US\$ 107.80	+Info
476367	Aconia TTML 98*14MM B4	US\$ 107.80	+Info
476368	Aconia TTML 98*14MM C1	US\$ 107.80	+Info
476369	Aconia TTML 98*14MM C2	US\$ 107.80	+Info
476370	Aconia TTML 98*14MM C3	US\$ 107.80	+Info
476371	Aconia TTML 98*14MM C4	US\$ 107.80	+Info
476372	Aconia TTML 98*14MM D1	US\$ 107.80	+Info
476373	Aconia TTML 98*14MM D2	US\$ 107.80	+Info
476374	Aconia TTML 98*14MM D3	US\$ 107.80	+Info
476375	Aconia TTML 98*14MM D4	US\$ 107.80	+Info

16mm

16mm

Code	Name	Price	+
476376	Aconia TTML 98*16MM A1	US\$ 112.00	+Info
476377	Aconia TTML 98*16MM A2	US\$ 112.00	+Info
476378	Aconia TTML 98*16MM A3	US\$ 112.00	+Info
476379	Aconia TTML 98*16MM A4	US\$ 112.00	+Info

476380	Aconia TTML 98*16MM B1	US\$ 112.00	+Info
476381	Aconia TTML 98*16MM B2	US\$ 112.00	+Info
476382	Aconia TTML 98*16MM B3	US\$ 112.00	+Info
476383	Aconia TTML 98*16MM B4	US\$ 112.00	+Info
476384	Aconia TTML 98*16MM C1	US\$ 112.00	+Info
476385	Aconia TTML 98*16MM C2	US\$ 112.00	+Info
476386	Aconia TTML 98*16MM C3	US\$ 112.00	+Info
476387	Aconia TTML 98*16MM C4	US\$ 112.00	+Info
476388	Aconia TTML 98*16MM D1	US\$ 112.00	+Info
476389	Aconia TTML 98*16MM D2	US\$ 112.00	+Info
476390	Aconia TTML 98*16MM D3	US\$ 112.00	+Info
476391	Aconia TTML 98*16MM D4	US\$ 112.00	+Info

18mm
18mm

Code	Name	Price	+
476392	Aconia TTML 98*18MM A1	US\$ 126.00	+Info
476393	Aconia TTML 98*18MM A2	US\$ 126.00	+Info
476394	Aconia TTML 98*18MM A3	US\$ 126.00	+Info
476395	Aconia TTML 98*18MM A4	US\$ 126.00	+Info
476396	Aconia TTML 98*18MM B1	US\$ 126.00	+Info
476397	Aconia TTML 98*18MM B2	US\$ 126.00	+Info
476398	Aconia TTML 98*18MM B3	US\$ 126.00	+Info
476399	Aconia TTML 98*18MM B4	US\$ 126.00	+Info
476400	Aconia TTML 98*18MM C1	US\$ 126.00	+Info
476401	Aconia TTML 98*18MM C2	US\$ 126.00	+Info
476402	Aconia TTML 98*18MM C3	US\$ 126.00	+Info
476403	Aconia TTML 98*18MM C4	US\$ 126.00	+Info
476404	Aconia TTML 98*18MM D1	US\$ 126.00	+Info
476405	Aconia TTML 98*18MM D2	US\$ 126.00	+Info
476406	Aconia TTML 98*18MM D3	US\$ 126.00	+Info
476407	Aconia TTML 98*18MM D4	US\$ 126.00	+Info

ZIRCONIA DISCS | AMANN GIRRBACH
H



MONO DISCS AMANN GIRRBACH

ZOLID FX WHITE monolayer Zirconia



AMANN GIRRACH

Code	Name	Price	+
468156	Ceramill ZOLID FX 98x12 N	US\$ 203.86	+Info
468157	Ceramill ZOLID FX 98x14 N	US\$ 203.86	+Info
468158	Ceramill ZOLID FX 98x16 N	US\$ 239.55	+Info
468159	Ceramill ZOLID FX 98x18 N	US\$ 240.81	+Info
468160	Ceramill ZOLID FX 98x20 N	US\$ 279.03	+Info

ZOLID HT+ monolayer Zirconia



AMANN GIRRACH

HT+ White 71mm

Code	Name	Price	+
468578	Ceramill Zolid HT+ white 71x14	US\$ 178.38	+Info
468579	Ceramill Zolid HT+ white 71x16	US\$ 172.93	+Info
468581	Ceramill Zolid HT+ white 71x20	US\$ 226.79	+Info
468582	Ceramill Zolid HT+ white 71x25	US\$ 265.01	+Info

HT+ WHITE 98mm

Code	Name	Price	+
468583	Ceramill ZOLID HT+ white 98x10	US\$ 152.90	+Info
468584	Ceramill ZOLID HT+ white 98x12	US\$ 178.38	+Info
468585	Ceramill ZOLID HT+ white 98x14	US\$ 183.47	+Info
468586	Ceramill ZOLID HT+ white 98x16	US\$ 202.58	+Info
468587	Ceramill ZOLID HT+ white 98x18	US\$ 215.33	+Info
468588	Ceramill ZOLID HT+ white 98x20	US\$ 228.07	+Info
468589	Ceramill ZOLID HT+ white 98x25	US\$ 263.75	+Info

ZOLID HT+ PreShaded 71mm

Code	Name	Price	+
468591	Ceramill Zolid HT+ PS BL 71x16	US\$ 198.77	+Info
468593	Ceramill Zolid HT+ PS A1 71x14	US\$ 179.66	+Info

468594	Ceramill Zolid HT+ PS A1 71x16	US\$ 198.77	+Info
468596	Ceramill Zolid HT+ PS A2 71x14	US\$ 179.66	+Info
468597	Ceramill Zolid HT+ PS A2 71x16	US\$ 198.77	+Info
468600	Ceramill Zolid HT+ PS A3 71x16	US\$ 198.77	+Info
468608	Ceramill Zolid HT+ PS B1 71x14	US\$ 179.66	+Info
468609	Ceramill Zolid HT+ PS B1 71x16	US\$ 198.77	+Info
468612	Ceramill Zolid HT+ PS B2 71x16	US\$ 198.77	+Info
468621	Ceramill Zolid HT+ PS C1 71x16	US\$ 198.77	+Info
468624	Ceramill Zolid HT+ PS C2 71x16	US\$ 198.77	+Info
468633	Ceramill Zolid HT+ PS D2 71x16	US\$ 198.77	+Info
468634	Ceramill Zolid HT+ PS D2 71x20	US\$ 226.79	+Info
468671	Ceramill ZOLID HT+ PS C1 98x14	US\$ 184.74	+Info

ZOLID HT+ PreShaded 98

Code	Name	Price	+
468680	Ceramill ZOLID HT+ PS C4 98x14	US\$ 184.74	+Info
468667	Ceramill ZOLID HT+ PS B3 98x20	US\$ 228.07	+Info
468668	Ceramill ZOLID HT+ PS B4 98x14	US\$ 184.74	+Info
468669	Ceramill ZOLID HT+ PS B4 98x16	US\$ 202.58	+Info
468670	Ceramill ZOLID HT+ PS B4 98x20	US\$ 228.07	+Info
468672	Ceramill ZOLID HT+ PS C1 98x16	US\$ 202.58	+Info
468673	Ceramill ZOLID HT+ PS C1 98x20	US\$ 228.07	+Info
468674	Ceramill ZOLID HT+ PS C2 98x14	US\$ 184.74	+Info
468675	Ceramill ZOLID HT+ PS C2 98x16	US\$ 202.58	+Info
468676	Ceramill ZOLID HT+ PS C2 98x20	US\$ 228.07	+Info
468677	Ceramill ZOLID HT+ PS C3 98x14	US\$ 184.74	+Info
468678	Ceramill ZOLID HT+ PS C3 98x16	US\$ 202.58	+Info
468679	Ceramill ZOLID HT+ PS C3 98x20	US\$ 228.07	+Info
468666	Ceramill ZOLID HT+ PS B3 98x16	US\$ 202.58	+Info
468681	Ceramill ZOLID HT+ PS C4 98x16	US\$ 202.58	+Info
468682	Ceramill ZOLID HT+ PS C4 98x20	US\$ 228.07	+Info
468683	Ceramill ZOLID HT+ PS D2 98x14	US\$ 184.74	+Info
468684	Ceramill ZOLID HT+ PS D2 98x16	US\$ 202.58	+Info
468685	Ceramill ZOLID HT+ PS D2 98x20	US\$ 228.07	+Info
468686	Ceramill ZOLID HT+ PS D3 98x14	US\$ 184.74	+Info
468687	Ceramill ZOLID HT+ PS D3 98x16	US\$ 202.58	+Info
468688	Ceramill ZOLID HT+ PS D3 98x20	US\$ 228.07	+Info
468689	Ceramill ZOLID HT+ PS D4 98x14	US\$ 184.74	+Info
468690	Ceramill ZOLID HT+ PS D4 98x16	US\$ 202.58	+Info
468691	Ceramill ZOLID HT+ PS D4 98x20	US\$ 228.07	+Info
468654	Ceramill ZOLID HT+ PS A3,5 98x16	US\$ 202.58	+Info
468642	Ceramill ZOLID HT+ PS BL 98x16	US\$ 202.58	+Info
468643	Ceramill ZOLID HT+ PS BL 98x20	US\$ 228.07	+Info
468644	Ceramill ZOLID HT+ PS A1 98x14	US\$ 184.74	+Info
468645	Ceramill ZOLID HT+ PS A1 98x16	US\$ 202.58	+Info
468646	Ceramill ZOLID HT+ PS A1 98x20	US\$ 168.77	+Info
468647	Ceramill ZOLID HT+ PS A2 98x14	US\$ 184.74	+Info
468648	Ceramill ZOLID HT+ PS A2 98x16	US\$ 149.91	+Info
468649	Ceramill ZOLID HT+ PS A2 98x20	US\$ 228.07	+Info

468650	Ceramill ZOLID HT+ PS A3 98x14	US\$ 184.74	+Info
468651	Ceramill ZOLID HT+ PS A3 98x16	US\$ 202.58	+Info
468652	Ceramill ZOLID HT+ PS A3 98x20	US\$ 228.07	+Info
468653	Ceramill ZOLID HT+ PS A3,5 98x14	US\$ 184.74	+Info
468641	Ceramill ZOLID HT+ PS BL 98x14	US\$ 184.74	+Info
468655	Ceramill ZOLID HT+ PS A3,5 98x20	US\$ 228.07	+Info
468656	Ceramill ZOLID HT+ PS A4 98x14	US\$ 184.74	+Info
468657	Ceramill ZOLID HT+ PS A4 98x16	US\$ 202.58	+Info
468658	Ceramill ZOLID HT+ PS A4 98x20	US\$ 228.07	+Info
468659	Ceramill ZOLID HT+ PS B1 98x14	US\$ 184.74	+Info
468660	Ceramill ZOLID HT+ PS B1 98x16	US\$ 202.58	+Info
468661	Ceramill ZOLID HT+ PS B1 98x20	US\$ 228.07	+Info
468662	Ceramill ZOLID HT+ PS B2 98x14	US\$ 184.74	+Info
468663	Ceramill ZOLID HT+ PS B2 98x16	US\$ 202.58	+Info
468664	Ceramill ZOLID HT+ PS B2 98x20	US\$ 228.07	+Info
468665	Ceramill ZOLID HT+ PS B3 98x14	US\$ 184.74	+Info

ZOLID HT+ PS

ZOLID HT+ PS

Code	Name	Price	+
468630	Ceramill Zolid HT+ PS C4 71x16	US\$ 198.77	+Info
468619	Ceramill Zolid HT+ PS B4 71x20	US\$ 226.79	+Info
468620	Ceramill Zolid HT+ PS C1 71x14	US\$ 179.66	+Info
468622	Ceramill Zolid HT+ PS C1 71x20	US\$ 226.79	+Info
468623	Ceramill Zolid HT+ PS C2 71x14	US\$ 179.66	+Info
468625	Ceramill Zolid HT+ PS C2 71x20	US\$ 226.79	+Info
468626	Ceramill Zolid HT+ PS C3 71x14	US\$ 179.66	+Info
468627	Ceramill Zolid HT+ PS C3 71x16	US\$ 198.77	+Info
468628	Ceramill Zolid HT+ PS C3 71x20	US\$ 226.79	+Info
468629	Ceramill Zolid HT+ PS C4 71x14	US\$ 179.66	+Info
468618	Ceramill Zolid HT+ PS B4 71x16	US\$ 198.77	+Info
468631	Ceramill Zolid HT+ PS C4 71x20	US\$ 226.79	+Info
468632	Ceramill Zolid HT+ PS D2 71x14	US\$ 179.66	+Info
468635	Ceramill Zolid HT+ PS D3 71x14	US\$ 179.66	+Info
468636	Ceramill Zolid HT+ PS D3 71x16	US\$ 198.77	+Info
468637	Ceramill Zolid HT+ PS D3 71x20	US\$ 226.79	+Info
468638	Ceramill Zolid HT+ PS D4 71x14	US\$ 179.66	+Info
468639	Ceramill Zolid HT+ PS D4 71x16	US\$ 198.77	+Info
468640	Ceramill Zolid HT+ PS D4 71x20	US\$ 226.79	+Info
468606	Ceramill Zolid HT+ PS A4 71x16	US\$ 198.77	+Info
468592	Ceramill Zolid HT+ PS BL 71x20	US\$ 226.79	+Info
468595	Ceramill Zolid HT+ PS A1 71x20	US\$ 226.79	+Info
468598	Ceramill Zolid HT+ PS A2 71x20	US\$ 226.79	+Info
468599	Ceramill Zolid HT+ PS A3 71x14	US\$ 179.66	+Info
468601	Ceramill Zolid HT+ PS A3 71x20	US\$ 226.79	+Info
468602	Ceramill Zolid HT+ PS A3,5 71x14	US\$ 179.66	+Info
468603	Ceramill Zolid HT+ PS A3,5 71x16	US\$ 198.77	+Info
468604	Ceramill Zolid HT+ PS A3,5 71x20	US\$ 226.79	+Info
468605	Ceramill Zolid HT+ PS A4 71x14	US\$ 179.66	+Info

468590	Ceramill Zolid HT+ PS BL 71x14	US\$ 179.66	+Info
468607	Ceramill Zolid HT+ PS A4 71x20	US\$ 226.79	+Info
468610	Ceramill Zolid HT+ PS B1 71x20	US\$ 226.79	+Info
468611	Ceramill Zolid HT+ PS B2 71x14	US\$ 179.66	+Info
468613	Ceramill Zolid HT+ PS B2 71x20	US\$ 226.79	+Info
468614	Ceramill Zolid HT+ PS B3 71x14	US\$ 179.66	+Info
468615	Ceramill Zolid HT+ PS B3 71x16	US\$ 198.77	+Info
468616	Ceramill Zolid HT+ PS B3 71x20	US\$ 226.79	+Info
468617	Ceramill Zolid HT+ PS B4 71x14	US\$ 179.66	+Info

ZOLID HT+ White

ZOLID HT+ White

Code	Name	Price	+
468576	Ceramill Zolid HT+ white 71x10	US\$ 145.24	+Info
468577	Ceramill Zolid HT+ white 71x12	US\$ 172.01	+Info
468580	Ceramill Zolid HT+ white 71x18	US\$ 212.79	+Info

MULTILAYER DISCS AMANN GIRRBACH

ZOLID FX MULTILAYER - AMANN G IRRBACH



AMANN GIRRBACH

Code	Name	Price	+
468481	Ceramill ZOLID FX ML B3/B4 98x14 N	US\$ 229.34	+Info
468495	Ceramill ZOLID FX ML D3/D4 98x20 N	US\$ 293.05	+Info
468494	Ceramill ZOLID FX ML D3/D4 98x16 N	US\$ 254.83	+Info
468493	Ceramill ZOLID FX ML D3/D4 98x14 N	US\$ 229.34	+Info
468492	Ceramill ZOLID FX ML D2/D3 98x20 N	US\$ 293.05	+Info
468491	Ceramill ZOLID FX ML D2/D3 98x16 N	US\$ 254.83	+Info
468490	Ceramill ZOLID FX ML D2/D3 98x14 N	US\$ 229.34	+Info
468489	Ceramill ZOLID FX ML C3/C4 98x20 N	US\$ 293.05	+Info
468488	Ceramill ZOLID FX ML C3/C4 98x16 N	US\$ 254.83	+Info
468487	Ceramill ZOLID FX ML C3/C4 98x14 N	US\$ 229.34	+Info
468486	Ceramill ZOLID FX ML C1/C2 98x20 N	US\$ 293.05	+Info
468485	Ceramill ZOLID FX ML C1/C2 98x16 N	US\$ 254.83	+Info
468484	Ceramill ZOLID FX ML C1/C2 98x14 N	US\$ 229.34	+Info
468483	Ceramill ZOLID FX ML B3/B4 98x20 N	US\$ 293.05	+Info
468482	Ceramill ZOLID FX ML B3/B4 98x16 N	US\$ 254.83	+Info
468466	Ceramill ZOLID FX ML 0/A1 98x14 N	US\$ 229.34	+Info
468480	Ceramill ZOLID FX ML B2/B3 98x20 N	US\$ 293.05	+Info
468479	Ceramill ZOLID FX ML B2/B3 98x16 N	US\$ 254.83	+Info
468478	Ceramill ZOLID FX ML B2/B3 98x14 N	US\$ 229.34	+Info
468477	Ceramill ZOLID FX ML 0/B1 98x20 N	US\$ 293.05	+Info

468476	Ceramill ZOLID FX ML 0/B1 98x16 N	US\$ 254.83	+Info
468475	Ceramill ZOLID FX ML 0/ B1 98x14 N	US\$ 229.34	+Info
468474	Ceramill ZOLID FX ML A3,5/A4 98x20	US\$ 293.05	+Info
468473	Ceramill ZOLID FX ML A3,5/A4 98x16	US\$ 254.83	+Info
468472	Ceramill ZOLID FX ML A3,5/A4 98x14	US\$ 229.34	+Info
468471	Ceramill ZOLID FX ML A2/A3 98x20 N	US\$ 293.05	+Info
468470	Ceramill ZOLID FX ML A2/A3 98x16 N	US\$ 254.83	+Info
468469	Ceramill ZOLID FX ML A2/A3 98x14 N	US\$ 229.34	+Info
468468	Ceramill ZOLID FX ML 0/A1 98x20 N	US\$ 293.05	+Info
468467	Ceramill ZOLID FX ML 0/A1 98x16 N	US\$ 254.83	+Info

ZOLID GEN+ GRADIENT - AMANN GIRRBACH



AMANNGIRRBACH

ZOLID GEN+ ZOLID GEN+

Code	Name	Price	+
468812	ZOLID GEN-X C3 98x16	US\$ 267.57	+Info
468773	ZOLID GEN-X B1 98x22	US\$ 331.27	+Info
468774	ZOLID GEN-X B1 98x25	US\$ 356.76	+Info
468777	ZOLID GEN-X B2 98x16	US\$ 267.57	+Info
468779	ZOLID GEN-X B2 98x20	US\$ 305.79	+Info
468786	ZOLID GEN-X B3 98x20	US\$ 305.79	+Info
468797	ZOLID GEN-X C1 98x14	US\$ 229.34	+Info
468798	ZOLID GEN-X C1 98x16	US\$ 267.57	+Info
468800	ZOLID GEN-X C1 98x20	US\$ 305.79	+Info
468805	ZOLID GEN-X C2 98x16	US\$ 267.57	+Info
468807	ZOLID GEN-X C2 98x20	US\$ 305.79	+Info
468811	ZOLID GEN-X C3 98x14	US\$ 229.34	+Info
468772	ZOLID GEN-X B1 98x20	US\$ 305.79	+Info
468814	ZOLID GEN-X C3 98x20	US\$ 305.79	+Info
468821	ZOLID GEN-X C4 98x20	US\$ 305.79	+Info
468825	ZOLID GEN-X D2 98x14	US\$ 229.34	+Info
468826	ZOLID GEN-X D2 98x16	US\$ 267.57	+Info
468828	ZOLID GEN-X D2 98x20	US\$ 305.79	+Info
468832	ZOLID GEN-X D3 98x14	US\$ 229.34	+Info
468833	ZOLID GEN-X D3 98x16	US\$ 267.57	+Info
468835	ZOLID GEN-X D3 98x20	US\$ 305.79	+Info
468839	ZOLID GEN-X D4 98x14	US\$ 229.34	+Info
468840	ZOLID GEN-X D4 98x16	US\$ 267.57	+Info
468842	ZOLID GEN-X D4 98x20	US\$ 305.79	+Info
468746	ZOLID GEN-X A2 98x25	US\$ 356.76	+Info

468723	ZOLID GEN-X BL1 98x20	US\$ 305.79	+Info
468725	ZOLID GEN-X BL1 98x25	US\$ 356.76	+Info
468727	ZOLID GEN-X BL3 98x14	US\$ 229.34	+Info
468728	ZOLID GEN-X BL3 98x16	US\$ 238.14	+Info
468730	ZOLID GEN-X BL3 98x20	US\$ 305.79	+Info
468734	ZOLID GEN-X A1 98x14	US\$ 229.34	+Info
468735	ZOLID GEN-X A1 98x16	US\$ 267.57	+Info
468737	ZOLID GEN-X A1 98x20	US\$ 305.79	+Info
468741	ZOLID GEN-X A2 98x14	US\$ 229.34	+Info
468742	ZOLID GEN-X A2 98x16	US\$ 238.14	+Info
468744	ZOLID GEN-X A2 98x20	US\$ 305.79	+Info
468721	ZOLID GEN-X BL1 98x16	US\$ 267.57	+Info
468748	ZOLID GEN-X A3 98x14	US\$ 229.34	+Info
468749	ZOLID GEN-X A3 98x16	US\$ 267.57	+Info
468751	ZOLID GEN-X A3 98x20	US\$ 305.79	+Info
468756	ZOLID GEN-X A3,5 98x16	US\$ 267.57	+Info
468758	ZOLID GEN-X A3,5 98x20	US\$ 305.79	+Info
468762	ZOLID GEN-X A4 98x14	US\$ 229.34	+Info
468763	ZOLID GEN-X A4 98x16	US\$ 267.57	+Info
468765	ZOLID GEN-X A4 98x20	US\$ 305.79	+Info
468769	ZOLID GEN-X B1 98x14	US\$ 229.34	+Info
468770	ZOLID GEN-X B1 98x16	US\$ 267.57	+Info

ZOLID GEN-X

- ZOLID GEN-X

Code	Name	Price	+
468799	ZOLID GEN-X C1 98x18	US\$ 280.31	+Info
468813	ZOLID GEN-X C3 98x18	US\$ 280.31	+Info
468810	ZOLID GEN-X C3 98x12	US\$ 216.60	+Info
468809	ZOLID GEN-X C2 98x25	US\$ 356.76	+Info
468808	ZOLID GEN-X C2 98x22	US\$ 331.27	+Info
468806	ZOLID GEN-X C2 98x18	US\$ 280.31	+Info
468804	ZOLID GEN-X C2 98x14	US\$ 229.34	+Info
468803	ZOLID GEN-X C2 98x12	US\$ 216.60	+Info
468802	ZOLID GEN-X C1 98x25	US\$ 356.76	+Info
468801	ZOLID GEN-X C1 98x22	US\$ 331.27	+Info
468815	ZOLID GEN-X C3 98x22	US\$ 331.27	+Info
468796	ZOLID GEN-X C1 98x12	US\$ 216.60	+Info
468795	ZOLID GEN-X B4 98x25	US\$ 356.76	+Info
468794	ZOLID GEN-X B4 98x22	US\$ 331.27	+Info
468793	ZOLID GEN-X B4 98x20	US\$ 305.79	+Info
468792	ZOLID GEN-X B4 98x18	US\$ 280.31	+Info
468791	ZOLID GEN-X B4 98x16	US\$ 267.57	+Info
468790	ZOLID GEN-X B4 98x14	US\$ 229.34	+Info
468789	ZOLID GEN-X B4 98x12	US\$ 216.60	+Info
468788	ZOLID GEN-X B3 98x25	US\$ 356.76	+Info

468829	ZOLID GEN-X D2 98x22	US\$ 331.27	+Info
468844	ZOLID GEN-X D4 98x25	US\$ 356.76	+Info
468843	ZOLID GEN-X D4 98x22	US\$ 331.27	+Info
468841	ZOLID GEN-X D4 98x18	US\$ 280.31	+Info
468838	ZOLID GEN-X D4 98x12	US\$ 216.60	+Info
468837	ZOLID GEN-X D3 98x25	US\$ 356.76	+Info
468836	ZOLID GEN-X D3 98x22	US\$ 331.27	+Info
468834	ZOLID GEN-X D3 98x18	US\$ 280.31	+Info
468831	ZOLID GEN-X D3 98x12	US\$ 216.60	+Info
468830	ZOLID GEN-X D2 98x25	US\$ 356.76	+Info
468787	ZOLID GEN-X B3 98x22	US\$ 331.27	+Info
468827	ZOLID GEN-X D2 98x18	US\$ 280.31	+Info
468824	ZOLID GEN-X D2 98x12	US\$ 216.60	+Info
468823	ZOLID GEN-X C4 98x25	US\$ 356.76	+Info
468822	ZOLID GEN-X C4 98x22	US\$ 331.27	+Info
468820	ZOLID GEN-X C4 98x18	US\$ 280.31	+Info
468819	ZOLID GEN-X C4 98x16	US\$ 267.57	+Info
468818	ZOLID GEN-X C4 98x14	US\$ 229.34	+Info
468817	ZOLID GEN-X C4 98x12	US\$ 216.60	+Info
468816	ZOLID GEN-X C3 98x25	US\$ 356.76	+Info
468738	ZOLID GEN-X A1 98x22	US\$ 331.27	+Info
468754	ZOLID GEN-X A3,5 98x12	US\$ 216.60	+Info
468753	ZOLID GEN-X A3 98x25	US\$ 356.76	+Info
468752	ZOLID GEN-X A3 98x22	US\$ 331.27	+Info
468750	ZOLID GEN-X A3 98x18	US\$ 280.31	+Info
468747	ZOLID GEN-X A3 98x12	US\$ 216.60	+Info
468745	ZOLID GEN-X A2 98x22	US\$ 331.27	+Info
468743	ZOLID GEN-X A2 98x18	US\$ 280.31	+Info
468740	ZOLID GEN-X A2 98x12	US\$ 216.60	+Info
468739	ZOLID GEN-X A1 98x25	US\$ 356.76	+Info
468755	ZOLID GEN-X A3,5 98x14	US\$ 229.34	+Info
468736	ZOLID GEN-X A1 98x18	US\$ 280.31	+Info
468733	ZOLID GEN-X A1 98x12	US\$ 216.60	+Info
468732	ZOLID GEN-X BL3 98x25	US\$ 356.76	+Info
468731	ZOLID GEN-X BL3 98x22	US\$ 331.27	+Info
468729	ZOLID GEN-X BL3 98x18	US\$ 280.31	+Info
468726	ZOLID GEN-X BL3 98x12	US\$ 216.60	+Info
468724	ZOLID GEN-X BL1 98x22	US\$ 331.27	+Info
468722	ZOLID GEN-X BL1 98x18	US\$ 280.31	+Info
468720	ZOLID GEN-X BL1 98x14	US\$ 229.34	+Info
468771	ZOLID GEN-X B1 98x18	US\$ 280.31	+Info
468785	ZOLID GEN-X B3 98x18	US\$ 280.31	+Info
468784	ZOLID GEN-X B3 98x16	US\$ 267.57	+Info
468783	ZOLID GEN-X B3 98x14	US\$ 229.34	+Info
468782	ZOLID GEN-X B3 98x12	US\$ 216.60	+Info
468781	ZOLID GEN-X B2 98x25	US\$ 356.76	+Info
468780	ZOLID GEN-X B2 98x22	US\$ 331.27	+Info
468778	ZOLID GEN-X B2 98x18	US\$ 280.31	+Info
468776	ZOLID GEN-X B2 98x14	US\$ 229.34	+Info

468775	ZOLID GEN-X B2 98x12	US\$ 216.60	+Info
468719	ZOLID GEN-X BL1 98x12	US\$ 216.60	+Info
468768	ZOLID GEN-X B1 98x12	US\$ 216.60	+Info
468767	ZOLID GEN-X A4 98x25	US\$ 356.76	+Info
468766	ZOLID GEN-X A4 98x22	US\$ 331.27	+Info
468764	ZOLID GEN-X A4 98x18	US\$ 280.31	+Info
468761	ZOLID GEN-X A4 98x12	US\$ 216.60	+Info
468760	ZOLID GEN-X A3,5 98x25	US\$ 356.76	+Info
468759	ZOLID GEN-X A3,5 98x22	US\$ 331.27	+Info
468757	ZOLID GEN-X A3,5 98x18	US\$ 280.31	+Info

ZIRCONIA BLOCKS



ZIRCONIA BLOCKS | 3M ESPE



Zirconia Blocks Cerec - 3M ESPE



Code	Name	Price	+
428368	Chairside Zirconia A3 - 39MM for CE	US\$ 162.50	+Info
428370	Chairside Zirconia A1 - 39MM for CE	US\$ 195.44	+Info
428373	Chairside Zirconia D2 - 39MM for CE	US\$ 143.00	+Info
428374	Chairside Zirconia BL - 39MM for CE	US\$ 222.94	+Info
428375	Chairside Zirconia B1 - 39MM for CE	US\$ 195.44	+Info
428378	Chairside Zirconia C1 - 39MM for CE	US\$ 222.94	+Info
428381	Chairside Zirconia A3.5 39MM 5pk	US\$ 156.00	+Info
428383	Chairside Zirconia A2 - 39MM for CE	US\$ 195.44	+Info

Zirconia Discs MonoLayer - 3M ESPE

Zirconia Discs MonoLayer - 3M ESPE



Code	Name	Price	+
428367	Chairside Zirconia B1 - 20MM for CE	US\$ 188.03	+Info
428369	Chairside Zirconia A3 - 20MM for CE	US\$ 176.76	+Info
428372	Chairside Zirconia D2 - 20MM for CE	US\$ 188.03	+Info
428376	Chairside Zirconia A1 - 20MM for CE	US\$ 188.03	+Info
428377	Chairside Zirconia C1 - 20MM for CE	US\$ 188.03	+Info
428379	Chairside Zirconia A3.5 20mm 5pk	US\$ 188.03	+Info
428380	Chairside Zirconia BL - 20MM for CE	US\$ 188.03	+Info
428382	Chairside Zirconia A2 - 20MM for CE	US\$ 183.73	+Info

ZIRCONIA BLOCKS | DIGITECH

CEREC ZIRCONIA BLOCKS

CEREC ZIRCONIA BLOCKS - DIGITECH

* Our Zirconia blocks, Standard Color

Zirconia Blocks Chemical compositions:

- * $ZrO_2 + HfO_2 + Y_2O_3$ 99.6 %
- * Y_2O_3 5.2 %
- * Al_2O_3 0.5 %
- * Bending Strength 1200 Mpa
- * Vickers Hardness 12.8 Gpa
- * COLD ISOSTATIC PRESSING with hydrostatic pressure.

Firing Chart:

- * Sinterization temperature
- HS and ST 1530 Celsius 2 hours
- AT 1450 Celcius 2 hours

CE & FDA CERTIFIED

#UsddCerec



CEREC ANTERIOR TRANSLUCENCY Zirconia Blocks



Code	Name	Price	+
123217	DIGITECH - MONO CEREC Anterior Translucency Dental Zirconia Blocks 55x19x15.5mm...	US\$ 16.37	+Info

CEREC HS HIGH STRENGTH MONOLAYER Zirconia Blocks



Code	Name	Price	+
101977	DIGITECH - CEREC MONO HS Dental Zirconia Blocks HS MONO 20x14x15mm for MC XL &...	US\$ 3.34	+Info
101980	DIGITECH - CEREC MONO HS Dental Zirconia Blocks 40x19x15mm for MC XL & InLab -...	US\$ 6.96	+Info

KAVO ZIRCONIA BLOCKS

Zirconia Discs MonoLayer - DIGITECH



HS MONO WHITE

HS MONO WHITE

Code	Name	Price	+
473542	5X LOT of DIGITECH - HS MONO Dental	US\$ 194.40	+Info

MANUAL MILLING ZIRCONIA BLO CKS

ZIRKONZAHN ZIRCONIA BLOCKS - FOR MANUAL MILLING - DIGITECH

FDA CERTIFIED

* Our Zirconia blocks, Standard Color

Zirconia Blocks Chemical compositions:

* $ZrO_2 + HfO_2 + Y_2O_3$ 99.6 %

* Y_2O_3 5.2 %

* Al_2O_3 0.5 %

* Bending Strength 1200 Mpa

* Vickers Hardness 12.8 Gpa

* COLD ISOSTATIC PRESSING with hydrostatic pressure.

Firing Chart:

* Sinterization temperature

HS and ST 1530 Celsius 2 hours

AT 1450 Celsius 2 hours

ZIRKONZAHN is a trademark of ZIRKONZAHN GMBH



16mm

Code	Name	Price	+
111992	DIGITECH - ZIRKONZHAN Dental Zirconia discs #16 Low D-Form Semiround 93x75x16...	US\$ 34.22	+Info

PMMA TOOTH COLOR DISC

S



DIGILUX PMMA DISCS

MILLING Discs & Blocks

Features DIGILUX PMMA MILLING DISCS DIGILUX PMMA Discs are made of cross linked Polymethylmethacrylate with great natural light transmission. Indicated for individual restorations & bridges up to 4 elements

DIGILUX is a multilayered PMMA dental disc material from Brazil. This ultra aesthetic temporary material combines 16 layers of colors into one disc, giving the restoration the most natural look possible. Top quality gradient in DIGILUX PMMA discs is very translucent (incisal) and gradually shades down to the most intense bottom (cervical) with no flags effect.

Designed for long-term provisional crowns and bridges of any span

DIGILUX scope of materials

PMMA Mono & Multilayer
NANO Nano silica Multilayer
SUPER CRYSTAL Hybrid Ceramic Multilayer

DIGILUX PMMA is AVAILABLE

Monolayer and Multilayer with 16 Layers
Monolayer CLEAR, TOOTH COLOR or DENTURE PINK
Multilayer TOOTH COLOR 16 vita shades

Available in 98mm 93mm (Aman), 95mm (ZZ)
and Cerec compatible C14 and B40
Various thicknesses between 12 to 30mm

Multi layer , Mono layer, Clear, Pattern red & Pink for dentures4

Technical Parameters DIGILUX PMMA DISCS

Flexural Strength: 100Mpa / 120 mpa
DIGI-NANO and 165 mpa SUPER-CRYSTAL
Fast milling
Dry or Wet milling 6000 9000 rpms
recommended

DIGI-NANOMULTILAYER Discs & Blocks

Features DIGI-NANO milling discs

DIGILUX NANO ORMOCER PMMA Disks (



Organic Modified Ceramic) are made of a unique combination of Polymethylmethacrylate ensuring flexural strength, abrasion resistance and exceptional natural light transmission, resulting in high durability for long term provisionals in full mouth restorations

Chemical Compositions DIGI-NANO milling discs

Inorganic filler: **Nano Silica 7 microns 2-3%**

Radiopacity DIGI-NANO milling discs

With the addition of barium, becomes visible in X-rays in post treatment controls

SUPER CRYSTAL

HYBRID CERAMIC Discs & Blocks.

MULTILAYER Discs & Blocks FEATURES SUPER-CRYSTAL Hybrid milling discs

DIGILUX SUPER-CRYSTAL Hybrid Ceramic Discs are made of a perfect combination of Ceramic in the size of nano particles & the molecule Bis-Gma (bisphenol A-glycidyl methacrylate) resulting in a chemical bonding increasing the resistance to abrasion and flexural strength.

These physical properties SUPER-CRYSTAL Hybrid Ceramic Discs also reduce some of the down sides of zirconia restorations like clacking, weight, sinterization, specialized equipment & overall costs.

165 Mpa flexural strength

SUPER-CRYSTAL Hybrid Ceramic Discs are excellent choice for permanent restorations and long term provisionals in full mouth rehabilitations

Chemical Compositions SUPER-CRYSTAL Hybrid Ceramic Discs

Inorganic filler: **Nano Silica 30-55%**

Radiopacity SUPER-CRYSTAL Hybrid Ceramic Discs

With the addition of barium, becomes visible in X-rays in post treatment controls

MILLING TECHNIQUE

Preferable should be milled wet.

As most of the discs milling units in the market are DRY milling, consider always a slower RPM protocol. PMMA protocol will work

for Zirkonzahn machines choose protocol for material THERAPON

DIGI-NANO MULTI-LAYER PM

MA

DIGILUX PMMA Milling discsDIGI-NANO MULTILAYER Discs & Blocks

Features DIGI-NANO milling discs

DIGILUX NANO PMMA Disks (Organic Modified Ceramic) are made of a unique combination of Polymethylmethacrylate ensuring flexural strength, abrasion resistance and exceptional natural light transmission, resulting in high durability for long term provisionals in full mouth restorations

Chemical Compositions DIGI-NANO milling discs

Inorganic filler: Nano Silica 7 microns 2-3% & 120 mpa

Radiopacity DIGI-NANO milling discs with the addition of barium, becomes more visible in X-rays in post treatment controls



98mm

DIGILUX PMMA Milling discsDIGI-NANO MULTILAYER Discs & Blocks

Features DIGI-NANO milling discs

DIGILUX NANO PMMA Disks (Organic Modified Ceramic) are made of a unique combination of Polymethylmethacrylate ensuring flexural strength, abrasion resistance and exceptional natural light transmission, resulting in high durability for long term provisionals in full mouth restorations

Chemical Compositions DIGI-NANO milling discs

Inorganic filler: Nano Silica 7 microns 2-3% & 120 mpa

Radiopacity DIGI-NANO milling discs with the addition of barium, becomes more visible in X-rays in post treatment controls



98.5mm x 20mm

Code	Name	Price	+
464526	PMMA MULTI NANO 98.5x20mm A1	US\$ 50.40	+Info
464527	PMMA MULTI 98.5x20mm A2	US\$ 47.88	+Info
464528	PMMA MULTI 98.5x20mm A3	US\$ 50.40	+Info
464529	PMMA MULTI 98.5x20mm B1	US\$ 50.40	+Info
464530	PMMA MULTI 98.5x20mm B2	US\$ 48.38	+Info
474493	PMMA MULTI 98.5x20mm BL1	US\$ 50.40	+Info
474494	PMMA MULTI 98.5x20mm BL2	US\$ 50.40	+Info
464531	PMMA MULTI 98.5x20mm B0 BLEACH	US\$ 50.40	+Info

95mm

DIGILUX PMMA Milling discsDIGI-NANO MULTILAYER Discs & Blocks

Features DIGI-NANO milling discs

DIGILUX NANO PMMA Disks (Organic Modified Ceramic) are made of a unique combination of Polymethylmethacrylate ensuring flexural strength, abrasion resistance and exceptional natural light transmission, resulting in high durability for long term provisionals in full mouth restorations

Chemical Compositions DIGI-NANO milling discs

Inorganic filler: Nano Silica 7 microns 2-3% & 120 mpa

Radiopacity DIGI-NANO milling discs with the addition of barium, becomes more visible in X-rays in post treatment controls

16mm

Code	Name	Price	+
475045	NANO PMMA Disc Multi-Layer 95mm x 16mm Color : B1 ML	US\$ 50.40	+Info
475046	NANO PMMA Disc Multi-Layer 95mm x 16mm Color : B2 ML	US\$ 50.40	+Info
475047	NANO PMMA Disc Multi-Layer 95mm x 16mm Color : BL1 ML	US\$ 50.40	+Info
475048	NANO PMMA Disc Multi-Layer 95mm x 16mm Color : A1 ML	US\$ 50.40	+Info
475049	NANO PMMA Disc Multi-Layer 95mm x 16mm Color : A2 ML	US\$ 50.40	+Info
475050	NANO PMMA Disc Multi-Layer 95mm x 16mm Color : A3.5 ML	US\$ 50.40	+Info
475053	NANO PMMA Disc Multi-Layer 95mm x 16mm Color : A3 ML	US\$ 50.40	+Info

20mm

Code	Name	Price	+
464532	PMMA MULTI 95x20mm B1	US\$ 50.40	+Info
470506	PMMA MULTI 95x20mm A2	US\$ 50.40	+Info
464533	PMMA MULTI 95x20mm A1	US\$ 50.40	+Info
475044	NANO PMMA ML 95x20 BL1 BLEACH	US\$ 50.40	+Info
464534	NANO PMMA ML 95x20 B0 BLEACH	US\$ 50.40	+Info

MULTI-LAYER PMMA



98.5mm

16mm

Code	Name	Price	+
124503	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 16mm Color : A1 - #124503 ML	US\$ 36.00	+Info
124504	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 16mm Color : A2 - #124504 ML	US\$ 36.00	+Info
124505	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 16mm Color : A3 - #124505 ML	US\$ 36.00	+Info
124506	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 16mm Color : A3.5 - #124506 ML	US\$ 36.00	+Info

124507	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 16mm Color : A4 - #124507 ML	US\$ 35.28	+Info
124508	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 16mm Color : B1 - #124508 ML	US\$ 36.00	+Info
124509	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 16mm Color : B2 - #124509 ML	US\$ 36.00	+Info
473720	DIGILUX PMMA Disc MULTILAYER B0 Bleach 98.5mm (98mm) x 16mm 1 Round block per ...	US\$ 36.00	+Info

20mm

Code	Name	Price	+
124171	DIGILUX PMMA Disc MULTILAYER (ML) 98.5mm (98mm) x 20mm Color : A1 - #124171...	US\$ 36.00	+Info
124172	DIGILUX PMMA Disc MULTILAYER (ML) 98.5mm (98mm) x 20mm Color : A2 - #124172...	US\$ 36.00	+Info
124173	DIGILUX PMMA Disc MULTILAYER (ML) 98.5mm (98mm) x 20mm Color : A3 - #124173...	US\$ 36.00	+Info
124174	DIGILUX PMMA Disc MULTILAYER (ML) 98.5mm (98mm) x 20mm Color : A3.5 - #1241...	US\$ 35.28	+Info
124175	DIGILUX PMMA Disc MULTILAYER (ML) 98.5mm (98mm) x 20mm Color : A4 - #124175...	US\$ 36.00	+Info
124176	DIGILUX PMMA Disc MULTILAYER (ML) 98.5mm (98mm) x 20mm Color : B1 - #124176...	US\$ 36.00	+Info
124177	DIGILUX PMMA Disc MULTILAYER (ML) 98.5mm (98mm) x 20mm Color : B2 - #124177...	US\$ 36.00	+Info
476433	DIGILUX PMMA Disc MULTILAYER (ML) 98.5mm (98mm) x 20mm Color : C3 - #476433 ML	US\$ 36.00	+Info
476434	DIGILUX PMMA Disc MULTILAYER (ML) 98.5mm (98mm) x 20mm Color : D3 - #476433 ML	US\$ 36.00	+Info
416570	DIGILUX PMMA Disc MULTILAYER (ML) 98.5mm (98mm) x 20mm Color : B0 (510K K10...	US\$ 36.00	+Info
476429	DIGILUX PMMA Disc MULTILAYER (ML) 98.5mm (98mm) x 20mm Color : BL1 BLEACH (...)	US\$ 36.00	+Info
476430	DIGILUX PMMA Disc MULTILAYER (ML) 98.5mm (98mm) x 20mm Color : BL2 (510K K1...	US\$ 36.00	+Info

25mm

Code	Name	Price	+
124178	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 25mm Color : A1 - #124178 ML	US\$ 48.60	+Info
124179	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 25mm Color : A2 - #124179 ML	US\$ 48.60	+Info
124180	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 25mm Color : A3 - #ML 98.5 25MM A3	US\$ 48.60	+Info
124181	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 25mm Color : A3.5 - #124181 ML	US\$ 42.12	+Info
124182	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 25mm Color : A4 - #124182 ML	US\$ 35.57	+Info
417347	DIGILUX PMMA Disc MULTILAYER BL4 98.5mm (98mm) x 25mm 1pk ML (BL on the box)	US\$ 45.68	+Info
124183	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 25mm Color : B1 - #124183 ML	US\$ 48.60	+Info

124184	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 25mm Color : B2 - #124184 ML	US\$ 48.60	+Info
474262	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 25mm Color : C1- # ML	US\$ 46.80	+Info
474389	DIGILUX PMMA Disc MULTILAYER 95mm x 25mm Color : C2- # ML	US\$ 45.86	+Info

30mm

Code	Name	Price	+
476418	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 30mm Color : A2 - # ML	US\$ 54.00	+Info
476419	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 30mm Color : BL1- # ML	US\$ 54.00	+Info
476420	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 30mm Color : BL2 - # ML	US\$ 54.00	+Info
476421	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 30mm Color : A3 - # ML	US\$ 54.00	+Info
476417	DIGILUX PMMA Disc MULTILAYER 98.5mm (98mm) x 30mm Color : A1 - # ML	US\$ 54.00	+Info

95mm

95x16mm

Code	Name	Price	+
416965	DIGILUX PMMA Disc MULTILAYER Disc A1 ZZ 95mm x 16mm - 1 block per box - for Zirk...	US\$ 36.00	+Info
416966	DIGILUX PMMA Disc MULTILAYER Disc A2 ZZ 95mm x 16mm - 1 block per box - for Zirk...	US\$ 36.00	+Info
416967	DIGILUX PMMA Disc MULTILAYER Disc A3 ZZ 95mm x 16mm - 1 block per box - for Zirk...	US\$ 36.00	+Info
416968	DIGILUX PMMA Disc MULTILAYER Disc A3.5 ZZ 95mm x 16mm - 1 block per box - for Zi...	US\$ 36.00	+Info
123901	DIGILUX PMMA Disc MULTILAYER Disc B0 Bleach ZZ 95mm x 16mm - 1 block per box - f...	US\$ 36.00	+Info
475051	DIGILUX PMMA Disc MULTILAYER Disc B1 Bleach ZZ 95mm x 16mm - 1 block per box - f...	US\$ 36.00	+Info
475052	DIGILUX PMMA Disc MULTILAYER Disc BL1 Bleach ZZ 95mm x 16mm - 1 block per box - ...	US\$ 36.00	+Info

95x20mm

Code	Name	Price	+
124213	DIGILUX PMMA Disc MULTILAYER (ML) ZZ 95mm x 20mm Color : A1 - #124213 ML	US\$ 30.60	+Info
124214	DIGILUX PMMA Disc MULTILAYER (ML) ZZ 95mm x 20mm Color : A2 - #124214 ML	US\$ 34.56	+Info
124215	DIGILUX PMMA Disc MULTILAYER (ML) ZZ 95mm x 20mm Color : A3 - #124215 ML	US\$ 34.56	+Info
124216	DIGILUX PMMA Disc MULTILAYER (ML) ZZ 95mm x 20mm Color : A3.5 - #124216 ML	US\$ 34.56	+Info
124217	DIGILUX PMMA Disc MULTILAYER (ML) ZZ 95mm x 20mm Color : A4 - #124217 ML	US\$ 22.32	+Info
416571	DIGILUX PMMA Disc MULTILAYER (ML) B0 Bleach ZZ 95mm x 20mm 1 block per box - ...	US\$ 35.28	+Info
124218	DIGILUX PMMA Disc MULTILAYER (ML) ZZ 95mm x 20mm Color : B1 - #124218 ML	US\$ 36.00	+Info
124219	DIGILUX PMMA Disc MULTILAYER (ML) ZZ 95mm x 20mm Color : B2 - #124219 ML	US\$ 34.92	+Info

416756

DIGILUX PMMA Disc MULTILAYER (ML) ZZ 95mm x 20mm Color : C2 - #416756 ML

US\$ 35.28

[+Info](#)

95x25mm

Code	Name	Price	+
124220	DIGILUX PMMA Disc MULTILAYER ZZ 95mm x 25mm Color : A1 - #124220 ML	US\$ 45.68	+Info
124221	DIGILUX PMMA Disc MULTILAYER ZZ 95mm x 25mm Color : A2 - #124221 ML	US\$ 45.40	+Info
124222	DIGILUX PMMA Disc MULTILAYER ZZ 95mm x 25mm Color : A3 - #124222 ML	US\$ 45.40	+Info
124223	DIGILUX PMMA Disc MULTILAYER ZZ 95mm x 25mm Color : A3.5 - #124223 ML	US\$ 31.36	+Info
124224	DIGILUX PMMA Disc MULTILAYER ZZ 95mm x 25mm Color : A4 - #124224 ML	US\$ 28.55	+Info
417348	DIGILUX PMMA Disc MULTILAYER B0 Bleach ZZ 95mm x 25mm 1 block per box - for Zirk...	US\$ 45.68	+Info
124225	DIGILUX PMMA Disc MULTILAYER ZZ 95mm x 25mm Color : B1 - #124225 ML	US\$ 48.60	+Info
124226	DIGILUX PMMA Disc MULTILAYER ZZ 95mm x 25mm Color : B2 - #124226 ML	US\$ 45.40	+Info

71mm

Code	Name	Price	+
124241	DIGILUX PMMA Disc AMANN GIRRBACH MULTILAYER (ML) 93mm / 71mm x 20mm AMANN GI...	US\$ 22.93	+Info
124242	DIGILUX PMMA Disc AMANN GIRRBACH MULTILAYER (ML) 93mm / 71mm x 20mm AMANN GI...	US\$ 23.40	+Info
124243	DIGILUX PMMA Disc AMANN GIRRBACH MULTILAYER (ML) 93mm / 71mm x 20mm AMANN GI...	US\$ 20.64	+Info
124244	DIGILUX PMMA Disc AMANN GIRRBACH MULTILAYER (ML) 93mm / 71mm x 20mm AMANN GI...	US\$ 19.80	+Info
124246	DIGILUX PMMA Disc AMANN GIRRBACH MULTILAYER (ML) 93mm / 71mm x 20mm AMANN GI...	US\$ 23.40	+Info
124247	DIGILUX PMMA Disc AMANN GIRRBACH MULTILAYER (ML) 93mm / 71mm x 20mm AMANN GI...	US\$ 19.80	+Info
416755	DIGILUX PMMA Disc AMANN GIRRBACH MULTILAYER (ML) 93mm / 71mm x 20mm AMANN GI...	US\$ 19.80	+Info
416995	DIGILUX PMMA Disc AMANN GIRRBACH MULTILAYER (ML) 93mm / 71mm x 20mm AMANN GI...	US\$ 20.64	+Info
416996	DIGILUX PMMA Disc AMANN GIRRBACH MULTILAYER (ML) 93mm / 71mm x 20mm AMANN GI...	US\$ 22.93	+Info

MONOLAYER PMMA

MONOLAYER PMMA DISCS - DIGILUX

Features:

Our Monolayer PMMA discs are perfect for dental laboratories who want to reduce labor costs and improve consistency by using modern milling systems. They offer excellent physical and biocompatibility characteristics due to highly cross-linked polymers and result in superior surface finish.



98.5mm MONO PMMA DIGILUX

98.5x16mm

Code	Name	Price	+
124150	DIGILUX PMMA Disc MONO 98.5mm (98mm) x 16mm Color : A1 - #124150	US\$ 23.40	+Info
124151	DIGILUX PMMA Disc MONO 98.5mm (98mm) x 16mm Color : A2 - #124151	US\$ 22.93	+Info
124152	DIGILUX PMMA Disc MONO 98.5mm (98mm) x 16mm Color : A3 - #124152	US\$ 22.46	+Info
124153	DIGILUX PMMA Disc MONO 98.5mm (98mm) x 16mm Color : A3.5 - #124153	US\$ 23.40	+Info
124154	DIGILUX PMMA Disc MONO 98.5mm (98mm) x 16mm Color : A4 - #124154	US\$ 23.40	+Info
124155	DIGILUX PMMA Disc MONO 98.5mm (98mm) x 16mm Color : B1 - #124155	US\$ 23.40	+Info
124156	DIGILUX PMMA Disc MONO 98.5mm (98mm) x 16mm Color : B2 - #124156	US\$ 19.19	+Info

98.5x20mm

Code	Name	Price	+
124157	DIGILUX PMMA Disc MONO 98.5mm (98mm) x 20mm Color : A1 - #124157	US\$ 23.40	+Info
124158	DIGILUX PMMA Disc MONO 98.5mm (98mm) x 20mm Color : A2 - #124158	US\$ 23.40	+Info
124159	DIGILUX PMMA Disc MONO 98.5mm (98mm) x 20mm Color : A3 - #124159	US\$ 23.40	+Info
124160	DIGILUX PMMA Disc MONO 98.5mm (98mm) x 20mm Color : A3.5 - #124160	US\$ 23.40	+Info
124161	DIGILUX PMMA Disc MONO 98.5mm (98mm) x 20mm Color : A4 - #124161	US\$ 23.40	+Info
124162	DIGILUX PMMA Disc MONO 98.5mm (98mm) x 20mm Color : B1 Bleach- #124162	US\$ 30.60	+Info
124163	DIGILUX PMMA Disc MONO 98.5mm (98mm) x 20mm Color : B2 - #124163	US\$ 23.40	+Info

98.5x25mm

DIGILUX

Code	Name	Price	+
124164	PMMA MONO 98.5x25mm A1	US\$ 30.60	+Info
124166	PMMA MONO 98.5x25mm A3	US\$ 38.41	+Info
124167	PMMA MONO 98.5x25mm A3.5	US\$ 28.25	+Info
124168	PMMA MONO 98.5x25mm A4	US\$ 18.01	+Info
124169	PMMA MONO 98.5x25mm B1	US\$ 30.60	+Info
124170	PMMA MONO 98.5x25mm B2	US\$ 35.31	+Info
402379	PMMA MONO 98.5x25mm A0 BLEACH 1pk	US\$ 30.29	+Info

DIGITECH

Code	Name	Price	+
124165	PMMA MONO 98.5x25mm A2	US\$ 28.05	+Info

95mm MONO PMMA DIGILUX

95x16mm

Code	Name	Price	+
124192	DIGILUX PMMA Disc MONO ZZ 95mm x 16mm Color : A1 - for ZirkonZahn #124192	US\$ 19.66	+Info
124193	DIGILUX PMMA Disc MONO ZZ 95mm x 16mm Color : A2 - #124193	US\$ 23.40	+Info
124194	DIGILUX PMMA Disc MONO ZZ 95mm x 16mm Color : A3 - for ZirkonZahn #124194	US\$ 16.38	+Info
124195	DIGILUX PMMA Disc MONO ZZ 95mm x 16mm Color : A3.5 - for ZirkonZahn #124195	US\$ 17.32	+Info
124196	DIGILUX PMMA Disc MONO ZZ 95mm x 16mm Color : A4 - for ZirkonZahn #124196	US\$ 14.51	+Info
124197	DIGILUX PMMA Disc MONO ZZ 95mm x 16mm Color : B1 - for ZirkonZahn #124197	US\$ 21.06	+Info

124198

DIGILUX PMMA Disc MONO ZZ 95mm x 16mm Color : B2 - for ZirkonZahn #124198

US\$ 14.27

[+Info](#)

95x20mm

Code	Name	Price	+
124200	DIGILUX PMMA Disc MONO ZZ 95mm x 20mm Color : A2 - for ZirkonZahn #124200	US\$ 23.40	+Info
124201	DIGILUX PMMA Disc MONO ZZ 95mm x 20mm Color : A3 - for ZirkonZahn #124201	US\$ 21.06	+Info
124202	DIGILUX PMMA Disc MONO ZZ 95mm x 20mm Color : A3.5 - for ZirkonZahn #124202	US\$ 14.27	+Info
124203	DIGILUX PMMA Disc MONO ZZ 95mm x 20mm Color : A4 - for ZirkonZahn #124203	US\$ 14.27	+Info
124204	DIGILUX PMMA Disc MONO ZZ 95mm x 20mm Color : B1 - for ZirkonZahn #124204	US\$ 22.46	+Info
124205	DIGILUX PMMA Disc MONO ZZ 95mm x 20mm Color : B2 - for ZirkonZahn #124205	US\$ 23.40	+Info

95x25mm

Code	Name	Price	+
124206	DIGILUX PMMA Disc MONO ZZ 95mm x 25mm Color : A1 - for ZirkonZahn #124206	US\$ 23.40	+Info
124207	DIGILUX PMMA Disc MONO ZZ 95mm x 25mm Color : A2 - for ZirkonZahn #124207	US\$ 16.38	+Info
124208	DIGILUX PMMA Disc MONO ZZ 95mm x 25mm Color : A3 - for ZirkonZahn #124208	US\$ 16.38	+Info
124209	DIGILUX PMMA Disc MONO ZZ 95mm x 25mm Color : A3.5 - for ZirkonZahn #124209	US\$ 16.38	+Info
124210	DIGILUX PMMA Disc MONO ZZ 95mm x 25mm Color : A4 - for ZirkonZahn #124210	US\$ 16.38	+Info
124211	DIGILUX PMMA Disc MONO ZZ 95mm x 25mm Color : B1 - for ZirkonZahn #124211	US\$ 29.99	+Info
124212	DIGILUX PMMA Disc MONO ZZ 95mm x 25mm Color : B2 - for ZirkonZahn #124212	US\$ 23.40	+Info

71mm Amann Girbach

16mm

DIGITECH

DIGITECH

Code	Name	Price	+
113883	PMMA Disc Amann 93x71x16 A2	US\$ 13.11	+Info
113884	PMMA Disc Amann 93x71x16 A2	US\$ 16.72	+Info

ECCO

ECCO

Code	Name	Price	+
113889	PMMA Disc Amann 93x75x16 B1	US\$ 14.82	+Info

20mm

Code	Name	Price	+
124234	DIGILUX PMMA Disc AMANN GIRRBACH A MONO 93mm x 20mm AMANN GIRRBACH Color : A1 - ...	US\$ 13.11	+Info
124235	DIGILUX PMMA Disc AMANN GIRRBACH A MONO 93mm x 20mm AMANN GIRRBACH Color : A2 - ...	US\$ 17.67	+Info
124236	DIGILUX PMMA Disc AMANN GIRRBACH A MONO 93mm x 20mm AMANN GIRRBACH Color : A3 - ...	US\$ 13.11	+Info
124237	DIGILUX PMMA Disc AMANN GIRRBACH A MONO 93mm x 20mm AMANN GIRRBACH Color : A3.5 ...	US\$ 13.11	+Info
124238	DIGILUX PMMA Disc AMANN GIRRBACH A MONO 93mm x 20mm AMANN GIRRBACH Color : A4 - ...	US\$ 17.10	+Info

124239	DIGILUX PMMA Disc AMANN GIRRBACH A MONO 93mm x 20mm AMANN GIRRBACH Color : B1 - ...	US\$ 13.11	+Info
124240	DIGILUX PMMA Disc AMANN GIRRBACH A MONO 93mm x 20mm AMANN GIRRBACH Color : B2 - ...	US\$ 13.11	+Info

UNIDESA PMMA DISCS

unidesa · odi

MONOLAYER PMMA DISCS | UNIDESA

UNIDESA UNILAYER PMMA Discs MonoColor

A product you can trust for the preparation of crowns and bridges. Made with PMMA resins of a high molecular weight which have been tested many times and which fulfil all the current legal requirements and certifications on biocompatibility. The same used by Unidesa-Odi for the preparation of its artificial teeth. A material that provides hardness and resistance, perfect for the preparation of provisional pieces. You can make the most of the product since it can be reused many times. It is manufactured in a wide range of colours and it helps the professional to every day make profitable its investment in CAD-CAM technology.



98.5mm

16mm Wieland

Code	Name	Price	+
123240	IDODENTINE - PMMA 98.5mm/16mm/A1 Mono-layer Disc for Regular/Wieland/Opend Mfg #...	US\$ 45.66	+Info
123243	IDODENTINE - PMMA 98.5mm/16mm/B1 Mono-layer Disc for Regular/Wieland/Opend Mfg #...	US\$ 45.66	+Info
123244	IDODENTINE - PMMA 98.5mm/16mm/BG3 (Bleach WHITE (A0)) Mono-layer Disc for Regu...	US\$ 45.66	+Info

20mm Wieland

Code	Name	Price	+
123233	IDODENTINE - PMMA 98.5mm/20mm/A2 Mono-layer Disc for Regular/Wieland/Opend Mfg #...	US\$ 30.24	+Info

123234

IDODENTINE - PMMA 98.5mm/20mm/A3 Mono-layer Disc for Regular/Wieland/Opend
Mfg #...

US\$ 36.00
[+Info](#)

25 mm Wieland

Code

Name

Price

+

420279

IDODENTINE - PMMA 98.5mm/25mm/A1 Mono-layer Disc for Regular/Wieland/Opend
Mfg #...

US\$ 39.00
[+Info](#)

95mm

95mm x 16mm

Code

Name

Price

+

123258

IDODENTINE - PMMA ZZ 95mm/16mm/A3 Mono-layer Disc for Mfg #4252Z/A3 - for
Zirkon...

US\$ 41.12
[+Info](#)

123260

IDODENTINE - PMMA ZZ 95mm/16mm/BG3 Bleach WHITE Mono-layer Disc for Mfg
#4252Z/B...

US\$ 48.57
[+Info](#)

95mm x 20mm

Code

Name

Price

+

123250

IDODENTINE - PMMA ZZ 95mm/20mm/A3 Mono-layer Disc for Mfg #4250Z/A3 - for
Zirkon...

US\$ 33.48
[+Info](#)

101mm Amann

Code

Name

Price

+

123264

IDODENTINE - PMMA 93mm/71mm/20mm/A1 Mono-layer Disc for Ceramill Mfg
#4250AG/A1....

US\$ 42.88
[+Info](#)

123265

IDODENTINE - PMMA 93mm/71mm/20mm/A2 Mono-layer Disc for Ceramill . MFG
#4250AG/A...

US\$ 23.40
[+Info](#)

123266

IDODENTINE - PMMA 93mm/71mm/20mm/A3 Mono-layer Disc for Ceramill . Mfg
#4250AG/A...

US\$ 23.40
[+Info](#)

123267

IDODENTINE - PMMA 93mm/71mm/20mm/B1 Mono-layer Disc for Ceramill . Mfg
#4250AG/B...

US\$ 23.40
[+Info](#)

MULTILAYER PMMA | UNIDESA

UNIDESA - MULTILAYER PMMA Disc

Multicolor

Features:

- * For CAD-CAM milling machines
- * 100% biocompatible and antiallergenic
- * Permanent natural coloration (does not change once it has been processed)
- * Residual monomer Free
- * High resistance to abrasion
- * Easy to finish and polish
- * Made with PMMA resins of a high molecular weight.
- * Super High Cross Linked
- * High Speed milling

#USDDUNIDENT


98.5mm

16mm

Code	Name	Price	+
123278	IDODENTINE - PMMA 98.5mm/16mm/A1 Multi-Layer (5 Layers) Disc for Regular/Wiel...	US\$ 47.00	+Info
123279	IDODENTINE - PMMA 98.5mm/16mm/A2 Multi-Layer (5 Layers) Disc for Regular/Wiel...	US\$ 47.00	+Info
123280	IDODENTINE - PMMA 98.5mm/16mm/A3 Multi-Layer (5 Layers) Disc for Regular/Wiel...	US\$ 47.00	+Info
420865	IDODENTINE - PMMA 98.5mm/16mm/B1 Multi-Layer (5 Layers) Disc for Regular/Wiel...	US\$ 47.00	+Info
123282	IDODENTINE - PMMA 98.5mm/16mm/B3 Multi-Layer (5 Layers) Disc for Regular/Wiel...	US\$ 46.53	+Info
123283	IDODENTINE - PMMA 98.5mm/16mm/C2 Multi-Layer (5 Layers) Disc for Regular/Wiel...	US\$ 47.00	+Info

20mm

Code	Name	Price	+
123272	IDODENTINE - PMMA 98.5mm/20mm/A1 Multi-Layer (5 Layers) Disc for Regular/Wiel...	US\$ 49.99	+Info
123273	IDODENTINE - PMMA 98.5mm/20mm/A2 Multi-Layer (5 Layers) Disc for Regular/Wiel...	US\$ 49.99	+Info
123274	IDODENTINE - PMMA 98.5mm/20mm/A3 Multi-Layer (5 Layers) Disc for Regular/Wiel...	US\$ 48.49	+Info
123275	IDODENTINE - PMMA 98.5mm/20mm/B0 Multi-Layer (5 Layers) Disc for Regular/Wiel...	US\$ 49.99	+Info
420867	IDODENTINE - PMMA 98.5mm/20mm/B1 Multi-Layer (5 Layers) Disc for Regular/Wiel...	US\$ 49.99	+Info
123276	IDODENTINE - PMMA 98.5mm/20mm/B3 Multi-Layer (5 Layers) Disc for Regular/Wiel...	US\$ 49.99	+Info
123277	IDODENTINE - PMMA 98.5mm/20mm/C2 Multi-Layer (5 Layers) Disc for Regular/Wiel...	US\$ 49.99	+Info
464880	UNIDESIA - PMMA 98.5mm/20mm/ BL2 Multi-Layer (5 Layers) - Casting Disc 1pk - M...	US\$ 49.99	+Info
464881	UNIDESIA - PMMA 98.5mm/20mm/ SCAN Casting Disc 1pk - Mfg #4250/SCAN	US\$ 53.96	+Info

25mm & 30mm

PMMA CAD/CAM MILLING Discs - UNIDESIA

Code	Name	Price	+
470755	MULTI 98.5mm/25mm/A1	US\$ 60.85	+Info
470756	MULTI 98.5mm/25mm/A2	US\$ 60.85	+Info
470757	MULTI 98.5mm/25mm/BL1	US\$ 60.85	+Info
470758	MULTI 98.5mm/30mm/A1	US\$ 64.99	+Info
470759	MULTI 98.5mm/30mm/A2	US\$ 64.99	+Info
470760	MULTI 98.5mm/30mm/BL1	US\$ 64.99	+Info

95mm

16mm

Code	Name	Price	+
123290	IDODENTINE - PMMA ZZ 95mm/16mm/A1 Multi-Layer (5 Layers) Disc for - Mfg #4516...	US\$ 41.36	+Info
123291	IDODENTINE - PMMA ZZ 95mm/16mm/A2 Multi-Layer (5 Layers) Disc for - Mfg #4516...	US\$ 45.59	+Info
123292	IDODENTINE - PMMA ZZ 95mm/16mm/A3 Multi-Layer (5 Layers) Disc for - Mfg #4516...	US\$ 43.24	+Info
123293	IDODENTINE - PMMA ZZ 95mm/16mm/B0 Multi-Layer (5 Layers) Disc for - Mfg #4516...	US\$ 46.06	+Info
123294	IDODENTINE - PMMA ZZ 95mm/16mm/B3 Multi-Layer (5 Layers) Disc for - Mfg #4300...	US\$ 47.00	+Info
123295	IDODENTINE - PMMA ZZ 95mm/16mm/C2 Multi-Layer (5 Layers) Disc for - Mfg #4300...	US\$ 47.00	+Info

20mm

Code	Name	Price	+
123284	IDODENTINE - PMMA ZZ 95mm/20mm/A1 Multi-Layer (5 Layers) Disc for - Mfg #4500...	US\$ 49.99	+Info
123285	IDODENTINE - PMMA ZZ 95mm/20mm/A2 Multi-Layer (5 Layers) Disc for - Mfg #4500...	US\$ 49.99	+Info
123286	IDODENTINE - PMMA ZZ 95mm/20mm/A3 Multi-Layer (5 Layers) Disc for - Mfg #4500...	US\$ 49.99	+Info
123287	IDODENTINE - PMMA ZZ 95mm/20mm/B1 Multi-Layer (5 Layers) Disc for - Mfg #4500...	US\$ 49.99	+Info
123288	IDODENTINE - PMMA ZZ 95mm/20mm/B3 Multi-Layer (5 Layers) Disc for - Mfg #4500...	US\$ 38.99	+Info
123289	IDODENTINE - PMMA ZZ 95mm/20mm/C2 Multi-Layer (5 Layers) Disc for - Mfg #4500...	US\$ 44.49	+Info

71mm Amann Girschbach

Amann Girschbach

Code	Name	Price	+
123296	IDODENTINE - PMMA 93mm / 71mm /20mm/A1 Multi-Layer (5 Layers) Disc for Cerami...	US\$ 32.49	+Info
123297	IDODENTINE - PMMA 93mm / 71mm /20mm/A2 Multi-Layer (5 Layers) Disc for Cerami...	US\$ 32.49	+Info
123298	IDODENTINE - PMMA 93mm / 71mm /20mm/A3 Multi-Layer (5 Layers) Disc for Cerami...	US\$ 32.49	+Info
123299	IDODENTINE - PMMA 93mm / 71mm /20mm/B0 Multi-Layer (5 Layers) Disc for Cerami...	US\$ 32.49	+Info
123300	IDODENTINE - PMMA 93mm / 71mm /20mm/B3 Multi-Layer (5 Layers) Disc for Cerami...	US\$ 32.49	+Info
123301	IDODENTINE - PMMA 93mm / 71mm /20mm/C2 Multi-Layer (5 Layers) Disc for Cerami...	US\$ 32.49	+Info

VIPI PMMA DISCS



TRILUX LAYERED TOOTH COLOR PMMA DISC ViPi

TRILUX LAYERED TOOTH COLOR PMMA DISC ViPi

The TRILUX LAYERED DISC is used in the milling technique of temporary crowns and bridges by CAD/CAM systems of various manufacturers. The recommendations of the manufacturer of the Software and the machining center have to be followed.

Advantages:

- * Excellent Milling Capacity
- * Biocompatibility
- * High Color Stability

Benefits:

- * Bending Resistance
- * Natural Fluorescence
- * Biocompatibility
- * Color and Translucency
- * High Color Stability

Typical Mechanical Properties of PMMA Blocks:

- * Modulus of Elasticity: 2,200 MPa
- * Bending Strength: 100 MPa
- * Toughness (Charpy): 10 KJ/ m2 (Izod): 1.3 KJ/m2
- * Tensile Strength: 75 MPa
- * Vickers Hardness HV (0,05/10): 135
- * Absorption of Water according to ISO 10477: 20 ug/mm3
- * Water solubility: 0.8 ug/mm3
- * Residual Monomer Content (24 h at 37 C and 30 min in air): < 0.7 %

The TRILUX LAYERED DISC are composed of acrylic resins with OMC (Organically Modified Ceramics) Nanotechnology materials that elaborate the three layerd tooth TRILUX, Giving a visual appearance close to the natural look.



VIPI BLOCK - CAD/CAM MATERIALS
Various forms. Several compositions. Infinite possibilities.

ZIRKON ZAHN

Code	Name	Price	+
115280	TRILUX - PMMA Disc Multi 95x20 A1 1 pk ML	US\$ 40.43	+Info
115284	TRILUX Tooth Color Multi PMMA Disc A-2 ZZ 95mm x 20mm - 1 block per box - for Zi...	US\$ 33.29	+Info
115290	TRILUX Tooth Color Multi PMMA Disc A3 ZZ 95mm x 20mm - 1 block per box - for Zir...	US\$ 40.60	+Info

95mm PMMA DISCS

95mm PMMA DISCS TOOTH COLOR ZIRKONZAHN - VIPI

VIPI BLOCK is used in the milling technique of temporary crowns and bridges by CAD/CAM systems of various manufacturers. The recommendations of the manufacturer of the Software and the machining center have to be followed.

Typical Mechanical Properties of PMMA Blocks:

- * Modulus of Elasticity: 2,200 MPa



- * Bending Strength: 100 MPa
- * Toughness (Charpy): 10 KJ/ m2 (Izod): 1.3 KJ/m2
- * Tensile Strength: 75 MPa
- * Vickers Hardness HV (0,05/10): 135
- * Absorption of Water according to ISO 10477: 20 ug/mm3
- * Water solubility: 0.8 ug/mm3
- * Residual Monomer Content (24 h at 37 C and 30 min in air): < 0.7 %

#Usdd95

#Zirkon

ZIRKONZAHN is a trademark of ZIRKONZAHN GMBH

A1 Vita

Code	Name	Price	+
114059	VIPI-PMMA Tooth Color Mono Disc #A1 ZZ 95mm x 16mm - 1 block per box - for Zirko...	US\$ 14.40	+Info

98.5mm PMMA DISCS

VIPI BLOCK is used in the milling technique of temporary crowns and bridges by **CAD/CAM** sytems of various manufacturers. The recommendations of the manufacturer of the Software and the machining center hate to be followed.



Composition:

- * Polymethylmethacrylate
- * Polymerized Ethylene Glycol Dimethacrylate (EDMA)
- * OMC (Organically Modified Ceramics)
- * Fluoresent
- * Biocompatible Pigments
- * DCL - Double Cross Linkage

Typical Mechanical Properties of PMMA Blocks:

- * Modulus of Elasticity: 2,200 MPa
- * Bending Stregth: 100 MPa
- * Toughness (Charpy) : 10 KJ/ m2 (Izod) : 1.3 KJ/m2
- * Tensile Strength: 75 MPa
- * Vickers Hardness HV (0,05/10) : 135
- * Absorption of Water according to ISO 10477: 20 ug/mm3
- * Water solubility: 0.8 ug/mm3

Residual Monomer Content (24 h at 37 C and 30 min in air) : < 0.7 %

#USDD98.5

A1 Vita

Code	Name	Price	+
113880	VIPI-PMMA Tooth Color Mono A1 Disc 98.5mm x 20mm - WIELAND - 1 Round block per b...	US\$ 27.29	+Info

A2 Vita

Code	Name	Price	+
115697	VIPI-PMMA Tooth Color Mono A2 Disc 98.5mm x 20mm - WIELAND - 1 Round block per b...	US\$ 39.46	+Info

B1 Vita

Code	Name	Price	+
113891	VIPI-PMMA Tooth Color Mono B1 Disc 98.5mm x 16mm - WIELAND - 1 Round block per b...	US\$ 14.40	+Info

PMMA Cerec Blocks - ALL BRANDS



BLOCKS PMMA - SHOFU

BLOCKS PMMA - SHOFU

The composite resin that was densely packed the filler is pressurizing and heating polymerized hybrid resin block. By using a dental cad / cam system, before molars crown, you can use in the fabrication of the final prosthetic devices, such as inlay. The shofu block HC to processing in dwx-50 is required jig adapter (sold separately).

Available:

LT / 6 colors W2-LT, A1-LT, A2-LT, A3-LT, A3.5-LT, B3-LT

HT / 3-color A1-HT, A2-HT, A3-HT



Code	Name	Price	+
116916	SHOFU Pmma Cerec Block HC M Oc 5pk #2172 (blow out while inventory lasts)	US\$ 60.00	+Info

DIGILUX

DIGILUX

Code	Name	Price	+
124250	PMMA MULTI CERECx A3	US\$ 23.50	+Info

UNIDESA

UNIDESA

Cerec

Cerec

Code	Name	Price	+
123302	15mm/19mm/39mm/A1 3pk	US\$ 24.00	+Info
123303	15mm/19mm/39mm/A2 3pk	US\$ 24.00	+Info
123304	15mm/19mm/39mm/A3 3pk	US\$ 24.00	+Info

DIGILUX MONO

DIGILUX MONO

Code	Name	Price	+
400101	MONO 15mm/19mm/39mm/A1 3pk	US\$ 13.00	+Info
400102	MONO 15mm/19mm/39mm/A3. 3pk	US\$ 13.00	+Info

MONO CLEAR

MONO CLEAR

Code	Name	Price	+
123305	15mm/19mm/39mm/Violet 3pk	US\$ 10.40	+Info

HYBRID MULTI-LAYER CERA MIC DISCS



SUPER CRYSTAL MULTILAYER H YBRID CERAMIC DISCS - DIGILUX



98.5mm

16mm

16mm

DIGILUX SUPER CRYSTAL

DIGILUX SUPER CRYSTAL

Code	Name	Price	+
416796	DIGILUX Disc HC 98.5mm LT A2 16mm	US\$ 200.00	+Info
416798	DIGILUX Disc HC 98.5mm LT A3.5 16mm	US\$ 200.00	+Info
416799	DIGILUX Disc HC 98.5mm LT A4 16mm	US\$ 200.00	+Info
416902	DIGILUX Disc HC 98.5mm LT B1 16mm	US\$ 200.00	+Info
416906	DIGILUX Disc HC 98.5mm LT C2 16mm	US\$ 182.00	+Info
416907	DIGILUX Disc HC 98.5mm LT D3 16mm	US\$ 182.00	+Info
416908	DIGILUX Disc HC 98.5mm LT C3 16mm	US\$ 182.00	+Info
474899	DIGILUX Disc HC 98.5mm LT A1 16mm	US\$ 200.00	+Info

SUPER CRYSTAL

SUPER CRYSTAL

Code	Name	Price	+
416903	DIGILUX Disc HC 98.5x16mm BL3	US\$ 200.00	+Info

20mm

20mm

DIGILUX SUPER CRYSTAL

DIGILUX SUPER CRYSTAL

Code	Name	Price	+
421123	Super Crystal HYBRID 98.5x20 LT BL	US\$ 200.00	+Info
464876	DIGILUX Disc HC 98.5mm LT B1 20mm	US\$ 200.00	+Info
464877	DIGILUX Disc HC 98.5mm LT A1 20mm	US\$ 200.00	+Info
474333	DIGILUX Disc HC 98.5mm LT A2 20mm	US\$ 200.00	+Info
474898	DIGILUX Disc HC 98.5mm LT A3 20mm	US\$ 194.00	+Info

SUPER CRYSTAL

SUPER CRYSTAL

Code	Name	Price	+
464878	DIGILUX Disc HC 98.5x20mm LT BL3	US\$ 200.00	+Info
475043	DIGILUX Disc HC 98.5x20mm LT BL2	US\$ 200.00	+Info

25mm

25mm

Code	Name	Price	+
473695	DIGILUX Disc HC 98.5x25mm BL3	US\$ 265.42	+Info

95mm

Code	Name	Price	+
416971	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc 95mm LT A1 x 16mm	US\$ 192.00	+Info
474846	Zirc-Na-n DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc 95mm LT A2 x 16mm	US\$ 192.00	+Info
474847	Zirc-Na-n DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc 95mm LT A3 x 16mm	US\$ 192.00	+Info
416974	Zirc-Na-n DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc 95mm LT A3.5 x 16mm	US\$ 200.00	+Info
416975	Zirc-Na-n DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc 95mm LT A4 x 16mm	US\$ 200.00	+Info
	Zirc-Na-n...		

416976	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc 95mm LT B1 x 16mm	US\$ 192.00	+Info
464542	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc 95mm LT BL1 x 16mm	US\$ 198.00	+Info
416981	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc 95mm LT C3 x 16mm	US\$ 200.00	+Info
416978	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc 95mm LT D2 x 16mm	US\$ 200.00	+Info
416980	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc 95mm LT D3 x 16mm	US\$ 200.00	+Info

93mm Super Crystal

Code	Name	Price	+
416909	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc AG 93mm x 71mm LT A1 x 16mm...	US\$ 144.00	+Info
416940	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc AG 93mm x 71mm LT A2 x 16mm...	US\$ 144.00	+Info
416943	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc AG 93mm x 71mm LT A3 x 16mm...	US\$ 144.00	+Info
416944	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc AG 93mm x 71mm LT A3.5 x 16...	US\$ 144.00	+Info
416945	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc AG 93mm x 71mm LT A4 x 16mm...	US\$ 144.00	+Info
416941	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc AG 93mm x 71mm LT B1 x 16mm...	US\$ 144.00	+Info
416942	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc AG 93mm x 71mm LT B2 x 16mm...	US\$ 144.00	+Info
416948	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc AG 93mm x 71mm LT C2 x 16mm...	US\$ 144.00	+Info
416949	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc AG 93mm x 71mm LT C3 x 16mm...	US\$ 144.00	+Info
416946	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc AG 93mm x 71mm LT D2 x 16mm...	US\$ 144.00	+Info
416947	DIGILUX SUPER CRYSTAL Hybrid Ceramic Multilayer Disc AG 93mm x 71mm LT D3 x 16mm...	US\$ 144.00	+Info

HC HYBRID CERAMIC MONO DISCS 98.5MM - SHOFU

HC DISC 98.5MM - SHOFU



Features:

Shofu Blocks & Disks HC are made of a unique hybrid ceramic material combination ensuring exceptional natural light transmission as well as high durability and flexural strength.

This new, highly aesthetic CAD/CAM material is suitable for a wide range of indications - including minimally invasive inlays, onlays, cosmetic veneers, full crowns for anterior and posterior teeth, and implant-supported restorations.

- * Force-absorbing hybrid ceramic
- * High wear resistance
- * Tooth-like light transmission and fluorescence
- * Excellent stain resistance
- * Precise milling-fast chairtime
- * Easy high gloss polishing

Size: 14mm

Available in Color: A1-A2- A3-A3.5-B3- Incisal

HT

Code	Name	Price	+
416578	SHOFU Disc HC 98.5mm HT A1 x 14mm Zirconium - silicate nanoceramic CAD/CAM HYBR...	US\$ 245.00	+Info
416579	SHOFU Disc HC 98.5mm HT A2 x 14mm Zirconium - silicate nanoceramic CAD/CAM HYBR...	US\$ 242.55	+Info
416580	SHOFU Disc HC 98.5mm HT A3 x 14mm Zirconium - silicate nanoceramic CAD/CAM HYBR...	US\$ 245.00	+Info

Incisal

Code	Name	Price	+
416581	SHOFU Disc HC 98.5mm 59 INCISAL x 14mm Zirconium - silicate nanoceramic CAD/CAM...	US\$ 242.55	+Info

LT

Code	Name	Price	+
416573	SHOFU Disc HC 98.5mm LT A1 x 14mm Zirconium - silicate nanoceramic CAD/CAM	US\$ 242.55	+Info
416574	SHOFU Disc HC 98.5mm LT A2 x 14mm Zirconium - silicate nanoceramic CAD/CAM	US\$ 245.00	+Info
416575	SHOFU Disc HC 98.5mm LT A3 x 14mm Zirconium - silicate nanoceramic CAD/CAM	US\$ 245.00	+Info
416576	SHOFU Disc HC 98.5mm LT A3.5 x 14mm Zirconium - silicate nanoceramic CAD/CAM HY...	US\$ 245.00	+Info
416577	SHOFU Disc HC 98.5mm LT B3 x 14mm Zirconium - silicate nanoceramic CAD/CAM	US\$ 245.00	+Info

LITHIUM DISILICATE



Lithium Disilicate - DIGITECH

Press Lithium Disilicate | DIGITECH
CH
Pressable Ingots - DIGITECH



HT High Translucen
HT High Translucen

Code	Name	Price	+
421165	DISIL. ingot HT BL1 5pk	US\$ 45.00	+Info
421166	DISIL. ingot HT BL2 5pk	US\$ 45.00	+Info
BB5836	DISIL. ingot HT A2 5pk	US\$ 45.00	+Info
BB6129	DISIL. ingot HT B4 5pk	US\$ 45.00	+Info
BB6470	DISIL. ingot HT A3 5pk	US\$ 45.00	+Info
BB6494	DISIL. ingot HT A1 5pk	US\$ 45.00	+Info
BB6495	DISIL. ingot HT B1 5pk	US\$ 45.00	+Info

LT Low Translucenc
LT Low Translucenc

Code	Name	Price	+
476422	DISIL. ingot LT BL1 5pk	US\$ 45.00	+Info
476423	DISIL. ingot LT BL2 5pk	US\$ 45.00	+Info
BB6427	DISIL. ingot LT A1 5pk	US\$ 45.00	+Info
BB6430	DISIL. ingot LT B1 5pk	US\$ 45.00	+Info
BB6491	DISIL. ingot LT A2 5pk	US\$ 45.00	+Info
BB6492	DISIL. ingot LT A3 5pk	US\$ 45.00	+Info
BB6493	DISIL. ingot LT C2 5pk	US\$ 45.00	+Info

MO Medium Opacity
PRESS TABS

Code	Name	Price	+
BB5122	DISIL. ingot MO1 5pk	US\$ 45.00	+Info

Milling Lithium Disilicate - DIGITECH

Lithium Disilicate - DIGITECH



HT (High Translucency)

HT (High Transluc

C14

Code	Name	Price	+
BB3178	DIGITECH - HT A1 C14 Lithium Disilicate Milling Block - #123178	US\$ 9.50	+Info
BB6126	DIGITECH - HT A2 C14 Lithium Disilicate Milling Block - #116126	US\$ 9.50	+Info
BB3179	DIGITECH - HT A3 C14 Lithium Disilicate Milling Block - #123179	US\$ 9.50	+Info
469860	DIGITECH - HT A3.5 C14 CAD -Lithium Disilicate Milling Block - #469860	US\$ 8.74	+Info
473705	DIGITECH - HT A4 C14 Lithium Disilicate Milling Block - HT A4 C14	US\$ 9.50	+Info
BB6463	DIGITECH - HT B1 C14 Lithium Disilicate Milling Block HT B1 C14 High Translucenc...	US\$ 9.50	+Info
469861	DIGITECH - HT B2 C14 Lithium Disilicate Milling Block - #469861	US\$ 9.50	+Info
473706	DIGITECH - HT B4 C14 Lithium Disilicate Milling Block - HT B4 C14	US\$ 9.50	+Info
469792	DIGITECH - HT BL1 C14 Lithium Disilicate Milling Block - HT BL1 C14	US\$ 9.31	+Info
473747	DIGITECH - HT BL3 C14 Lithium Disilicate Milling Block - HT BL3 C14	US\$ 9.50	+Info
473707	DIGITECH - HT C3 C14 Lithium Disilicate Milling Block - HT C3 C14	US\$ 9.50	+Info

LT (Low Translucency)

LT (Low Translucency)

B40

B40

Code	Name	Price	+
BB2982	GLASS CER BLOCK LT B1 40mm	US\$ 33.78	+Info
BB2983	GLASS CER BLOCK LT A3 40mm	US\$ 33.78	+Info
BB2984	GLASS CER BLOCK LT A2 40mm B40	US\$ 33.78	+Info

C14

C14

Code	Name	Price	+
BB3180	GLASS CER BLOCK LT BL1 C14	US\$ 8.31	+Info
BB6122	GLASS CER BLOCK LT B1 C14 e	US\$ 9.50	+Info

BB6123	GLASS CER BLOCK LT A3 C14	US\$ 9.50	+Info
BB5349	GLASS CER BLOCK LT A1 C14 ea	US\$ 9.50	+Info
BB5350	GLASS CER BLOCK LT A2 C14	US\$ 9.50	+Info
BB7456	GLASS CER BLOCK LT C1 C14	US\$ 9.50	+Info
BB7457	GLASS CER BLOCK LT C2 C14 e	US\$ 9.50	+Info
BB7458	GLASS CER BLOCK LT B2 C14 e	US\$ 9.50	+Info
BB7459	GLASS CER BLOCK LT B3 C14	US\$ 9.50	+Info
469866	GLASS CERAMIC BLOCK LT A3.5 C14	US\$ 9.50	+Info
473882	GLASS CERAM BLOCK LT BL2 C14	US\$ 8.31	+Info
473883	GLASS CER BLOCK LT BL3 C14	US\$ 9.50	+Info

MO (Medium Opacity)

MO (Medium Opacity)

Code	Name	Price	+
402366	GLASS CER BLOCK MO1 C14 e	US\$ 9.50	+Info

Lithium Disilicate - AMBER HASS

MILLING BLOCKS - HASS



AMBER MILL LITHIUM DISILICATE - HASS

AMBER MILL LITHIUM DISILICATE - HASS

Features:

Amber Mill is the first dental glass ceramic milling block to be released in lithium disilicate. Excellent processability, enhanced physical properties and high aesthetics provide benefits for both patients and clinics.

Natural milky, fluorescent and multi-chromatic gradient effects.

Prostheses finished with Amber Mill can see the difference in quality as a result. Amber Mill's excellent whiteness and fluorescence create a natural color spectrum from the cervical area to the cut / occlusal surface of the prosthesis, complete with glazing without staining.

Denser and crosslinked crystal structure. Excellent product properties.

Amber Mill's compact, cross-linked crystal structure provides excellent product properties. After final crystallization heat treatment, Amber Mill's biaxial flexural strength is approximately 10% higher than other products.



I decide for myself what lightness I want .

With just one Amber Mill, you can create prostheses with varying light transmission. After selecting the shades you need, heat treatment the light at the temperature corresponding to the desired light transmittance.

HT

Code	Name	Price	+
420381	HASS - Amber Mill A1 x 5 ingots 12x14x18mm #AMB C14 A1	US\$ 101.00	+Info

LT

Code	Name	Price	+
420372	HASS - Amber Mill A1 x 5 ingots 12x14x18mm #AMB C14 A1	US\$ 101.00	+Info
420373	HASS - Amber Mill A2 x 5 ingots 12x14x18mm #AMB C14 A2	US\$ 101.00	+Info
420374	HASS - Amber Mill A3 x 5 ingots 12x14x18mm #AMB C14 A3	US\$ 101.00	+Info
420375	HASS - Amber Mill A3.5 x 5 ingots 12x14x18mm #AMB C14 A3.5	US\$ 101.00	+Info
420376	HASS - Amber Mill B1 x 5 ingots 12x14x18mm #AMB C14 B1	US\$ 101.00	+Info
420377	HASS - Amber Mill B2 x 5 ingots 12x14x18mm #AMB C14 B2	US\$ 101.00	+Info
420378	HASS - Amber Mill B3 x 5 ingots 12x14x18mm #AMB C14 B3	US\$ 101.00	+Info
420379	HASS - Amber Mill C2 x 5 ingots 12x14x18mm #AMB C14 C2	US\$ 101.00	+Info
420380	HASS - Amber Mill D3 x 5 ingots 12x14x18mm #AMB C14 D3	US\$ 101.00	+Info
476443	Amber press master MT (Medium Traslucency) A2.5	US\$ 83.05	+Info
476442	Amber press master MT (Medium Traslucency) A1.5	US\$ 83.05	+Info
476441	Amber press master MT (Medium Traslucency) A0.5	US\$ 83.05	+Info
476440	Amber press master MT (Medium Traslucency) W0.5	US\$ 83.05	+Info
476439	Amber press master HT+ (High Traslucency Plus) A2.5	US\$ 83.05	+Info
476438	Amber press master HT+ (High Traslucency Plus) A1.5	US\$ 83.05	+Info
476437	Amber press master HT+ (High Traslucency Plus) A0.5	US\$ 83.05	+Info
476436	Amber press master HT+ (High Traslucency Plus) W0.5	US\$ 90.00	+Info
476444	Amber press master LO (Low Opacity) 0	US\$ 83.05	+Info
476445	Amber press master LO (Low Opacity) 2	US\$ 83.05	+Info
420370	HASS - Amber Mill W1 x 5 ingots 12x14x18mm #AMB C14 W1	US\$ 101.00	+Info
420371	HASS - Amber Mill W2 x 5 ingots 12x14x18mm #AMB C14 W2	US\$ 101.00	+Info
420979	HASS - Amber Mill W3 x 5 ingots 12x14x18mm #AMB C14 W3	US\$ 101.00	+Info
421124	HASS - Amber Mill W4 x 5 ingots 12x14x18mm #AMB C14 W4	US\$ 101.00	+Info

ROSETTA BM LEUCITE - HASS

ROSETTA BM LEUCITE - HASS

Leucite Glass-Ceramic Blocks for Dental CAD/CAM System

Rosetta BM is the leucite glass ceramic CAD/CAM block (SAME AS EMP-RESS CAD IVO CLAR) with minimized shrinkage for simple and efficient restorations, stable and durable aesthetics

with high precision in the outcomes.

Simple & Fast Process

Efficiency and precision. Succeed in fulfilling two targets at once.

Rosetta BM enables post-milling process to be free from additional heat application. Enjoy the highest level of efficiency and satisfaction by achieving simplicity, fastness and efficiency with the

utmost accur

Durability in Aesthetics

True beauty.

The value of durability.

Lifelike and natural look of the resulted restorations without any failures such as discoloration, debonded case, etc. is one of the most important criteria and Rosetta BM is completely fulfilling

those standards. Excellent opalescence, fluorescence and brightness of Rosetta BM provides you with highly integrated aesthetic value.

High Edge Stability

Purpose-oriented stability and versatility.

uperior edge stability and long-lasting state of restoration after cementation are quintessential to evaluate the quality of dental work. In particular, those points become even more prominent

when considering smaller sized works such as inlays and onlays. Rosetta BM is your unquestionable solution for such cases.



Quantity: 5 ingots per pack

HT (High translucency)

Code	Name	Price	+
124089	HASS - ROSETTA BM HT A1 x 5 ingots 12x14x18mm #C14 HT A1	US\$ 82.00	+Info
124090	HASS - ROSETTA BM HT A2 x 5 ingots 12x14x18mm #C14 L HT A2	US\$ 82.00	+Info
124091	HASS - ROSETTA BM HT A3 x 5 ingots 12x14x18mm #C14 HT A3	US\$ 82.00	+Info
124092	HASS - ROSETTA BM HT A3.5 x 5 ingots 12x14x18mm #C14 HT A3.5	US\$ 82.00	+Info
124094	HASS - ROSETTA BM HT B2 x 5 ingots 12x14x18mm #C14 HT B2	US\$ 82.00	+Info
124095	HASS - ROSETTA BM HT B3 x 5 ingots 12x14x18mm #C14 HT B3	US\$ 82.00	+Info
124096	HASS - ROSETTA BM HT C2 x 5 ingots 12x14x18mm #C14 HT C2	US\$ 82.00	+Info
124097	HASS - ROSETTA BM HT D3 x 5 ingots 12x14x18mm #C14 HT D3	US\$ 82.00	+Info

LT (Low translucency)

Code	Name	Price	+
124080	HASS - ROSETTA BM LT A1 x 5 ingots 12x14x18mm #C14 LT A1	US\$ 82.00	+Info
124081	HASS - ROSETTA BM LT A2 x 5 ingots 12x14x18mm #C14 LT A2	US\$ 82.00	+Info
124082	HASS - ROSETTA BM LT A3 x 5 ingots 12x14x18mm #C14 LT A3	US\$ 82.00	+Info
124083	HASS - ROSETTA BM LT A3.5 x 5 ingots 12x14x18mm #C14 LT A3.5	US\$ 82.00	+Info
124084	HASS - ROSETTA BM LT B1 x 5 ingots 12x14x18mm #C14 LT B1	US\$ 65.60	+Info
124085	HASS - ROSETTA BM LT B2 x 5 ingots 12x14x18mm #C14 LT B2	US\$ 82.00	+Info
124086	HASS - ROSETTA BM LT B3 x 5 ingots 12x14x18mm #C14 LT B3	US\$ 82.00	+Info
124087	HASS - ROSETTA BM LT C2 x 5 ingots 12x14x18mm #C14 LT C2	US\$ 82.00	+Info
124088	HASS - ROSETTA BM LT D3 x 5 ingots 12x14x18mm #C14 LT D3	US\$ 82.00	+Info
124066	HASS - ROSETTA BM LT W1 x 5 ingots 12x14x18mm #C14 LT W1	US\$ 82.00	+Info
124079	HASS - ROSETTA BM LT W2 x 5 ingots 12x14x18mm #C14 LT W2	US\$ 82.00	+Info

Press Lithium Disilicate | AMBER HA SS



HT (High translucency)
HT

Code	Name	Price	+
474663	Amber Press:HT W3 - Press Ingots 5p	US\$ 75.00	+Info
474696	Amber Press:HT D4 (R10) - PPressing	US\$ 75.00	+Info
474695	Amber Press:HT D3 (R10) - Pressing	US\$ 75.00	+Info
474694	Amber Press:HT D2 (R10) - Pressing	US\$ 75.00	+Info
474693	Amber Press:HT C4 (R10) - PPressing	US\$ 75.00	+Info
474692	Amber Press:HT C3 (R10) - Pressing	US\$ 75.00	+Info
474691	Amber Press:HT C2 (R10) - Pressing	US\$ 75.00	+Info
474690	Amber Press:HT C1 (R10) - Pressing	US\$ 75.00	+Info
474689	Amber Press:HT B4 (R10) - Pressing	US\$ 75.00	+Info
474664	Amber Press:HT W4 - Press Ingots 5p	US\$ 75.00	+Info
474654	Amber Press:HT A1 - Press Ingots 5p	US\$ 75.00	+Info
474662	Amber Press:HT W2 - Press Ingots 5p	US\$ 75.00	+Info
474661	Amber Press:HT W1 - Press Ingots 5p	US\$ 75.00	+Info
474660	Amber Press:HT B3 (R10) - Pressing	US\$ 75.00	+Info
474659	Amber Press:HT B2 - Press Ingots 5p	US\$ 75.00	+Info
474658	Amber Press:HT B1 - Press Ingots 5p	US\$ 75.00	+Info
474657	Amber Press:HT A3.5 - Press Ingots	US\$ 75.00	+Info
474656	Amber Press:HT A3 - Press Ingots 5p	US\$ 75.00	+Info
474655	Amber Press:HT A2 - Press Ingots 5p	US\$ 75.00	+Info
474893	Amber Press:W0 HT - Press 5pk	US\$ 75.00	+Info

LT (Low translucency)

LT (Low translucency)

Code	Name	Price	+
474674	Amber Press:LT C2 - Press Ingots 5p	US\$ 75.00	+Info
474683	Amber Press:LT W4 - Press Ingots 5p	US\$ 75.00	+Info
474682	Amber Press:LT W3 - Press Ingots 5p	US\$ 75.00	+Info
474681	Amber Press:LT W2 - Press Ingots 5p	US\$ 75.00	+Info
474680	Amber Press:LT W1 - Press Ingots 5p	US\$ 75.00	+Info
474679	Amber Press:LT D4 - Press Ingots 5p	US\$ 75.00	+Info
474678	Amber Press:LT D3 - Press Ingots 5p	US\$ 75.00	+Info
474677	Amber Press:LT D2 - Press Ingots 5p	US\$ 75.00	+Info
474676	Amber Press:LT C4 - Press Ingots 5p	US\$ 75.00	+Info
474675	Amber Press:LT C3 - Press Ingots 5p	US\$ 75.00	+Info
474665	Amber Press:LT A1 - Press Ingots 5p	US\$ 75.00	+Info
474673	Amber Press:LT C1 - Press Ingots 5p	US\$ 75.00	+Info
474672	Amber Press:LT B4 - Press Ingots 5p	US\$ 75.00	+Info
474671	Amber Press:LT B3 - Press Ingots 5p	US\$ 75.00	+Info
474670	Amber Press:LT B2 - Press Ingots 5p	US\$ 75.00	+Info
474669	Amber Press:LT B1 - Press Ingots 5p	US\$ 75.00	+Info
474668	Amber Press:LT A3.5 - Press Ingots	US\$ 75.00	+Info
474667	Amber Press:LT A3 - Press Ingots 5p	US\$ 75.00	+Info
474666	Amber Press:LT A2 - Press Ingots 5p	US\$ 75.00	+Info
474894	Amber Press:W0 LT - Press 5pk	US\$ 75.00	+Info

Medium Opacity MO

Other Features

Code	Name	Price	+
474684	Amber Press:MO0 - Press Ingots 5pk	US\$ 75.00	+Info
474685	Amber Press:MO1 - Press Ingots 5pk	US\$ 75.00	+Info
474686	Amber Press:MO2 - Press Ingots 5pk	US\$ 75.00	+Info
474687	Amber Press:MO3 - Press Ingots 5pk	US\$ 75.00	+Info
474688	Amber Press:MO4 - Press Ingots 5pk	US\$ 75.00	+Info

Shade Guides - AMBER HAAS

Shade Guides - AMBER HAAS

Amber MILL

Amber MILL

Code	Name	Price	+
474506	Selection Tip Amber MILL	US\$ 98.01	+Info

AMBER PRESS

AMBER PRESS

Code	Name	Price	+
474505	Selection Tip Rossetta Amber Press	US\$ 98.01	+Info

INITIAL LISI MILL BLOCKS - GC AMERICA



INITIAL LISI BLOCK

INITIAL LISI BLOCK

Code	Name	Price	+
474394	LiSi Block 14 Cerec A2	US\$ 171.06	+Info

Other Models

Other Models

Code	Name	Price	+
474401	Initial LiSi Block 14, Cerec A1 HT	US\$ 171.06	+Info
474402	Initial LiSi Block 14, Cerec A3 HT	US\$ 171.06	+Info
474403	Initial LiSi Block 14, Cerec B1 HT	US\$ 169.35	+Info
474404	Initial LiSi Block 14, Cerec A1 LT	US\$ 171.06	+Info
474405	Initial LiSi Block 14, Cerec B1 LT	US\$ 171.06	+Info

Lithium Disilicate - IVOCLAR

IPS E.MAX Lithium Disilicate Blocks - IVOCLAR VIVADENT



HT

B40

Code	Name	Price	+
442438	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab HT Block, Shade A3 Size B40 ~...	US\$ 428.39	+Info
444272	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab HT Block, Shade BL2 Size B40,...	US\$ 0.00	+Info

C14

Code	Name	Price	+
442207	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar HT Block, Shade A1 Size C14, Pack...	US\$ 222.37	+Info
442196	IVOCLAR VIVADENT - IPS e.max CAD HT Block, Shade A2 Size C14, 5/Pk. is a system ...	US\$ 233.99	+Info
442203	IVOCLAR VIVADENT - IPS e.max CAD HT Block, Shade A3 Size C14, 5/Pk. is a system ...	US\$ 238.28	+Info
442230	IVOCLAR VIVADENT - IPS e.max CAD HT Block, Shade A3.5 Size C14, 5/Pk. is a syste...	US\$ 236.08	+Info
442262	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar HT Block, Shade A4 Size C14, Pack...	US\$ 228.22	+Info
442240	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar HT Block, Shade B1 Size C14, Pack...	US\$ 210.20	+Info
442310	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar HT Block, Shade B2 Size C14, Pack...	US\$ 207.99	+Info
442420	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar HT Block, Shade B3 Size C14, Pack...	US\$ 235.51	+Info
442445	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab HT Block, Shade BL1 Size C14,...	US\$ 206.67	+Info
442352	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab HT Block, Shade BL2 Size C14,...	US\$ 289.41	+Info
442353	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab HT Block, Shade BL3 Size C14,...	US\$ 215.12	+Info
442270	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar HT Block, Shade C1 Size C14, Pack...	US\$ 240.96	+Info
442259	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar HT Block, Shade C2 Size C14, Pack...	US\$ 225.06	+Info
442274	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar HT Block, Shade C3 Size C14, Pack...	US\$ 224.18	+Info

442390	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar HT Block, Shade C4 Size C14, Pack...	US\$ 242.64	+Info
442269	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar HT Block, Shade D2 Size C14, Pack...	US\$ 234.97	+Info
442282	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar HT Block, Shade D3 Size C14, Pack...	US\$ 234.81	+Info

INLAB HT

Code	Name	Price	+
442256	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab HT Block, Shade A1 Size I12, ...	US\$ 232.78	+Info
442219	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab HT Block, Shade A2 Size I12, ...	US\$ 234.44	+Info
442229	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab HT Block, Shade A3 Size I12, ...	US\$ 234.47	+Info
442322	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab HT Block, Shade A3.5 Size I12...	US\$ 234.25	+Info
442365	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab HT Block, Shade B1 Size I12, ...	US\$ 209.53	+Info
444291	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab HT Block, Shade B2 Size I12, ...	US\$ 245.04	+Info

LT

A16(L)

Code	Name	Price	+
442313	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade A1 Size A16(L) 5/...	US\$ 520.74	+Info
442242	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade A2 Size A16(L) 5/...	US\$ 504.14	+Info
442247	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade A3 Size A16(L) 5/...	US\$ 449.15	+Info
442294	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade A3.5 Size A16(L) ...	US\$ 476.46	+Info

A16(S)

Code	Name	Price	+
442458	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade A1 Size A16(S) 5/...	US\$ 450.10	+Info
442293	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade A2 Size A16(S) 5/...	US\$ 337.99	+Info
442517	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade A3 Size A16(S) 5/...	US\$ 487.01	+Info
442705	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade A3.5 Size A16(S) ...	US\$ 488.05	+Info
442684	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade B1 Size A16(S) 5/...	US\$ 508.74	+Info
444264	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade B2 Size A16(S) 5/...	US\$ 570.83	+Info
444273	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade BL2 Size A16(S) 5...	US\$ 388.26	+Info

444261	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade C1Size A16(S) 5/P...	US\$ 481.70	+Info
444287	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade D2 Size A16(S) 5/...	US\$ 570.75	+Info

B32

Code	Name	Price	+
442241	IVOCLAR VIVADENT - IPS e.max CAD Blocks LT, Size B32 Shade A1 3/Pk. For CAD/CAM ...	US\$ 417.81	+Info
442218	IVOCLAR VIVADENT - IPS e.max CAD Blocks LT, Size B32 Shade A2 3/Pk. For CAD/CAM ...	US\$ 339.74	+Info
442232	IVOCLAR VIVADENT - IPS e.max CAD Blocks LT, Size B32 Shade A3 3/Pk. For CAD/CAM ...	US\$ 429.49	+Info
442301	IVOCLAR VIVADENT - IPS e.max CAD Blocks LT, Size B32 Shade A3.5 3/Pk. For CAD/CA...	US\$ 417.79	+Info
442324	IVOCLAR VIVADENT - IPS e.max CAD Blocks LT, Size B32 Shade B1 3/Pk. For CAD/CAM ...	US\$ 435.27	+Info
442315	IVOCLAR VIVADENT - IPS e.max CAD Blocks LT, Size B32 Shade D2 3/Pk. For CAD/CAM ...	US\$ 359.26	+Info

BL32

Code	Name	Price	+
442527	IVOCLAR VIVADENT - IPS e.max CAD Blocks LT, Size BL2 Shade B1 3/Pk. For CAD/CAM ...	US\$ 492.58	+Info

C14

Code	Name	Price	+
442199	IVOCLAR VIVADENT - IPS e.max CAD LT Block, Shade A1 Size C14, 5/Pk. is a system ...	US\$ 235.64	+Info
442191	IVOCLAR VIVADENT - IPS e.max CAD LT Block, Shade A2 Size C14, 5/Pk. is a system ...	US\$ 208.04	+Info
442194	IVOCLAR VIVADENT - IPS e.max CAD LT Block, Shade A3 Size C14, 5/Pk. is a system ...	US\$ 233.96	+Info
442211	IVOCLAR VIVADENT - IPS e.max CAD LT Block, Shade A3.5 Size C14, 5/Pk. is a syste...	US\$ 237.59	+Info
442252	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade A4 Size C14, Pack...	US\$ 232.74	+Info
442221	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade B1 Size C14, Pack...	US\$ 230.68	+Info
442268	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade B2 Size C14, Pack...	US\$ 212.59	+Info
442374	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade B3 Size C14, Pack...	US\$ 271.82	+Info
442285	IVOCLAR VIVADENT - IPS e.max CAD for cerec LT Block, Shade BL2 Size C14 5/pk - #...	US\$ 192.57	+Info
442300	IVOCLAR VIVADENT - IPS e.max CAD for cerec LT Block, Shade BL3 Size C14 5/pk - #...	US\$ 216.12	+Info
442321	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade BL4 Size C14, Pac...	US\$ 229.29	+Info
442236	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade C1 Size C14, Pack...	US\$ 252.25	+Info

444265	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade C2 Size A16(S) 5/...	US\$ 422.97	+Info
442228	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade C2 Size C14, Pack...	US\$ 229.11	+Info
442238	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade C3 Size C14, Pack...	US\$ 236.09	+Info
442292	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade C4 Size C14, Pack...	US\$ 238.15	+Info
442233	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade D2 Size C14, Pack...	US\$ 235.65	+Info
442243	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade D3 Size C14, Pack...	US\$ 223.91	+Info
442337	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar LT Block, Shade D4 Size C14, Pack...	US\$ 230.73	+Info

INLAB LT

Code	Name	Price	+
442209	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab LT Block, Shade A2 Size I12, ...	US\$ 218.76	+Info
442226	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab LT Block, Shade A3 Size I12, ...	US\$ 234.64	+Info
442363	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab LT Block, Shade A3.5 Size I12...	US\$ 237.04	+Info
442443	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab LT Block, Shade A4 Size I12, ...	US\$ 265.20	+Info
442340	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab LT Block, Shade B1 Size I12, ...	US\$ 279.50	+Info
442442	IVOCLAR VIVADENT - IPS e.max CAD for CEREC and inLab - LT Block, Shade BL1 Size ...	US\$ 229.79	+Info
442478	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab LT Block, Shade C1 Size I12, ...	US\$ 285.93	+Info
442586	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab LT Block, Shade C3 Size I12, ...	US\$ 270.70	+Info
442585	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab LT Block, Shade D2 Size I12, ...	US\$ 233.93	+Info
444300	IVOCLAR VIVADENT - IPS e.max CAD for CEREC / inLab LT Block, Shade D4 Size I12, ...	US\$ 0.00	+Info

MO

A14(L)

Code	Name	Price	+
442331	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar MO Block, Shade 1 Size A14(L) 5/P...	US\$ 331.49	+Info
442475	IVOCLAR VIVADENT - IPS e.max CAD Ivoclar MO Shade 3, Size A14(L) 5/Pk. Medium Op...	US\$ 473.26	+Info

A14(S)

Code	Name	Price	+
444290	IVOCLAR VIVADENT - IPS e.max CAD CEREC Ivoclar MO Block, Shade 2 Size A14(S) 5/P...	US\$ 296.06	+Info

MT

Code	Name	Price	+
442216	IVOCLAR - IPS e.max CAD MT(Medium Translucency) ceramic milling Blocks 5pk - Sh...	US\$ 223.81	+Info
442201	IVOCLAR VIVADENT - IPS e.max CAD MT(Medium Translucency) ceramic milling Blocks,...	US\$ 234.73	+Info
442206	IVOCLAR VIVADENT - IPS e.max CAD MT(Medium Translucency) ceramic milling Blocks,...	US\$ 220.61	+Info
442588	IVOCLAR VIVADENT - IPS e.max CAD MT(Medium Translucency) ceramic milling Blocks,...	US\$ 270.32	+Info
442344	IVOCLAR VIVADENT - IPS e.max CAD MT(Medium Translucency) ceramic milling Blocks,...	US\$ 235.30	+Info
442368	IVOCLAR VIVADENT - IPS e.max CAD MT(Medium Translucency) ceramic milling Blocks,...	US\$ 253.41	+Info

IPS Empress CAD Leucite - IVOCLAR

IPS Empress CAD are leucite-reinforced glass-ceramic blocks for the CAD/CAM application. They are suitable for the fabrication of full-contour single-tooth restorations.



HT (High Trans)

HT (High Trans)

Code	Name	Price	+
442328	IPS Empress CAD HT Block, Shade A3	US\$ 119.63	+Info
442333	IPS Empress CAD HT Block, Shade A2	US\$ 144.23	+Info
442377	IPS Empress CAD HT Block, Shade A1	US\$ 112.98	+Info
442465	IPS Empress CAD LT Block, Shade B1	US\$ 95.71	+Info
442543	IPS Empress CAD HT (High Translucen	US\$ 148.12	+Info
442562	IPS Empress CAD CEREC / inLab HT bl	US\$ 0.00	+Info
442617	IPS Empress CAD HT (High Translucen	US\$ 220.99	+Info
444257	IPS Empress CAD CEREC / inLab HT bl	US\$ 188.50	+Info
444258	IPS Empress CAD HT(High Translucenc	US\$ 179.40	+Info
444280	IPS Empress CAD HT Block, Shade C2	US\$ 129.79	+Info
444283	IPS Empress CAD CEREC / inLab HT bl	US\$ 146.83	+Info

LT (Low Trans)

LT (Low Trans)

Code	Name	Price	+
442295	IPS Empress CAD LT Block, Shade A2	US\$ 138.18	+Info
442362	IPS Empress CAD LT Block, Shade A3	US\$ 149.16	+Info
442422	IPS Empress CAD LT(Low Translucency	US\$ 126.61	+Info
442439	IPS Empress CAD LT Block, Shade A1	US\$ 77.86	+Info
442451	IPS Empress CAD LT Block, Shade A2	US\$ 130.10	+Info
442502	IPS Empress CAD CEREC / inLab LT bl	US\$ 127.05	+Info
442513	IPS Empress CAD LT Block, Shade A3	US\$ 122.23	+Info

442525	IPS Empress CAD CEREC / inLab LT bl	US\$ 142.93	+Info
444274	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info

Multi Multi

Code	Name	Price	+
442224	IPS Empress CAD CEREC Ivoclar Multi	US\$ 161.30	+Info
442239	IPS Empress CAD CEREC Ivoclar Multi	US\$ 152.87	+Info
442271	IPS Empress CAD CEREC Ivoclar Multi	US\$ 188.31	+Info
442302	IPS Empress CAD CEREC Ivoclar Multi	US\$ 182.90	+Info
442327	IPS Empress CAD CEREC Ivoclar Multi	US\$ 184.98	+Info
442348	IPS Empress CAD Multi Block, Shade	US\$ 170.51	+Info
442364	IPS Empress CAD CEREC / inLab Multi	US\$ 168.54	+Info
442382	IPS Empress CAD Multi Block, Shade	US\$ 141.57	+Info
442392	IPS Empress CAD CEREC Ivoclar Multi	US\$ 150.33	+Info
442411	IPS Empress CAD CEREC / inLab Multi	US\$ 176.54	+Info
442413	IPS Empress CAD CEREC / inLab Multi	US\$ 126.10	+Info
442608	IPS Empress CAD CEREC / inLab Multi	US\$ 181.91	+Info

HYBRID CAD/CAM MILLING BLOCKS



SUPER CRYSTAL HYBRID CERAMIC BLOCKS - DIGITECH

SUPER CRYSTAL HYBRID CERAMIC CEREC BLOCKS - DIGITECH

Features:

- * Suitable for CAD/CAM system & chairside system
- * Available in multi-layer shade



HYBRID MULTILAYER BLOCKS C14 (CEREC) - SUPER CRYSTAL

Code	Name	Price	+
124259	CRYSTAL 18x15x13 ML A1	US\$ 9.50	+Info
124260	HYBRID 18x15x13 ML C14 A2	US\$ 9.50	+Info

124261	HYBRID 18x15x13 ML A3	US\$ 9.50	+Info
124262	HYBRID 18x15x13 ML B2	US\$ 9.50	+Info
464380	CRYSTAL 18x15x13 ML C14 B1	US\$ 9.50	+Info
464381	HYBRID SC 18x15x13 C14 ML BL1	US\$ 9.50	+Info
464882	HYBRID 18x15x13 ML C14 B2	US\$ 9.21	+Info
464883	HYBRID 18x15x13 ML A3.5	US\$ 9.50	+Info
473862	HYBRID SC 18x15x13 C14 ML BL2	US\$ 9.50	+Info
473863	HYBRID SC 18x15x13 C14 ML BL3	US\$ 9.50	+Info

HYBRID MULTILAYER BLOCKS B40 (CEREC) - SUPER CRYSTAL

Code	Name	Price	+
464382	HYBRID 40mm (B40)ML BL1	US\$ 21.00	+Info
464386	HYBRID 40mm B40 x Cerec ML A2	US\$ 21.00	+Info
464387	HYBRID SC 40mm -B40-for Cerec ML B1	US\$ 21.00	+Info

CERASMART UNIVERSAL - GC AMERICA



HT (High translucency)

SIZE 12

Code	Name	Price	+
421733	GC CERASMART Universal Size 12, Shade A1, HT block, pack of 5 blocks - #008653	US\$ 158.93	+Info
421592	GC CERASMART Universal size #12, Shade A2, HT block, pack of 5 blocks - #008654	US\$ 158.93	+Info
421591	GC CERASMART Universal size #12, Shade A3, HT block, pack of 5 blocks - #008655	US\$ 158.93	+Info
421722	GC CERASMART Universal Size 12, Shade A3.5, HT block, pack of 5 blocks - #008656	US\$ 158.93	+Info
421966	GC CERASMART Universal Size 12, Shade B1, HT block, pack of 5 blocks - #008657	US\$ 158.93	+Info
422236	GC CERASMART Universal Size 12, Shade BL (Bleach) block, pack of 5 blocks - #008663	US\$ 158.93	+Info

SIZE 14

Code	Name	Price	+
421660	GC CERASMART Universal size #14, Shade A1, HT block, pack of 5 blocks - #008664	US\$ 169.45	+Info
422241	GC CERASMART Universal Size 14, Shade A3.5, HT block, pack of 5 blocks - #008667	US\$ 169.45	+Info
422019	GC CERASMART Universal Size 14, Shade BL (Bleach), pack of 5 blocks with Univers...	US\$ 169.45	+Info
421738	GC CERASMART Universal Size 14, Shade B1, HT block, pack of 5 blocks - #008668	US\$ 163.40	+Info

SIZE 14L

Code	Name	Price	+
------	------	-------	---

421954	GC CERASMART Universal Size 14L, Shade A1, HT block, pack of 5 blocks - #008675	US\$ 188.93	+Info
421915	GC CERASMART Universal Size 14L, Shade A2, HT block, pack of 5 blocks - #008676	US\$ 188.93	+Info
421603	GC CERASMART Universal size #14L, Shade A3, HT block, pack of 5 blocks - #008677	US\$ 188.93	+Info
421704	GC CERASMART Universal size #14L, Shade A3.5, HT block, pack of 5 blocks - #008678	US\$ 188.93	+Info
422239	GC CERASMART Universal Size 14L, Shade B1, HT block, pack of 5 blocks - #008679	US\$ 188.93	+Info

LT (Low translucency)

SIZE 12

Code	Name	Price	+
421849	GC CERASMART Universal Size 12, Shade A1, LT block, pack of 5 blocks - #008658	US\$ 158.93	+Info
422242	GC CERASMART Universal Size 12, Shade A2, LT block, pack of 5 blocks - #008659	US\$ 158.93	+Info
422237	GC CERASMART Universal Size 12, Shade A3, LT block, pack of 5 blocks - #008660	US\$ 158.93	+Info
422246	GC CERASMART Universal Size 12, Shade A3.5, LT block, pack of 5 blocks - #008661	US\$ 158.93	+Info
422024	GC CERASMART Universal Size 12, Shade B1, LT block, pack of 5 blocks - #008662	US\$ 158.93	+Info

SIZE 14

Code	Name	Price	+
421712	GC CERASMART Universal Size 14, Shade A1, LT block, pack of 5 blocks - #008669	US\$ 169.45	+Info
421710	GC CERASMART Universal size #14, Shade A2, LT block, pack of 5 blocks - #008670	US\$ 169.45	+Info
421711	GC CERASMART Universal Size 14, Shade A3, LT block, pack of 5 blocks - #008671	US\$ 169.45	+Info
421899	GC CERASMART Universal Size 14, Shade A3.5, LT block, pack of 5 blocks - #008672	US\$ 163.40	+Info
422238	GC CERASMART Universal Size 14, Shade B1, LT block, pack of 5 blocks - #008673	US\$ 169.45	+Info

SIZE 14L

Code	Name	Price	+
422243	GC CERASMART Universal Size 14L, Shade A1, LT block, pack of 5 blocks - #008680	US\$ 188.93	+Info
421946	GC CERASMART Universal Size 14L, Shade A2, LT block, pack of 5 blocks - #008681	US\$ 182.18	+Info
422245	GC CERASMART Universal Size 14L, Shade A3, LT block, pack of 5 blocks - #008682	US\$ 188.93	+Info
422240	GC CERASMART Universal Size 14L, Shade A3.5, LT block, pack of 5 blocks - #008683	US\$ 188.93	+Info
422014	GC CERASMART Universal Size 14L, Shade B1, LT block, pack of 5 blocks - #008684	US\$ 188.93	+Info
422244	GC CERASMART Universal Size 14L, Shade BL (bleach), pack of 5 blocks - #008685	US\$ 188.93	+Info

HYBRID Cad/Cam Discs - AMANN GIRRBACH

HYBRID Cad/Cam Discs - AMANN GIRRBACH



Code	Name	Price	+
------	------	-------	---

468337	VITA ENAMIC, EM-14, 0M1-T (5-Stk)	US\$ 219.16	+Info
468338	VITA ENAMIC, EM-14, 1M1-T (5-Stk)	US\$ 219.16	+Info
468339	VITA ENAMIC, EM-14, 1M2-T (5-Stk)	US\$ 219.16	+Info
468340	VITA ENAMIC, EM-14, 2M2-T (5-Stk)	US\$ 219.16	+Info
468341	VITA ENAMIC, EM-14, 3M2-T (5-Stk)	US\$ 219.16	+Info
468342	VITA ENAMIC, EM-14, 0M1-HT (5-Stk)	US\$ 219.16	+Info
468343	VITA ENAMIC, EM-14, 1M1-HT (5-Stk)	US\$ 219.16	+Info
468344	VITA ENAMIC, EM-14, 1M2-HT (5-Stk)	US\$ 219.16	+Info
468345	VITA ENAMIC, EM-14, 2M2-HT (5-Stk)	US\$ 219.16	+Info
468346	VITA ENAMIC, EM-14, 3M2-HT (5-Stk)	US\$ 219.16	+Info
468359	VITA ENAMIC multiColor, 14, 1M1-HT(	US\$ 242.08	+Info
468360	VITA ENAMIC multiColor, 14, 1M2-HT(	US\$ 242.08	+Info
468361	VITA ENAMIC multiColor, 14, 2M2-HT(	US\$ 242.08	+Info
468362	VITA ENAMIC multiColor, 14, 3M2-HT(	US\$ 242.08	+Info
468363	VITA ENAMIC multiColor, 14, 4M2-HT(	US\$ 242.08	+Info

GUM / GINGIVA PMMA DISCS



UNIDESA - PMMA PINK DISC

UNIDESA - PMMA PINK DISC

Features:

- * For CAD-CAM milling machines
- * 100% biocompatible and antiallergenic
- * Permanent natural coloration (does not change once it has been processed)
- * Residual monomer Free
- * High resistance to abrasion
- * Easy to finish and polish
- * Made with PMMA resins of a high molecular weight.
- * Super High Cross Linked
- * High Speed milling

#USDDUNIDENT



98.5mm

Code	Name	Price	+
123237	UNIDESIA - PMMA 98.5mm / 20mm / GUM PINK Disc .	US\$ 49.39	+Info
123444	UNIDESIA - PMMA 98.5mm / 20mm / GUM Veined PINK Disc . 4250 / PV	US\$ 48.89	+Info
123245	UNIDESIA - PMMA 98.5mm / 25mm / GUM PINK Disc . 4225 / P	US\$ 50.14	+Info
123445	UNIDESIA - PMMA 98.5mm / 25mm GUM Veined PINK Disc . 4225 / PV	US\$ 46.94	+Info
417324	UNIDESIA - PMMA 98.5mm / 30mm GUM Veined PINK Disc . 4230 / PV	US\$ 34.03	+Info

95mm

Code	Name	Price	+
123261	UNIDESIA - PMMA ZZ 95mm / 16mm / GUM PINK Disc for - for ZirkonZahn	US\$ 30.94	+Info
123262	UNIDESIA - PMMA ZZ 95mm / 16mm / GUM PINK veined Disc for - for ZirkonZahn	US\$ 35.94	+Info
123253	UNIDESIA - PMMA ZZ 95mm / 20mm / GUM PINK Disc for - for ZirkonZahn	US\$ 50.40	+Info
123254	UNIDESIA - PMMA ZZ 95mm / 20mm / GUM PINK veined Disc for - for ZirkonZahn	US\$ 50.40	+Info
123442	UNIDESIA - PMMA ZZ 95mm / 25mm / GUM PINK Disc for - for ZirkonZahn	US\$ 37.34	+Info
123443	UNIDESIA - PMMA ZZ 95mm / 25mm / GUM Veined PINK Disc for - for ZirkonZahn 4225Z/PV	US\$ 50.67	+Info
417325	UNIDESIA - PMMA ZZ 95mm / 30mm / GUM Veined PINK Disc for .4230Z / PV - for Zirko...	US\$ 51.87	+Info

93 mm Amann

Code	Name	Price	+
123269	UNIDESIA - PMMA 93mm / 71mm / 20mm / GUM PINK Disc for Ceramill . Amann Girrbach	US\$ 31.90	+Info
123270	UNIDESIA - PMMA 93mm / 71mm / 20mm / GUM PINK Veined Disc for Ceramill . Amann Gi...	US\$ 34.67	+Info

GUM PMMA DISCS - DIGILUX

GUM PMMA DISCS - DIGILUX

Features:

Our GUM PMMA discs are perfect for dental laboratories who want to reduce labor costs and improve consistency by using modern milling systems. They offer excellent physical and biocompatibility characteristics due to highly cross-linked polymers and result in superior surface finish.



GUM

98.5mm

Code	Name	Price	+
124185	DIGILUX PMMA Disc 98.5mm (98mm) x 20mm Color : GUM - REGULAR PINK Color #124185	US\$ 38.00	+Info
124186	DIGILUX PMMA Disc 98.5mm (98mm) x 25mm Color : GUM - REGULAR PINK -Color #124186	US\$ 38.00	+Info
401074	DIGILUX PMMA Disc 98.5mm (98mm) x 25mm Color : GUM - LIGHT PINK Color #401074	US\$ 29.68	+Info
416525	DIGILUX PMMA Disc 98.5mm /30mm GUM - REGULAR PINK #416525	US\$ 44.16	+Info
419992	DIGILUX PMMA Disc 98.5mm /30mm GUM LIGHT PINK Blank () for /Wieland/Opend . #419992	US\$ 45.08	+Info

95mm

GUM COLOR

GUM COLOR

Code	Name	Price	+
124227	PMMA 95x20mm GUM LIGHT	US\$ 24.00	+Info
124228	PMMA 95x25mm GUM REGULAR	US\$ 21.36	+Info
419994	PMMA 95x25mm GUM LIGHT PINK	US\$ 23.04	+Info

Other Models

Other Models

Code	Name	Price	+
473712	PMMA 95x25mm GUM REGULAR PINK	US\$ 22.80	+Info

SPLINT CAD/CAM DISCS



BEGO PMMA Splint E



Code	Name	Price	+
474642	BEGO PMMA Splint E 20 x 98.5	US\$ 95.00	+Info
474643	BEGO PMMA Splint E 16 x 98.5	US\$ 95.00	+Info

Splint Cad/Cam Discs - ASTRON

Splint Cad/Cam Discs - ASTRON

Code	Name	Price	+
420878	Clear Splint Disc 98x20	US\$ 76.08	+Info
420879	95mm SPLINT /20mm/Clear Clear	US\$ 72.45	+Info

Splint Cad/Cam Discs - UNIDESA

Splint Cad/Cam Discs - UNIDESA

Code	Name	Price	+
124330	98.5mm/20mm/Clear GREEN	US\$ 41.09	+Info
417321	98.5mm/25mm/Clear VIOLET	US\$ 45.66	+Info

CASTING PMMA MILLING DISCS



CASTING PMMA DISC - UNIDESA

CLEAR & CASTING PMMA DISC - UNIDESA

Features:

- * For CAD-CAM milling machines
- * 100% biocompatible and antiallergenic
- * Permanent natural coloration (does not change once it has been processed)
- * Residual monomer Free
- * High resistance to abrasion
- * Easy to finish and polish
- * Made with PMMA resins of a high molecular weight.
- * Super High Cross Linked
- * High Speed milling



#USDDUNIDENT

98.5mm

Code	Name	Price	+
123247	UNIDESA - PMMA 98.5mm/16mm/CLEAR Casting Disc 1pk - Mfg #4252/C	US\$ 42.52	+Info
123239	UNIDESA - PMMA 98.5mm/20mm/CLEAR Casting Disc 1pk - Mfg #4250/C	US\$ 42.52	+Info
400110	UNIDESA - PMMA 98.5mm/25mm/CLEAR Casting Disc 1pk . 4225/C	US\$ 42.52	+Info

95mm

Code	Name	Price	+
123263	UNIDESA - PMMA ZZ 95mm/16mm/CLEAR Casting / Splint Disc 1pk - Mfg #4252Z/C - fo...	US\$ 43.78	+Info
123255	UNIDESA - PMMA ZZ 95mm/20mm/CLEAR Casting / Splint Disc 1pk - Mfg #4250Z/C - fo...	US\$ 45.60	+Info

400007

UNIDESA - PMMA ZZ 95mm/25mm/CLEAR Casting / Splint Disc 1pk- Mfg #4225Z/C - for...

US\$ 45.60

[+Info](#)

101mm Ammann Girschbach

Code	Name	Price	+
123271	UNIDESA - PMMA 93mm / 71mm /20mm/CLEAR Casting Disc 1pk for Ceramill - Mfg #4250...	US\$ 24.88	+Info

CASTING PMMA MILLING DISCS - DIGILUX

CLEAR PMMA DISCS - DIGILUX

Features:

Our Clear PMMA discs are perfect for dental laboratories who want to reduce labor costs and improve consistency by using modern milling systems. They offer excellent physical and biocompatibility characteristics due to highly cross-linked polymers and result in superior surface finish.



DIGILUX



CASTING

CLEAR 93mm AMANN

MONO CLEAR

MONO CLEAR

Code	Name	Price	+
113859	PMMA 93x16mm CLEAR	US\$ 10.92	+Info
124248	PMMA 93x20mm CLEAR	US\$ 22.86	+Info

Other Models

Other Models

Code	Name	Price	+
474286	PMMA 93x13mm CLEAR	US\$ 10.40	+Info

CLEAR 95mm

Code	Name	Price	+
124229	DIGILUX PMMA Disc ZZ 95mm x 16mm Color : CLEAR - 1pk #124229	US\$ 21.91	+Info
124230	DIGILUX PMMA Disc ZZ 95mm x 20mm Color : CLEAR - 1pk #124230	US\$ 23.78	+Info
124231	DIGILUX PMMA Disc ZZ 95mm x 25mm Color : CLEAR - #124231	US\$ 32.39	+Info

CLEAR 98.5mm

Code	Name	Price	+
124187	DIGILUX PMMA Disc 98.5mm (98mm) x 16mm Color : CLEAR - 1pk #124187	US\$ 24.77	+Info
124188	DIGILUX PMMA Disc 98.5mm (98mm) x 20mm Color : CLEAR - 1pk #124188	US\$ 24.77	+Info
124189	DIGILUX PMMA Disc 98.5mm (98mm) x 25mm Color : CLEAR - 1pk #124189	US\$ 32.39	+Info

RED PATTERN PMMA DISCS- DIGILUX

RED PATTERN PMMA DISCS - DIGILUX

Features:

Our Red Pattern PMMA discs are perfect for dental laboratories who want to reduce labor costs and improve consistency by using modern milling systems. They offer excellent physical and biocompatibility characteristics due to highly cross-linked polymers and result in superior surface finish.



95mm

Code	Name	Price	+
124232	DIGILUX PMMA Disc PATTERN ZZ 95mm x 16mm Color : RED - for ZirkonZahn #124232	US\$ 19.72	+Info
124233	DIGILUX PMMA Disc PATTERN ZZ 95mm x 20mm Color : RED - for ZirkonZahn #124233	US\$ 16.00	+Info

98.5mm

Code	Name	Price	+
124190	DIGILUX PMMA Disc PATTERN 98.5mm (98mm) x 16mm Color : RED 1pk - #124190	US\$ 15.84	+Info
124191	DIGILUX PMMA Disc PATTERN 98.5mm (98mm) x 20mm Color : RED 1pk - #124191	US\$ 15.84	+Info

CASTING PMMA MILLING DISCS - VIP

PMMA Casting & Gum Color MILLING DISCS - VIP



Code	Name	Price	+
113857	PMMA Cast Clear Roland 98x18	US\$ 12.80	+Info
114055	PMMA Cast Zirkon 95x16mm Clear	US\$ 22.40	+Info

MILLING DISCS & BLOCKS - DIGITECH

MILLING DISCS & BLOCKS - DIGITECH

Code	Name	Price	+
113535	Zirconia Disc HS 98mm x 16mm	US\$ 59.00	+Info

113536	Zirconia Disc HS 98mm x 18mm	US\$ 58.91	+Info
115488	Zirconia Disc HS 98mm x 12mm	US\$ 39.10	+Info

FIBER GLASS CAD/CAM DISCS



TRILOR HPP

TRILOR HPP

Bone

Bone

95mm

95mm

16mm

16mm

Code	Name	Price	+
475312	TriLor HPP, 95 x 16mm, Bone	US\$ 243.46	+Info

20mm

20mm

Code	Name	Price	+
475314	TriLor HPP, 95 x 20mm, Bone	US\$ 253.02	+Info

25mm

25mm

Code	Name	Price	+
475316	TriLor HPP, 95 x 25mm, Bone	US\$ 291.29	+Info

98.5mm

98.5mm

14mm

14mm

Code	Name	Price	+
475302	TriLor HPP, 98.5 x 14mm, Bone	US\$ 243.46	+Info

16mm

16mm

Code	Name	Price	+
475303	TriLor HPP, 98.5 x 16mm, Bone	US\$ 243.46	+Info

20mm

20mm

Code	Name	Price	+
475305	TriLor HPP, 98.5 x 20mm, Bone	US\$ 253.02	+Info

25mm

25mm

Code	Name	Price	+
475307	TriLor HPP, 98.5 x 25mm, Bone	US\$ 291.29	+Info

Pink

Pink

95mm

95mm

16mm

16mm

Code	Name	Price	+
475313	TriLor HPP, 95 x 16mm, Pink	US\$ 247.80	+Info

20mm

20mm

Code	Name	Price	+
475315	TriLor HPP, 95 x 20mm, Pink	US\$ 253.02	+Info

98.5mm

98.5mm

16mm

16mm

Code	Name	Price	+
475304	TriLor HPP, 98.5 x 16mm, Pink	US\$ 247.80	+Info

20mm

20mm

Code	Name	Price	+
475306	TriLor HPP, 98.5 x 20mm, Pink	US\$ 261.71	+Info

TRILOR HPP Arch

TRILOR HPP Arch

Code	Name	Price	+
475308	TriLor HPP Arch, 3.5mm (3pk)	US\$ 254.13	+Info
475309	TriLor HPP Arch, 5.5mm (3pk)	US\$ 294.40	+Info
475310	TriLor HPP Arch, 7.5mm (3pk)	US\$ 321.24	+Info
475311	TriLor HPP Arch Kit 7.5mm)	US\$ 265.76	+Info

FLEXIBLE MILLING DISCS

Unbreakable™



iFlex Milling Discs

iFlex

Light Dark Pink

Light Dark Pink

98x20mm

98x20mm

Code	Name	Price	+
477156	iFlex by TCS - Disc / 98mmx20mm / S	US\$ 64.68	+Info

98x25mm

98x25mm

Code	Name	Price	+
477158	iFlex by TCS - Disc / 98mmx25mm / S	US\$ 64.68	+Info

Natural

Natural

98x20mm

98x20mm

Code	Name	Price	+
477155	iFlex by TCS - Disc / 98mmx20mm / S	US\$ 64.68	+Info

98x25mm

98x25mm

Code	Name	Price	+
477159	iFlex by TCS - Disc / 98mmx25mm / S	US\$ 64.68	+Info

Standard Pink

Standard Pink

98x20mm

98x20mm

Code	Name	Price	+
477154	iFlex by TCS - Disc / 98mmx20mm / S	US\$ 64.68	+Info

98x25mm

98x25mm

Code	Name	Price	+
477157	iFlex by TCS - Disc / 98mmx25mm / S	US\$ 64.68	+Info

Unbreakable Milling Discs

Unbreakable

Light Dark Pink

Light Dark Pink

98x20mm

98x20mm

Code	Name	Price	+
477150	TCS Unbreakable - Disc / 98mmx20mm	US\$ 64.68	+Info

98x25mm

98x25mm

Code	Name	Price	+
477153	TCS Unbreakable - Disc / 98mmx25mm	US\$ 64.68	+Info

Natural

Natural

98x20mm

98x20mm

Code	Name	Price	+
477149	TCS Unbreakable - Disc / 98mmx20mm	US\$ 64.68	+Info

98x25mm

98x25mm

Code	Name	Price	+
477151	TCS Unbreakable - Disc / 98mmx25mm	US\$ 64.68	+Info

Standard Pink

Standard Pink

98x20mm

98x20mm

Code	Name	Price	+
477148	TCS Unbreakable - Disc / 98mmx20mm	US\$ 64.68	+Info

98x25mm

98x25mm

Code	Name	Price	+
477152	TCS Unbreakable - Disc / 98mmx25mm	US\$ 64.68	+Info

MILLING DISCS & BLOCKS - IVOCAR VIVADENT

MILLING DISCS & BLOCKS - IVOCAR VIVADENT

E.Max

e.Max

Code	Name	Price	+
442457	IPS e.max CAD CEREC Ivoclar LT Bloc	US\$ 451.43	+Info
442479	IPS e.max CAD CEREC Ivoclar HT Bloc	US\$ 230.49	+Info
442516	IPS e.max CAD CEREC Ivoclar LT Bloc	US\$ 515.35	+Info
442673	IPS e.max CAD CEREC Ivoclar MO Bloc	US\$ 477.04	+Info
444259	IPS e.max CAD CEREC Ivoclar LT Bloc	US\$ 473.03	+Info
444278	IPS e.max CAD CEREC Ivoclar LT Bloc	US\$ 450.79	+Info

IPS E-MAX

IPS E-MAX

Code	Name	Price	+
442249	IPS e.max CAD for CEREC / inLab LT	US\$ 246.39	+Info
444250	IPS e.max CAD for CEREC / inLab LT	US\$ 0.00	+Info

IPS Empress CAD

IPS Empress CAD

Code	Name	Price	+
444294	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444276	IPS Empress CAD CEREC / inLab HT bl	US\$ 0.00	+Info
444277	IPS Empress CAD CEREC / inLab HT bl	US\$ 0.00	+Info
444279	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444282	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444284	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info

444285	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444286	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444288	IPS Empress CAD CEREC / inLab HT bl	US\$ 0.00	+Info
444289	IPS Empress CAD CEREC / inLab HT bl	US\$ 145.53	+Info
444292	IPS Empress CAD CEREC / inLab HT bl	US\$ 0.00	+Info
444275	IPS Empress CAD CEREC / inLab HT bl	US\$ 0.00	+Info
444296	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444297	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444298	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444299	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444301	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444302	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444303	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444304	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444305	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444254	IPS Empress CAD CEREC / inLab HT bl	US\$ 0.00	+Info
442560	IPS Empress CAD CEREC / inLab HT bl	US\$ 0.00	+Info
442619	IPS Empress CAD CEREC / inLab HT bl	US\$ 149.43	+Info
442620	IPS Empress CAD LT Block, Shade B1	US\$ 118.43	+Info
444247	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444248	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444249	IPS Empress CAD CEREC / inLab HT bl	US\$ 0.00	+Info
444251	IPS Empress CAD CEREC / inLab HT bl	US\$ 0.00	+Info
444252	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444253	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
442416	IPS Empress CAD CEREC / inLab Multi	US\$ 185.63	+Info
444255	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444256	IPS Empress CAD LT(Low Translucency	US\$ 179.40	+Info
444260	IPS Empress CAD CEREC / inLab HT bl	US\$ 0.00	+Info
444263	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444266	IPS Empress CAD CEREC / inLab HT bl	US\$ 0.00	+Info
444268	IPS Empress CAD CEREC / inLab HT bl	US\$ 0.00	+Info
444269	IPS Empress CAD CEREC / inLab HT bl	US\$ 0.00	+Info
444270	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info
444271	IPS Empress CAD CEREC / inLab LT bl	US\$ 0.00	+Info

Other Models

Other Models

Code	Name	Price	+
442587	IPS e.max CAD for CEREC / inLab HT	US\$ 0.00	+Info
444178	IPS e.max CAD for CEREC / inLab HT	US\$ 279.50	+Info
444262	IPS e.max CAD for CEREC / inLab HT	US\$ 0.00	+Info
444295	IPS e.max CAD for CEREC / inLab HT	US\$ 0.00	+Info

PEEK CAD/CAM DISCS



PEEK - DIGITECH PEEK - DIGITECH



93mm

Natural

Code	Name	Price	+
416703	DIGITECH - PEEK AMANN GIRRBACH 93mm x 13mm NATURAL #PEEK AG 13MM	US\$ 76.86	+Info
416704	NAT DIGITECH - PEEK AMANN GIRRBACH 93mm x 20mm NATURAL #PEEK AG 20MM	US\$ 139.20	+Info
	NAT		

95 mm

Natural

Code	Name	Price	+
420153	DIGITECH - PEEK 95mm x 14mm NATURAL #PEEK 95MM 14NAT	US\$ 108.80	+Info
420154	DIGITECH - PEEK 95mm x 20mm NATURAL #PEEK 95MM 20MM NAT	US\$ 189.00	+Info

Pink

Code	Name	Price	+
420155	DIGITECH - PEEK 95mm x 20mm PINK #PEEK 95MM 20MM PINK	US\$ 199.50	+Info

98.5mm


Natural

Code	Name	Price	+
416710	DIGITECH - PEEK 98.5mm x 14mm NATURAL #PEEK 98.5MM 14 NAT	US\$ 159.60	+Info
416711	DIGITECH - PEEK 98.5mm x 20mm NATURAL #PEEK 98.5MM 20MM NAT	US\$ 136.00	+Info

Pink


Code	Name	Price	+
416712	DIGITECH - PEEK 98.5mm x 20mm PINK #PEEK 98.5MM 20MM PINK	US\$ 199.50	+Info

TITANIUM CAD/CAM DISCS

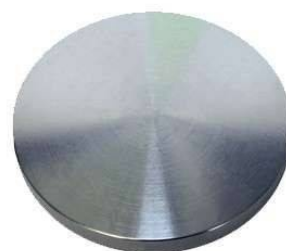

Cad/Cam Titanium Discs - ULTIMATE

Cad/Cam Titanium Discs - ULTIMATE

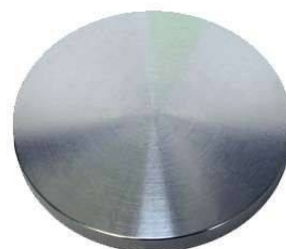

Code	Name	Price	+
471469	Soft Metal Blck 98.5mm x 12mm	US\$ 156.00	+Info
474835	Soft Metal Blck 98.5mm x 14mm	US\$ 156.00	+Info
420332	Soft Metal Blck 98.5mm x 20mm	US\$ 192.00	+Info

MEDILOY M-Ti4 Titanium Discs - BEG



Code	Name	Price	+
474629	Mediloy M-Ti4 12 mm mit Schulter	US\$ 151.37	+Info
474630	Mediloy M-Ti4 14 mm mit Schulter	US\$ 184.36	+Info
474631	Mediloy M-Ti4 16 mm mit Schulter	US\$ 209.51	+Info
474632	Mediloy M-Ti4 18 mm mit Schulter	US\$ 197.79	+Info
474633	Mediloy M-Ti4 20 mm mit Schulter	US\$ 234.68	+Info
474634	Mediloy M-Ti4 22 mm mit Schulter	US\$ 234.70	+Info
474635	Mediloy M-Ti4 25 mm mit Schulter	US\$ 273.02	+Info

MEDILOY M-Ti5 Titanium Discs - BEG



Code	Name	Price	+
474636	Mediloy M-Ti5 12 mm mit Schulter	US\$ 162.01	+Info
474637	Mediloy M-Ti5 14 mm mit Schulter	US\$ 184.36	+Info
474639	Mediloy M-Ti5 18 mm mit Schulter	US\$ 197.79	+Info
474640	Mediloy M-Ti5 20 mm mit Schulter	US\$ 234.68	+Info
474641	Mediloy M-Ti5 22 mm mit Schulter	US\$ 221.27	+Info
474638	Mediloy M-Ti 5 16 mm mit Schulter	US\$ 209.51	+Info

HATHO POLISHING KIT

Thermoplastics
like *never* before



Packing Diamond D®
Denture Acrylic with
Tecnoflask



Keystone Industries: Laboratory



WAX CAD/CAM DISCS



DIGILUX



WAXING

WAX CAD/CAM DISCS - DIGILUX

WAX CAD/CAM DISCS - DIGILUX

DIGILUX WAX DISC is a special high quality wax for milling, completely calcinable.

The WAX Disc Material allows the reduction in milling times and easy extraction. Ideal for casting and injection process.

Features:

- * Excellent Milling Capacity
- * Low Wear out of tools in CAM millings
- * Biocompatibility
- * Good Contrast

Available in 93mm/95mm/98mm/98.5mm



98.5mm

White Color

Code	Name	Price	+
115413	DIGILUX-Wax Casting Disc 98.5mm x 10mm - WIELAND - 1pk	US\$ 13.86	+Info
113894	DIGILUX Wax Casting Disc 98.5mm (98mm) x 20mm - WIELAND - 1pk	US\$ 14.00	+Info
115277	DIGILUX Wax Tooth Color Castable Disc 98.5mm (98mm) x 25mm - WIELAND - 1pk	US\$ 14.00	+Info

93mm

Tooth Color

Code	Name	Price	+
113918	DIGILUX Wax WHITE Color & Casting Disc 93mm x 71mm x 20mm - AMANN GIRRBACH 1 blo...	US\$ 13.86	+Info

95mm

Tooth Color

Code	Name	Price	+
115276	DIGILUX WAX Tooth Color Disc ZZ 95mm x 25mm - 1 block per box	US\$ 14.00	+Info

VERICORE WAX DISC - WHIP MIX

VERICORE WAX DISC - WHIP MIX

The Vericore Millwax discs are specially formulated to provide exceptional detail for full contour and coping wax forms, but are rigid enough for long span bridge patterns that will not warp. It mills cleanly and adheres to most waxes for adjustments and spruing.



The wax burns out cleanly and leaves no ash residue. Suitable for both casting or pressing techniques.

100 MM X 12 MM

Code	Name	Price	+
120364	WHIP MIX Vericore Wax Disc 100x12mm w/collar - #72432 (Cad/Cam wax discs)	US\$ 36.78	+Info
120363	WHIP MIX Vericore Wax Disc 100x12mm w/collar - #72432	US\$ 36.78	+Info

100 MM X 16 MM

Code	Name	Price	+
120361	WHIP MIX Vericore Wax Disc 100x16mm w/collar - #72434	US\$ 36.78	+Info
120362	WHIP MIX Vericore Wax Disc 100x16mm w/collar - #72434 (Cad/Cam wax discs)	US\$ 36.78	+Info

100 MM X 20 MM

Code	Name	Price	+
120360	WHIP MIX Vericore Wax Disc 100x20mm w/collar - #72436 (Cad/Cam wax discs)	US\$ 36.78	+Info
120359	WHIP MIX Vericore Wax Disc 100x20mm w/collar - #72436	US\$ 36.78	+Info

98 MM X 12 MM

Code	Name	Price	+
120358	WHIP MIX Vericore Wax Disc 98x12mm (98.5) w/collar - #72420 (Cad/Cam wax disc...	US\$ 36.78	+Info

98 MM X 16 MM

Code	Name	Price	+
120356	WHIP MIX Vericore Wax Disc 98x16mm (98.5) w/collar - #72422 (Cad/Cam wax disc...	US\$ 36.78	+Info
120327	WHIP MIX Vericore Wax Disc 98x16mm (98.5) w/collar - #72422	US\$ 36.78	+Info

98 MM X 20 MM

Code	Name	Price	+
120357	WHIP MIX Vericore Wax Disc 98x20mm (98.5) w/collar - #72424 (Cad/Cam wax disc...	US\$ 38.42	+Info
120325	WHIP MIX Vericore Wax Disc 98x20mm (98.5) w/collar - #72424	US\$ 38.42	+Info

WAX DISCS - BESQUAL



CAD/CAM Wax
SIZE: 98 x 18mm

Code	Name	Price	+
470847	Milling Wax Disc 98 x 12mm	US\$ 14.06	+Info
470848	Milling Wax Disc 98 x 10mm	US\$ 13.92	+Info
470849	Milling Wax Disc 98 x 14mm	US\$ 14.50	+Info
470850	Milling Wax Disc 98 x 16mm	US\$ 14.50	+Info
470851	Milling Wax Disc 98 x 18mm	US\$ 14.06	+Info
470852	Milling Wax Disc 98 x 20mm	US\$ 14.50	+Info
470853	Milling Wax Disc 98 x 25mm	US\$ 17.64	+Info

CAD/CAM

CAD/CAM EQUIPMENT

CAD/CAM MILLING UNITS



CAD/CAM Milling Units - UP3D



DRY MILLING UNITS | UP3D

DRY

DISC CHANGER

DISC CHANGER

Code	Name	Price	+
474123	P53DC Smart 5-axis Milling Machine	US\$ 0.00	+Info

NO- DISC CHANGER

NO- DISC CHANGER

Code	Name	Price	+
473287	P52 5 Axis Dental Milling Machine	US\$ 0.00	+Info
473288	P53 Smart 5-axis Milling Machine	US\$ 0.00	+Info
476416	P53 Smart 5-axis Milling Machine	US\$ 0.00	+Info

WET MILLING UNITS | UP3D

WET

Code	Name	Price	+
473290	P42 4-axis Wet Milling Machine	US\$ 0.00	+Info

Milling Units - ROLAND

Dental Milling Units - ROLAND



42W WET MILLINGS - ROLAND



DWX-42W

DWX-42W

Code	Name	Price	+
464630	DWX-42W Wet Mill Only	US\$ 0.00	+Info

DWX-42W PLUS

DWX-42W PLUS

Code	Name	Price	+
420224	DWX-42W Wet Plus Mill Only	US\$ 0.00	+Info

52D DRY MILLINGS

52D DRY MILLING - ROLAND

Features

The next generation of the industry's best-selling dental mill has kept everything users love, and then made it even better. Reliable, precise and affordable, the DWX-52D 5-Axis Dental Milling Machine now offers better disc handling, more material capabilities, improved tool management and other enhanced features, to deliver even more efficiency.

A new Snap-In Clamp System for disc changing includes six material adapters and offers a fast and secure set-up?eliminating mounting errors and minimizing damage during disc insertion. With the Snap-In Clamp System, reinserted material adapters keep discs indexed to their original position for production that requires secondary milling or two-stage processing. In addition, the DWX-52D features built-in storage drawers and an included disc adapter rack to ensure a highly organized and efficient workspace.

The Snap-In Clamp System enables multi-device adapter integration. The same material adapters can be used on the DWX-52D, DWX-52DC and DWX-52DCi for a convenient and scalable milling solution.



Reliable Dental Milling Now Comes with Added Efficiency

DGSHAPE by Roland takes thirty-years of reliability and innovation and combines it with next-level production and efficiency to meet the needs of today's growing lab businesses.

Code	Name	Price	+
420811	ROLAND DWX-52D OPEN BOX / REFURBISHED 5-Axis Dental Milling Machine ONLY 110vol ...	US\$ 0.00	+Info
420210	ROLAND DWX-52D PLUS 5-Axis Dental Milling unit ONLY 110vol - with air tube conne...	US\$ 0.00	+Info

52DCI - DRY MILLINGS + DISC CH ANGER

52DCI - DRY MILLING DISC CHANGER- ROLAND

Features

The ultimate in ease, material flexibility and unattended production, the DWX-52DCi dental milling machine with automatic disc changer expands your lab production and profit. With a 6-slot Automatic Disc Changer, 15-station Automatic Tool Changer and other automated features, it gives users of all skill levels an exceptionally powerful, efficient and affordable dental restoration solution.



Intelligent, Unattended Milling

The DWX-52DCi dental milling machine now comes with performance visualization software and other intelligent updates. Power it on, let it go, and automatically mill numerous dental restoration jobs with complete unattended confidence.

Reliable Milling Now Comes with Upgraded Efficiency

DGSHAPE by Roland takes thirty-years of reliability and innovation and combines it with next-level production to meet the needs of today's growing lab businesses. The new, upgraded DWX-52DCi dental mill features DWINDEX performance visualization software and an advanced dust collection system for a smart and clean operation.

Code	Name	Price	+
420399	ROLAND DWX-52DCI - OPEN BOX / REFURBISHED - Dry Milling unit 110vol with	US\$ 0.00	+Info
420211	ROLAND DWX-52DCi Dry Milling unit 110vol with Automatic Disc Changer NEW #DWX-...	US\$ 0.00	+Info

DXW-53DC DRY MILLINGS DXW-53DC

Code	Name	Price	+
474483	DWX-53DC Mill Auto Disc Ch ONLY	US\$ 0.00	+Info

DXW-53DC+RENFERT

DXW-53DC DRY MILLING + DUST COLLECTOR RENFERT

Code	Name	Price	+
474486	DWX-53DC Mill + Dust collector	US\$ 0.00	+Info
476660	DWX-52D plus Mill + Dust collector	US\$ 0.00	+Info
475026	DWX-53DC Mill + Dust collector	US\$ 0.00	+Info

Milling Units - CEREC

Milling Units - CEREC

Code	Name	Price	+
419965	CEREC MCXL USED	US\$ 10,500.00	+Info

Milling Units - ECCO

Milling Units - ECCO

Code	Name	Price	+
103899	Manual Zircon Milling Sys 110v	US\$ 10,832.40	+Info
112380	MILL PANTOGRAPH 110v	US\$ 4,436.80	+Info

SCANNERS CAD/CAM



DESKTOP SCANNERS

LS100 Desktop Scanners - BLZ

LS100

Code	Name	Price	+
476865	Lab Scanner LS100	US\$ 7,500.00	+Info

Desktop Scanners - UP3D

Desktop Scanners - UP3D



UP1000 OPTICAL SCANNER

UP1000 OPTICAL SCANNER

Code	Name	Price	+
474837	UP1000 Scanner	US\$ 0.00	+Info

UP400 OPTICAL SCANNER

UP400 OPTICAL SCANNER

Code	Name	Price	+
124252	UP400 Scanner	US\$ 0.00	+Info

UP560 OPTICAL SCANNER

UP560 OPTICAL SCANNER

Code	Name	Price	+
470645	UP560 Scanner	US\$ 0.00	+Info
473951	OPEN BOX UP560 Scanner	US\$ 0.00	+Info

Desktop Scanners - AMANN GIRRBACH

Desktop Scanners - AMANN GIRRBACH

Code	Name	Price	+
467354	Calibr. model Ceramill Map 300/400/	US\$ 115.16	+Info
467507	dust cover Map 600	US\$ 295.26	+Info
467677	side panel right Map 600	US\$ 220.54	+Info
467678	side panel left Map 600	US\$ 220.54	+Info
467679	top cover Map 600	US\$ 208.27	+Info
467680	interior Map 600	US\$ 260.96	+Info
467681	front panel incl. glass Map 600	US\$ 1,139.37	+Info
467682	axis cover Map 600	US\$ 156.82	+Info
467683	sponge rubber strip Map 600	US\$ 74.74	+Info
467684	Controller Card Map 600	US\$ 1,417.48	+Info
467685	Motordriver Card Map 600	US\$ 1,417.48	+Info

467686	rocker switch small double-pole Map	US\$ 60.04	+Info
467687	USB grommet Map 600	US\$ 49.00	+Info
467688	power supply MeanWell LRS-50-24 Map	US\$ 149.47	+Info
473727	Ceramill Mind	US\$ 6,607.88	+Info

Desktop Scanners - SHINING

Desktop Scanners - SHINING



DS-EX PRO

DS-EX PRO

Code	Name	Price	+
416559	SHINING 3D DS-EX PRO Scanner C	US\$ 6,320.00	+Info

INTRAORAL SCANNERS



Intraoral Scanners - SHINING

Intraoral Scanner - SHINING



AORALSCAN

AORALSCAN

Code	Name	Price	+
466321	MetiSmile Facial Scanner	US\$ 6,279.00	+Info

AORALSCAN 3

AORALSCAN 3

Code	Name	Price	+
470817	SHINING AORALSCAN 3 Intraoral	US\$ 11,623.50	+Info

IntraOral Scanners | BLZ - DENTAL

BLZ - DENTAL

Code	Name	Price	+
476431	INO200 intraoral-scanner	US\$ 7,500.00	+Info

BLOOMDENT

BLOOMDENT

Code	Name	Price	+
475027	BLOOMDENT AORALSCAN intraoral	US\$ 8,330.00	+Info

Intraoral Scanners - CEREC

Intraoral Scanners - CEREC

Code	Name	Price	+
124545	BLUE CAM CEREC USADA	US\$ 2,775.00	+Info
421106	OMNI CAM CEREC USADA	US\$ 9,075.00	+Info

Intraoral Scanners - DENTERPRISE

Intraoral Scanners - DENTERPRISE

Code	Name	Price	+
420170	QuickScan iOs 3D Instroral scan	US\$ 15,499.90	+Info

3D PRINTERS



3D Printers - ACKURETTA

3D Printers - ACKURETTA

Ready, Set, Print

All-Inclusive SOL LCD 3D Printer



SOL LCD 3D PRINTER

SOL LCD 3D

Code	Name	Price	+
470275	SOL LCD 3D Printer Bundle	US\$ 9,662.64	+Info
474935	SOL LCD 3D Printer	US\$ 5,192.78	+Info
476899	SOL LCD 3D Printer	US\$ 6,509.03	+Info

VARSEO 3D PRINTER - BEGO



Varseo Accessories



Varseo

Varseo

Code	Name	Price	+
419471	Foil - Varseo (Pckg. a 5 St.)	US\$ 189.18	+Info
419480	FILTER MATS (2 PCS) - VARSEO	US\$ 27.65	+Info
419504	Varseo S cartridge A (for CAD/Cast-	US\$ 927.32	+Info
419505	Varseo S cartridge B (for Surgical	US\$ 927.32	+Info
419825	Varseo Accessory Kit	US\$ 298.64	+Info
419826	FILTER FOR VARSEO MATERIALS (with y	US\$ 50.93	+Info
464390	Build Platform	US\$ 52.74	+Info
464391	Flash Module	US\$ 1,368.15	+Info
464392	Polymerization	US\$ 106.32	+Info
464393	Varseo XS Resin Tank	US\$ 170.74	+Info
464394	CAMcreator	US\$ 507.00	+Info
464395	Quick Guide	US\$ 26.00	+Info
464396	Operating Instructions	US\$ 26.00	+Info
464409	Varseo X2 2 Year Warranty Plus	US\$ 1,542.13	+Info
464887	Resin Tank Cover 3pk	US\$ 53.14	+Info

VarseoSmile

VarseoSmile

Code	Name	Price	+
464397	Activation of VarseoSmile Temp	US\$ 295.75	+Info

Varseo Printer unit

Varseo & Otoflash



Code	Name	Price	+
464398	VARSEO XS 100-240 V AC 50/60Hz whit	US\$ 4,300.34	+Info
464886	Otoflash Light Curing Device	US\$ 5,135.13	+Info

3D Printers - RENFERT

3D Printers - RENFERT

Code	Name	Price	+
474855	SIMPLEX Nozzle 0.4	US\$ 124.85	+Info
473956	SIMPLEX 3D-Filament-Printer	US\$ 3,890.25	+Info

3D PRINTERS - RITON 3D



3d Printer Accessories RITON

ACCESSORIES

ACCESSORIES

Code	Name	Price	+
474878	LCD exposure screen	US\$ 310.00	+Info
474879	Release film x1	US\$ 14.00	+Info

Other Models

Other Models

Code	Name	Price	+
474877	RXDent- W90 UV oven	US\$ 920.00	+Info

Other Models

Other Models

Code	Name	Price	+
474865	Ti-150 Laser 3D Printer one set	US\$ 0.00	+Info
474866	D-150 3D Laser Printer one set	US\$ 0.00	+Info
474867	Dual-150 3D Laser Printer one set	US\$ 0.00	+Info
474868	Dual-200 3D Laser Printer one set	US\$ 0.00	+Info
474869	design for dental processing	US\$ 0.00	+Info
474870	Sinter Furnace	US\$ 0.00	+Info
474871	design for dental processing	US\$ 0.00	+Info
474872	design for dental processing	US\$ 0.00	+Info
474875	RXDent- L 230 LCD Resin 3D Printer	US\$ 2,896.00	+Info
474876	RXDent- D190 DLP Resin 3D Printer	US\$ 0.00	+Info

3D Printers - FORMLABS

3D Printers - FORMLABS



ACCESSORIES for FORM 2 & 3 RESIN TANK

FORM 2 & FORM 3 RESIN TANK - FORMLAB

Features:

The Resin Tank blocks light and is made from shatter-resistant polycarbonate. It comes with a reusable lid and wiper, making it easy to stack tanks and store resin outside of the machine.



FORM 3

Code	Name	Price	+
510510	FORMLABS - Form 3 Resin Tank V2.1 #RT-F3-02-01	US\$ 219.03	+Info

FORM 2

Code	Name	Price	+
416586	FORMLABS - FORM-2 Resin Tank 2x #RT-F2-01 - #RTF201	US\$ 70.80	+Info
416587	FORMLABS - FORM-2 Resin Tank LT 20x #RT-F2-02 - #RTF202	US\$ 168.00	+Info

3D Printers - YESCOM

3D Printers - YESCOM

Code	Name	Price	+
420342	LIGHT BOX for 3D Resins 12lt	US\$ 123.04	+Info

CURING UNITS 4 3D Printers



Code	Name	Price	+
474958	ACKURETTA - CURIE PLUS Post processing UV Oven for 3D Printer - 110v + N2 Cover...	US\$ 2,646.00	+Info
474957	ACKURETTA - CURIE Post processing UV Oven for 3D Printer - 110v 1 year Warranty	US\$ 1,416.56	+Info

ACCESSORIES 4 3D Printers

3D PRINTING ACCESSORY EQUIPMENT



SOL 3D Printer ACCESSORIES

Building Platforms

Code	Name	Price	+
474976	ACKURETTA - SOL - Build Platform - Large - # KT210010106	US\$ 236.25	+Info
474977	ACKURETTA - SOL - Build Platform - Medium - # KT210050106	US\$ 236.25	+Info
474978	ACKURETTA - SOL - Build Platform - Small - # KT210040106	US\$ 236.25	+Info

Other Accessories

Other Models

Code	Name	Price	+
474979	Starter Tool Kit	US\$ 47.25	+Info
474980	SOL - Print Screen Pack with 2	US\$ 378.00	+Info
475017	SOL 3D Printer CONCIERGE VIP	US\$ 1,191.45	+Info

WASHING UNITS 4 3D Printing

SOL LCD 3D



ACKURETTA

ACKURETTA

Code	Name	Price	+
474815	SOL 3D Printer Washing Station	US\$ 724.50	+Info

RITON 3D

RITON 3D

Code	Name	Price	+
474882	Dental model washing	US\$ 96.00	+Info

CAD/CAM Burs



CAD/CAM Burs | ATTRITOR

CAD/CAM Burs - ATTRITOR



ROLAND

Roland CAD CAM Burs

DC (Diamond Coated)

DC (Diamond Coate

Code	Name	Price	+
124076	ROLANDCad Bur D round 0.6mm ea	US\$ 32.89	+Info
124077	ROLAND Cad Bur DC Round 1mm	US\$ 32.89	+Info
124078	ROLANDCad Bur D round 2mm ea	US\$ 32.89	+Info

DC for Titanium

Diamond Coated for Titanium

Code	Name	Price	+
------	------	-------	---

474906	ROLAND Cad Bur TI 1mm	US\$ 28.03	+Info
474907	ROLAND Cad Bur TI 2mm	US\$ 28.60	+Info
474908	ROLAND Cad Bur TI 3mm	US\$ 28.03	+Info

DLC



Kits /Set

Other Models

Code	Name	Price	+
420987	ROLAND Cad Bur DLC kit 13pcs	US\$ 228.80	+Info

Roland

Roland

Code	Name	Price	+
115566	ROLAND D50 - D51 Cad Bur 2mm	US\$ 20.24	+Info
115572	ROLAND D50 - D51 Cad Bur 1mm	US\$ 20.24	+Info
115578	ROLAND D50 - D51 Cad Bur 0.8mm	US\$ 25.05	+Info
115583	ROLAND D50 - D51 Cad Bur 0.6mm	US\$ 20.24	+Info

Wet Milling

Diamond Grinding burs for Wet Milling

Code	Name	Price	+
124121	Roland Diam Cad Bur 1mm	US\$ 36.43	+Info
124122	Roland Diam Cad Bur 0.6mm	US\$ 36.43	+Info
124123	Roland Diam Cad Bur 2.5mm	US\$ 36.43	+Info

AMANN GIRRBACH

Amann G.

DC (Diamond Coated)

DC (Diamond Coate

Code	Name	Price	+
124070	AG Cad Bur DC round 0.6mm ea	US\$ 34.32	+Info
124071	AG Cad Bur DC round 1mm ea	US\$ 34.32	+Info
124072	AG Cad Bur DC round 2.5mm ea	US\$ 34.32	+Info

DLC

DLC

Code	Name	Price	+
115565	AMANN GIRRBACH Cad Bur 2.5mm	US\$ 21.12	+Info
115571	AMANN GIRRBACH Cad Bur 1mm	US\$ 21.12	+Info
115577	AMANN GIRRBACH Cad Bur 0.6mm	US\$ 21.12	+Info

Others

Others

Code	Name	Price	+
416754	AMANN GIRRBACH Cad Bur 0.3mm	US\$ 26.40	+Info

Wet Milling

Wet Milling (Disi

Code	Name	Price	+
124124	AG Diam Cad Bur 0.4mm Glass	US\$ 42.24	+Info
124125	AG Diam Cad Bur 1mm Glass	US\$ 42.24	+Info
124126	AG Diam Cad Bur 1.4mm Glass	US\$ 42.24	+Info
124127	AG Diam Cad Bur 1.8mm Glass	US\$ 42.24	+Info

CEREC

Cerec

Compact

Compact

Code	Name	Price	+
113474	Cyl Ptd Bur Compact-6pk	US\$ 149.30	+Info
113475	Step Bur 12 for Compact-6pk	US\$ 139.64	+Info

MCXL

MCXL

Code	Name	Price	+
112250	Cylinder Ptd Bur #12S for MCXL	US\$ 24.89	+Info
112251	Step Bur #12S for MCXL ea	US\$ 24.89	+Info
116130	Cylinder Ptd Bur #20 for MCXL	US\$ 24.89	+Info
116131	Step Bur #20 for MCXL ea	US\$ 24.89	+Info

VHF

VHF

Code	Name	Price	+
124063	VHF Cad Bur DC coated 0.6mm ea	US\$ 28.60	+Info
124064	VHF Cad Bur Dc coated 1mm ea	US\$ 28.60	+Info
124065	VHF Cad Bur DC coated 2mm ea	US\$ 28.60	+Info

DLC

VHF & YETI DLC CAD/CAM COATED MILLING BUR - ATTRITOR

Features:

- * Ultra-precision.
- * Hard, wear-resistant material.
- * Last longer.
- * Can be applied to milling the Zirconia disc, PMMA disc, Acrylic disc for CAD/CAM

Quantity: 1 pc

Sizes: 0.6 mm - 1mm -2mm

Made in Germany



Set

Code	Name	Price	+
416999	VHF Cad/Cam DLC Coated Milling Bur - 40mm long 15 pcs set = 7 x 2mm - 5 1mm & 3 ...	US\$ 747.65	+Info

35mm 4 Axes

Code	Name	Price	+
115582	VHF & YETI Cad/Cam DLC Coated Milling Bur 0.6mm - 35mm long Each - Made in Germa...	US\$ 30.49	+Info
115575	VHF & YETI Cad/Cam DLC Coated Milling Bur 1mm - 35mm long Each - Made in Germany...	US\$ 30.49	+Info
115569	VHF & YETI Cad/Cam DLC Coated Milling Bur 2mm - 35mm long Each - Made in Germany...	US\$ 30.80	+Info

40mm 5 Axes

Code	Name	Price	+
124307	VHF & YETI Cad/Cam DLC Coated Milling Bur 0.6mm - 40mm long Each - Made in Germa...	US\$ 24.64	+Info
124306	VHF & YETI Cad/Cam DLC Coated Milling Bur 1mm - 40mm long Each - Made in Germany...	US\$ 30.18	+Info
124305	VHF & YETI Cad/Cam DLC Coated Milling Bur 2.5mm - 40mm long Each - Made in Germa...	US\$ 29.26	+Info

WET Milling



Code	Name	Price	+
417087	ATTRITOR - VHF & YETI Disilicate Sinter Diamond CadCam Bur 2.6mm FLAT END for di...	US\$ 34.50	+Info
417088	ATTRITOR - VHF & YETI Disilicate Sinter Diamond CadCam Bur 2.4mm BALL END for di...	US\$ 30.27	+Info
417081	ATTRITOR - VHF & YETI Disilicate Sinter Diamond CadCam Bur 1mm BALL END for disi...	US\$ 35.20	+Info
417086	ATTRITOR - VHF & YETI Disilicate Sinter Diamond CadCam Bur 0.6mm BALL END for di...	US\$ 35.20	+Info

WIELAND
WIELAND

DC (Diamond Coate

DC (Diamond Coate

Code	Name	Price	+
124067	WIEL Cad Bur D coated 0.7mm ea	US\$ 68.03	+Info
124068	WIEL Cad Bur D coated 1mm ea	US\$ 68.03	+Info
124069	WIEL Cad Bur D coated 2-5mm ea	US\$ 68.03	+Info

DLC

DLC

Code	Name	Price	+
115568	WLND Cad Bur 2.0mm x 40mm	US\$ 15.66	+Info
115574	WLND Cad Bur 1mm x 40mm	US\$ 16.54	+Info
115580	WLND Cad Bur 0.6mm x 40mm	US\$ 10.91	+Info

Wieland Mini

Wieland Mini

DLC

DLC

Code	Name	Price	+
115573	WLND MINI Cad Bur 1mm x 35mm	US\$ 17.42	+Info
115579	WLND MINI Cad Bur 0.6mm x 35mm	US\$ 17.42	+Info
420161	WLND MINI Cad Bur 0.6mm x 40mm	US\$ 17.60	+Info
420162	WLND MINI Cad Bur 1mm x 40mm	US\$ 17.07	+Info
420163	WLND MINI Cad Bur 2.5mm x 40mm	US\$ 17.60	+Info
115567	WLND MINI Cad Bur 2.5mm x 35mm	US\$ 17.60	+Info

ZIRKONZAHN

ZirkonZahn CAD CAM Burs

DC (Diamond Coated)

DC (Diamond Coate

Code	Name	Price	+
124073	ZIRK Cad Bur D coated 0.6mm ea	US\$ 19.36	+Info
124074	ZIRK Cad Bur D coated 1mm ea	US\$ 35.20	+Info
124075	ZIRK Cad Bur D coated 2mm ea	US\$ 20.24	+Info

DLC

DLC

Code	Name	Price	+
115570	ZirkonZahn M5 Cad Bur 2mm	US\$ 30.80	+Info
115576	ZirkonZahn M5 Cad Bur 1mm	US\$ 30.80	+Info
115581	ZirkonZahn M5 Cad Bur 0.5mm	US\$ 30.80	+Info

CAD/CAM Burs | ROLAND (original)



CARBIDE (DLC)

BALL END

Code	Name	Price	+
420215	ROLAND 0.3mm Carbide Ball End Mill - High Quality ZCB-15D-US	US\$ 75.00	+Info
420214	ROLAND 0.6mm Carbide Ball End Mill ZCB-30D-US	US\$ 75.00	+Info
420213	ROLAND 1mm Carbide Ball End Mill - High Quality ZCB-50D-US	US\$ 75.00	+Info
420212	ROLAND 2mm Carbide Ball End Mill - High Quality ZCB-100D-US	US\$ 75.00	+Info

FLAT END

Code	Name	Price	+
420300	ROLAND 0.6mm Carbide Flat End Mill #PT06-10 Extremo Chato	US\$ 38.51	+Info
420301	ROLAND 1mm Carbide Flat End Mill #PT1- 20 Extremo Chato	US\$ 38.51	+Info
420302	ROLAND 2mm Carbide Flat End Mill #PT2- 20 Extremo Chato	US\$ 38.51	+Info
420306	ROLAND 3mm Flat End Carbide Milling Bur #ZRB-100D Extremo Chato	US\$ 116.95	+Info

TORIC

Code	Name	Price	+
420334	ROLAND Rapid Razor-Sharp 3mm Carbide Toric End Bur SKU#RSB-150D-US Made in Germany	US\$ 121.71	+Info

DIAMOND COATED

Code	Name	Price	+
420218	ROLAND 0.6mm DIAMOND COATED DC Ball End Mill #ZDB-30D-US	US\$ 144.46	+Info
420217	ROLAND 1mm DIAMOND COATED DC Ball End Mill ZDB-50D-US	US\$ 144.46	+Info
420216	ROLAND 2mm DIAMOND COATED DC Ball End Mill ZDB-100D-US	US\$ 144.46	+Info

DIAMOND GRINDING BURS (WET)

For Glass Ceramics

LONG DIAMOND GRINDING BURS

Code	Name	Price	+
420316	ROLAND 0.5mm Long Diamond Grinding Bur #ZGB2- 25D ZBG2- 25D Wet Milling Cad Burs...	US\$ 94.81	+Info

420315	ROLAND 1mm Long Diamond Grinding Bur ZGB2-50D ZBG2-50D Wet Milling Cad Burs for ...	US\$ 94.81	+Info
420314	ROLAND 2.5mm Long Diamond Grinding Bur #ZGB2-125D ZBG2-125D Wet Milling Cad Burs...	US\$ 94.81	+Info
474834	ROLAND 1.5mm LongDiamond Grinding Bur #ZGB2-75D Wet Milling Cad Burs for Disilic...	US\$ 94.81	+Info

REGULAR DIAMOND GRINDING BURS

Code	Name	Price	+
420312	ROLAND 0.5mm Regular Diamond Grinding Bur #ZGB- 25D Wet Milling Cad Burs for Dis...	US\$ 67.78	+Info
420313	ROLAND 1mm Regular Diamond Grinding Bur #ZGB-50D Wet Milling Cad Burs for Disili...	US\$ 68.46	+Info
420311	ROLAND 2.5mm Regular Diamond Grinding Bur #ZGB-125D Wet Milling Cad Burs for Dis...	US\$ 68.46	+Info

SHORT DIAMOND GRINDING BURS

Code	Name	Price	+
420318	ROLAND 2.5mm Short Diamond Grinding Bur #ZGB2-125DS ZBG2-125DS Wet Milling Cad B...	US\$ 94.81	+Info

Diamond-Like

Code	Name	Price	+
475028	ROLAND 2.0mm Diamond-Like-Coated Milling Tool for PMMA Item # ZPB-100D	US\$ 83.51	+Info
475029	ROLAND 1.0mm Diamond-Like-Coated Milling Tool for PMMA Item # ZPB-50D	US\$ 83.51	+Info
475030	ROLAND 0.6mm Diamond-Like-Coated Milling Tool for PMMA Item # ZPB-30D	US\$ 83.51	+Info

LOLLIPOP

Code	Name	Price	+
420333	ROLAND 2mm Carbide Lollipop End Mill SKU#LOLLIPOP Made in Germany	US\$ 34.70	+Info

SPECIAL COATED

Code	Name	Price	+
420221	ROLAND 0.6mm Special Coated Ball End Mill-Ccomposites Hybrid #ZDB-30D-HY	US\$ 169.96	+Info
420220	ROLAND 1mm Special Coated Ball End Mill-Ccomposites Hybrid #ZDB-50D-HY	US\$ 169.96	+Info
420219	ROLAND 2mm Special Coated Ball End Mill-Ccomposites Hybrid #ZDB-100D-HY	US\$ 169.96	+Info

CAD/CAM Burs | UP3D

CAD/CAM Burs - UP3D



DC (Diam Coated)

DC (Diam Coated)



Code	Name	Price	+
474321	DC Diamond Coated Cad Bur 2mm	US\$ 30.71	+Info
474322	DC P52/53 Cad Bur 1mm	US\$ 30.71	+Info
474323	UP3D - P52/53 Cad Bur 0.6mm	US\$ 30.71	+Info
474324	UP3D - P52/53 Cad Bur 0.3mm	US\$ 30.71	+Info

DLC

DLC



Code	Name	Price	+
474339	UP3D - P52/53 Cad Bur 0.3mm	US\$ 20.73	+Info
474342	UP3D - P52/53 Cad Bur 0.6mm	US\$ 20.73	+Info
474343	UP3D - P52/53 Cad Bur 1.0mm	US\$ 20.73	+Info
474344	UP3D - P52/53 Cad Bur 2.0mm	US\$ 20.73	+Info

DLC (PMMA)

DLC (PMMA)

Code	Name	Price	+
474325	UP3D - PMMA Bur 2.0mm	US\$ 26.18	+Info
474326	UP3D - PMMA Bur 1.0mm	US\$ 26.18	+Info
474327	UP3D - PMMA Bur 0.6mm	US\$ 26.18	+Info

CAD/CAM Burs | AMANN GIRRBACH

CAD/CAM Burs - AMANN GIRRBACH



Other Models

Other Models

Code	Name	Price	+
420319	F3 Ceramill Roto Bur	US\$ 100.46	+Info

RFID

RFID

Code	Name	Price	+
468264	Roto RFID 1,2 Drill	US\$ 66.26	+Info
468265	Roto RFID 1,0 Wax Denture	US\$ 72.39	+Info
468266	Roto RFID 3,0 Wax Denture	US\$ 72.39	+Info

468267	Roto RFID 2,5 Model	US\$ 75.17	+Info
468268	Roto RFID 2,0 Model	US\$ 75.17	+Info

RFID CoCr

RFID CoCr

Code	Name	Price	+
468259	Roto RFID 2,5 CoCr	US\$ 72.39	+Info
468260	Roto RFID 1,5 CoCr	US\$ 72.39	+Info
468261	Roto RFID 1,0 CoCr	US\$ 72.39	+Info

RFID DMB DC

RFID DMB DC

Code	Name	Price	+
468255	Roto RFID 2,5 DMB DC	US\$ 219.16	+Info
468256	Roto RFID 1,0 DMB DC	US\$ 219.16	+Info
468257	Roto RFID 0,6 DMB DC	US\$ 219.16	+Info
468258	Roto RFID 0,3 DMB DC	US\$ 219.16	+Info

RFID PMMA

RFID PMMA

Code	Name	Price	+
468243	Roto RFID 0,6 PMMA/Wax	US\$ 92.03	+Info
468244	Roto RFID 0,3 PMMA/Wax	US\$ 92.03	+Info
468262	Roto RFID 2,5 PMMA Denture	US\$ 203.86	+Info
468263	Roto RFID 1,5 PMMA Denture	US\$ 203.86	+Info

RFID SC

RFID SC

Code	Name	Price	+
468251	Roto RFID 2,5 SC	US\$ 107.97	+Info
468252	Roto RFID 1,0 SC	US\$ 107.97	+Info
468253	Roto RFID 0,6 SC	US\$ 107.97	+Info
468254	Roto RFID 0,3 SC	US\$ 112.88	+Info

RFID Sintron

RFID Sintron

Code	Name	Price	+
468239	Roto RFID 2,5 Sintron	US\$ 93.01	+Info
468240	Roto RFID 1,0 Sintron	US\$ 93.01	+Info
468241	Roto RFID 0,6 Sintron	US\$ 93.01	+Info
468242	Roto RFID 0,3 Sintron	US\$ 93.01	+Info

RFID Titanio

RFID Titanio

Code	Name	Price	+
467774	Calibration pin Roto RFID	US\$ 137.42	+Info
468245	Roto RFID 2,0 Ti	US\$ 72.39	+Info
468246	Roto RFID 1,0 Ti	US\$ 72.39	+Info

469827	Roto RFID Calibration Pin	US\$ 142.70	+Info
--------	---------------------------	-------------	-----------------------

RFID Zirconia

RFID Zirconia

Code	Name	Price	+
468235	Roto RFID 2,5 ZI	US\$ 89.56	+Info
468236	Roto RFID 1,0 ZI	US\$ 89.56	+Info
468237	Roto RFID 0,6 ZI	US\$ 89.56	+Info
468238	Roto RFID 0,3 ZI	US\$ 89.56	+Info

ROTO MOTION

ROTO MOTION

Code	Name	Price	+
114391	Ceramill Roto Motion 1.0mm 1pk	US\$ 95.56	+Info
114392	Ceramill Roto Motion 2.5mm 1pk	US\$ 95.56	+Info

WET Roto Diamond

WET Roto Diamond

Code	Name	Price	+
420951	Roto Diamond Burs 1.4mm	US\$ 59.90	+Info
420952	Roto Diamond Burs 1.0mm	US\$ 59.90	+Info
464820	Roto Diamond Burs 1.8mm	US\$ 63.71	+Info

Zirconia KF

Zirconia KF

Code	Name	Price	+
760792	Zirconia Bur KF1.8mm #18	US\$ 100.46	+Info
760892	Zirconia Reader KT1.8mm #17	US\$ 100.46	+Info

Zirconia T

Zirconia T

Code	Name	Price	+
760709	Zirconia Bur F0.9mm #10	US\$ 100.46	+Info
760721	Zirconia Bur F2	US\$ 100.46	+Info
760809	Zirconia Reader T0.9mm #9	US\$ 100.46	+Info
760811	Zirconia Reader T1.2	US\$ 100.46	+Info
760821	Zirconia Reader T2	US\$ 74.34	+Info
760831	Zirconia Reader T3	US\$ 100.46	+Info

HANSONG CAD CAM BURS



Amann Girrbach Burs

Code	Name	Price	+
436910	Hansong Milling Burs Amann Girrbach	US\$ 34.85	+Info
436913	Hansong Milling Burs Amann Girrbach	US\$ 36.09	+Info
436926	Hansong Milling Burs Amann Girrbach	US\$ 39.88	+Info

Roland Burs

Code	Name	Price	+
436911	Hansong Milling Burs Roland Type-2	US\$ 35.85	+Info
436916	Hansong Milling Burs Roland Type- 2	US\$ 37.88	+Info
436919	Hansong Milling Burs Roland Type- 2	US\$ 37.88	+Info
436920	Hansong Milling Burs Roland Type- 2	US\$ 37.88	+Info
436922	Hansong Milling Burs Roland Type- 2	US\$ 37.88	+Info
436931	Hansong Milling Burs Roland Type-2O	US\$ 34.71	+Info
436939	Hansong Milling Burs Roland Type- 2	US\$ 37.88	+Info
436941	Hansong Milling Burs Roland Type-2	US\$ 41.58	+Info

Sirona Type Burs

Code	Name	Price	+
436927	Hansong Milling Burs Sirona Type- 2	US\$ 37.88	+Info
436930	Hansong Milling Burs Sirona Type- 2	US\$ 42.74	+Info
436934	Hansong Milling Burs Sirona Type- 2	US\$ 40.43	+Info
436936	Hansong Milling Burs Sirona Type- 2	US\$ 42.74	+Info
436937	Hansong Milling Burs Sirona Type- 2	US\$ 42.74	+Info
436943	Hansong Milling Burs Sirona Type- 2	US\$ 38.11	+Info

VHF Burs

Code	Name	Price	+
436912	Hansong Milling Burs VHF Type-2 VTR	US\$ 37.88	+Info

436914	Hansong Milling Burs VHF Type-2 VTR	US\$ 37.88	+Info
436923	Hansong Milling Burs VHF Type-2 VDR	US\$ 38.11	+Info
436924	Hansong Milling Burs VHF Type-2 VDR	US\$ 38.11	+Info
436929	Hansong Milling Burs VHF Type-2 VDR	US\$ 40.43	+Info
436932	Hansong Milling Burs VHF Type-2 VDR	US\$ 42.74	+Info
436933	Hansong Milling Burs VHF Type-2 VDR	US\$ 34.65	+Info
436935	Hansong Milling Burs VHF Type-2 VDR	US\$ 41.58	+Info
436938	Hansong Milling Burs VHF Type-2 VDR	US\$ 38.11	+Info
436944	Hansong Milling Burs VHF Type-2 VDR	US\$ 42.74	+Info

Wieland Type Burs


Code	Name	Price	+
436915	Hansong Milling Burs Wieland Type-	US\$ 41.58	+Info
436918	Hansong Milling Burs Wieland Type-	US\$ 42.74	+Info
436928	Hansong Milling Burs Wieland Type-	US\$ 38.11	+Info
436940	Hansong Milling Burs Wieland Type-	US\$ 41.58	+Info
436942	Hansong Milling Burs Wieland Type-	US\$ 45.05	+Info

ZirkonZahn Burs

Code	Name	Price	+
436917	Hansong Milling Burs Zirkonzahn Typ	US\$ 45.05	+Info
436921	Hansong Milling Burs Zirkonzahn Typ	US\$ 50.82	+Info
436925	Hansong Milling Burs Zirkonzahn Typ	US\$ 45.05	+Info

ZIRCONIA & EMAX Abrasives

ZIRCONIA & EMAX ABRASIVES - ATT RITOR

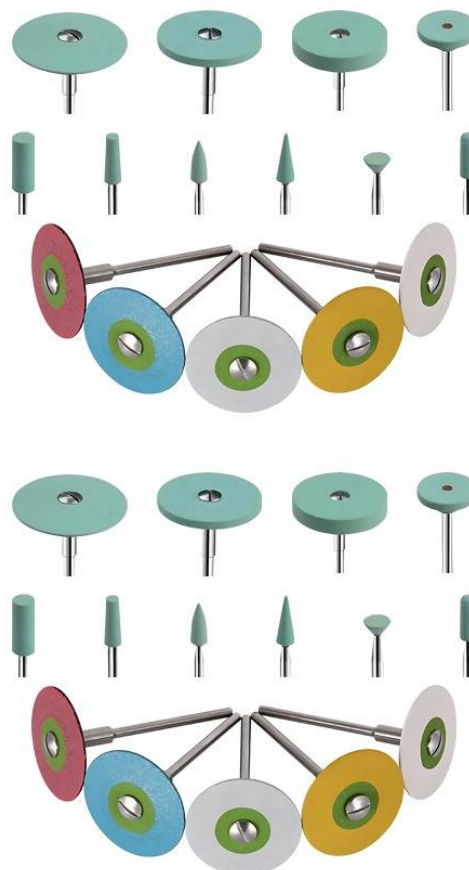

Zirconia & Emax GRINDERS

Zirconia & Emax Grinders - ATTRITOR - Hard
Stones made of SINTERED DIAMOND in ALUMINA
Matrix



Kits / Sets

Kits / Sets

Code	Name	Price	+
416991	ZIRCONIA GRINDER KIT 8pcs	US\$ 187.02	+Info

Z-GRINDER

Z-Grinder

Code	Name	Price	+
112947	Z-Grindr Zirconia Wheel 15x3.5	US\$ 22.00	+Info
112948	Z-Grindr Zirconia Taper 3.5x11	US\$ 13.75	+Info
112228	Z-Grindr Zirconia Wheel 22x2	US\$ 22.00	+Info
100754	Z-Grindr Zirconia Wheel 22x4.5	US\$ 22.00	+Info
100755	Z-Grindr Zirconia Med Cyl 5x13	US\$ 13.75	+Info
100756	Z-Grindr Zirconia Lg Barr 6x13	US\$ 16.50	+Info

Z-ZIRCON

Z-Zircon

Code	Name	Price	+
115225	Zirconia Z-Zircon Med Inv Cone	US\$ 24.75	+Info
114356	Zirconia Z-Zircon Disc 22x0.7	US\$ 22.00	+Info
114175	Zirconia Z-ZirconWheel 15x3.5	US\$ 22.00	+Info
114176	Zirconia Z-Zircon Wheel 22x2	US\$ 22.00	+Info
114177	Zirconia Z-Zircon Wheel 22x4.5	US\$ 26.95	+Info

Zirconia & Emax Polishers

Zirconia & Emax Polishers - ATTRITOR



DURA-DISCS Silicone

Dura-Discs

Code	Name	Price	+
112956	High Shine Bicolor Yell S-Fine	US\$ 16.90	+Info
112957	High Shine Bicolor Grey Fine	US\$ 19.01	+Info
112958	High Shine Bicolor Blue Medium	US\$ 16.90	+Info
112959	High Shine Bicolor red Coarse	US\$ 19.01	+Info

DURA-POINTS Silicone

Dura-Point

Code	Name	Price	+
100932	ATTR Dura Point - Fine - Pink	US\$ 15.84	+Info
100933	ATTR Dura Point- X Fine Blue	US\$ 19.80	+Info

DURA-WHEELS Silicone

Dura-Wheel

Code	Name	Price	+
100930	ATTR Dura Wheel- X-Coarse-Br	US\$ 16.90	+Info
100722	Dura Wheel - Medium-Blue	US\$ 16.90	+Info
100723	ATTR Dura Wheel -Coarse-Gray	US\$ 16.90	+Info
100763	ATTR Dura Wheel - Fine - Pink	US\$ 16.90	+Info
115563	ATTR Dura Wheel X Fine -White	US\$ 16.90	+Info

KITS / Sets

Kits / Sets

Code	Name	Price	+
416558	ZIRCONIA POLISHER KIT 6pcs	US\$ 133.37	+Info

RUBBER DISCS

Rubber Discs

Code	Name	Price	+
112949	Zirconia Polish Gray Sup Corse	US\$ 19.01	+Info
112950	Zirconia Polish Green Coarse	US\$ 16.90	+Info
112951	Zirconia Polish Blue Medium	US\$ 16.90	+Info
112952	Zirconia Polish red Fine	US\$ 16.90	+Info
112953	Zirconia Polish Lt Gray S-Fine	US\$ 19.01	+Info

ZIRCONIA & EMAX ABRASIVES - SHOFU



CERAMISTE HP | SHOFU

SHOFU - HP Ceramiste Polisher Standard - WH6 - 12pk
Porcelain Polishers

Advantage!

Ceramiste is specifically made for polishing **porcelains and enamel**.

Ceramiste Soft is designed for **contouring and polishing** **porcelain veneers**, inlays and on lays in the laboratory.

- * **Aluminum Oxide impregnated Silicon polishers**
- * Three-step system
- * Use Ceramiste Standard for pre-polishing. Ultra for polishing and Ultra II for super-polishing
- * Use Ceramiste Soft PA and Ceramiste Soft PB for smoothing rough surfaces and fragile porcelain margins
- * Use Ceramiste for: porcelains and enamel



Coral Stones (Shaping)

Code	Name	Price	+
420360	SHOFU Hp Coral Stones #1-12pk - #481A	US\$ 21.73	+Info
420361	SHOFU Hp Coral Stones #2-12pk - #482A	US\$ 21.73	+Info
420362	SHOFU Hp Coral Stones #3-12pk - #483A	US\$ 21.73	+Info
420363	SHOFU Hp Coral Stones #4-12pk - #484A	US\$ 21.08	+Info
420364	SHOFU Hp Coral Stones #5-12pk - #485A	US\$ 21.08	+Info
420365	SHOFU Hp Coral Stones #6-12pk - #486A	US\$ 21.73	+Info

Standard (Pre Polishing)

Code	Name	Price	+
106777	SHOFU HP Ceramiste Polisher Standard - KN7-12pk - #0261	US\$ 26.78	+Info
106790	SHOFU HP Ceramiste Polisher Standard - PC2-12pk - #0262	US\$ 26.78	+Info
106791	SHOFU HP Ceramiste Polisher Standard - WH6-12pk - #0263	US\$ 26.78	+Info

Ultra Polish (Polish A)

Code	Name	Price	+
106792	SHOFU HP Ceramiste Polisher UL - KN7-12pk - #0264	US\$ 46.00	+Info
106794	SHOFU HP Ceramiste Polisher UL - PC2-12pk - #0265	US\$ 23.62	+Info
106795	SHOFU HP Ceramiste Polisher UL - WH6-12pk - #0266	US\$ 26.51	+Info

Ultra Polish II (Polish B)

Code	Name	Price	+
106796	SHOFU HP Ceramiste Polisher UL2 - KN7-12pk - #0267	US\$ 36.80	+Info

106797	SHOFU HP Ceramiste Polisher UL2 - PC2-12pk - #0268	US\$ 26.78	+Info
106799	SHOFU HP Ceramiste Polisher UL2 - WH6-12pk - #0269	US\$ 26.78	+Info

Ceramiste Soft (Marginal Borders)

Code	Name	Price	+
111726	SHOFU HP Ceramiste Soft Pb Point PC2-12pk - #0563	US\$ 31.32	+Info
111730	SHOFU HP Ceramiste Soft Pa Wheel WH6-12pk - #0564	US\$ 32.99	+Info
111732	SHOFU HP Ceramiste Soft Cut Wheel WH6-12pk - #0565	US\$ 32.65	+Info

CERAMASTER HP POLISHERS

CERAMASTER HP POLISHERS - SHOFU

CeraMaster is designed to finish, **polish and super-polish porcelains and enamel**.

A carefully balanced blend of diamond particles provides polishers & fnishers that produce the smoothest finish and most glaze-like surfaces. Excellent for setting, retouching or adjusting, Porcelain, Zirconia, Full controur zirconia restorations, inlay-onlay, Full ceranmics like e.Max, or Enamel either Chairside or intraoral.



From Coarse to Super finishing, a 2 step finishing Diamond based system that makes easy set up restorations

Advantage

- * **Diamond impregnated silicone polishers** special for chairside
- * Saves time by providing a simple, two-step system
- * Use CeraMaster Coarse to polish to a satin finish and CeraMaster to super-polish
- * Use CeraMaster for: porcelains and enamel
- * Faster and simpler than CERAMISTE (aluminum oxide)

Kit

SHOFU - Hp Ceramaster Assorted Silicon Composite Polishers

111516 - Mfg # 130

Advantaje!

Diamond impregnated polishers with a unique blend of silicone and even distribution of diamond particles make these instruments the best choice for polishing porcelain and ceramic polymer to a glazed-like surface.

Features:

A simple 2 step finishing and polishing system, ideal for intra-oral, chairside or laboratory use. CeraMaster Coarse (blue polisher) - for finishing / polishing and adjusting margins of alumina copings
CeraMaster (blue and white bands) - for superpolishing

Durable, cost-effective and available in HP and CA shanks to suit all polishing surfaces.



Available:

PA grit (green band) - for contouring

PB grit (blue band) - for polishing

Contains: Bullet - Mini Point- Kn7 - Wh6

Code	Name	Price	+
111516	SHOFU HP Ceramaster ASSORTED - #0130	US\$ 50.74	+Info
111976	SHOFU Ceramaster Coarse ASSORTED Hp - #130C	US\$ 52.16	+Info

Ceramaster Coarse



Code	Name	Price	+
111966	SHOFU Ceramaster Coarse Hp Mini 3pk 3pk - #129C	US\$ 46.00	+Info
111963	SHOFU Ceramaster Coarse Hp Wheel WH6 3pk - #128C	US\$ 43.24	+Info
111962	SHOFU Ceramaster Coarse Hp Kn7 - #127C	US\$ 46.00	+Info
111961	SHOFU Ceramaster Coarse Hp Bullet 3pk - #126C	US\$ 46.00	+Info
111978	SHOFU Ceramaster Coarse Unmounted 1 Per Pack - #132C	US\$ 44.62	+Info

Ceramaster Regular



Code	Name	Price	+
111512	SHOFU HP Ceramaster Bullet 3pk - #0126	US\$ 46.00	+Info
111514	SHOFU HP Ceramaster Wheel WH6 3pk - #0128	US\$ 43.24	+Info
111513	SHOFU HP Ceramaster Knife Kn7 3pk - #0127	US\$ 46.00	+Info
111515	SHOFU HP Ceramaster MINI 3pk - #0129	US\$ 46.00	+Info

COREMASTER

SHOFU - Coremaster
Zirconia Copins
Features:

CoreMaster is designed for adjusting marginal areas and the thickness of aluminum oxide and **zirconia copings**.

- * Diamond impregnated polishers eliminate microfractures
- * Two grit sizes and three shapes
- * Contouring and finishing
- * No water cooling needed

Use CoreMaster for: zirconia frameworks, alumina frameworks and porcelains


Code

Name

Price

+

106772	Coremaster WH6 HP - 2pk - #0140C	US\$ 42.32	+Info
106767	Coremaster KN7 HP 2pk	US\$ 36.80	+Info
106773	Coremaster Stnd UM Wheel HP 1p	US\$ 45.25	+Info

ZIRCONIA & EMAX ABRASIVES - MPF



Zirconia & Emax Grinders

Zirconia & Emax Grinders - MPF BRUSH

WHITE PEARL GRINDERS

WHITE PEARL

Code	Name	Price	+
123202	WHITEPEARL Large Disc 22*4.5	US\$ 41.22	+Info
123203	WHITEPEARL Small Disc 15*3.5	US\$ 36.78	+Info
123204	WHITEPEARL Disc 22*2	US\$ 60.52	+Info
123205	WHITEPEARL Thin Disc 22*0.7	US\$ 35.03	+Info
123206	WHITEPEARL Small Barrel 4.8*13	US\$ 27.80	+Info
123207	WHITEPEARL Large Barrel 6.5*13	US\$ 31.59	+Info
123208	WHITEPEARL Cone 3.5*11	US\$ 33.17	+Info
123209	WHITEPEARL Invertd Cone 12*6	US\$ 20.10	+Info
123210	WHITEPEARL Invertd S Cone 6*8	US\$ 28.71	+Info

ZIRCO CERA GRINDERS

ZIRCO CERA

Code	Name	Price	+
123192	ZIRCO CERA Large Disc 22*4.5	US\$ 38.93	+Info
123194	ZIRCO CERA Small Disc 15*3.5	US\$ 32.04	+Info
123195	ZIRCO CERA Disc 22*2	US\$ 50.90	+Info
123196	ZIRCO CERA Thin Disc 22*0.7	US\$ 49.95	+Info
123197	ZIRCO CERA Small Barrel 4.8*13	US\$ 28.60	+Info
123198	ZIRCO CERA Large Barrel 6.5*13	US\$ 28.60	+Info
123199	ZIRCO CERA Cone 3.5*11	US\$ 30.03	+Info
123200	ZIRCO CERA Inverted Cone 12*6	US\$ 26.98	+Info
123201	ZIRCO CERA Inverted S Cone 6*8	US\$ 26.98	+Info

ZMAX GRINDERS

ZMAX

Code	Name	Price	+
123183	ZMAX Large Disc 22*4.5	US\$ 26.10	+Info
123184	ZMAX Small Disc 15*3.5	US\$ 32.04	+Info
123185	ZMAX Disc 22*2	US\$ 31.72	+Info
123186	ZMAX Thin Disc 22*0.7	US\$ 31.72	+Info
123187	ZMAX Small Barrel 4.8*13	US\$ 28.60	+Info
123188	ZMAX Large Barrel 6.5*13	US\$ 17.45	+Info
123189	ZMAX Cone 3.5*11	US\$ 28.60	+Info
123190	ZMAX Inverted Cone 12*6	US\$ 26.98	+Info
123191	ZMAX Inverted S Cone 6*8	US\$ 26.98	+Info

Zirconia & Emax Polishers

Zirconia & Emax Polishers - MPF BRUSH



Innovation Silicones

Innovation

Code	Name	Price	+
470835	Zirc / Disil Polish Adjust 5 pcs	US\$ 63.00	+Info
470836	Zirc / Disil Polish Adjust 5 pcs	US\$ 60.48	+Info
470837	Zirc / Disil Polish Adjust 5 pcs	US\$ 44.10	+Info
470838	Zirc / Disil Polish Adjust 5 pcs	US\$ 60.48	+Info
470839	Zirc / Disil Polish Adjust 5 pcs	US\$ 60.48	+Info


Rainbow Silicones

Rainbow

Code	Name	Price	+
123929	Rainbow Coarse red	US\$ 30.89	+Info
123930	Rainbow Medium Blue	US\$ 30.89	+Info
123931	Rainbow Fine Gray	US\$ 24.71	+Info
123932	Rainbow High shine Yellow	US\$ 29.96	+Info

Zirconia & Emax Polishers - UMBRELLA A CONCEPT

Zirconia & Emax Polishers - UMBRELLA
CONCEPT

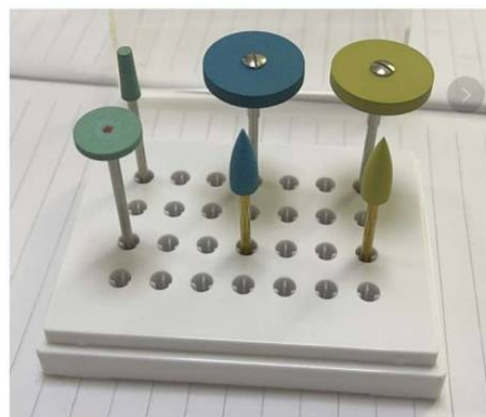

Code	Name	Price	+
470735	NANO Ceramic POLISHER KIT 5pcs	US\$ 24.50	+Info
470762	Zirconia Luster Polishers 6pcs	US\$ 78.00	+Info

Zirconia & Emax Abrasives - Zir.Care SMILE LINE



Code	Name	Price	+
473534	Zir.Care Assort Set (3 pcs)	US\$ 341.04	+Info
473535	Zir.Care CLEANER (1 pc)	US\$ 58.80	+Info
473536	Zir Care Soft, Set of 3 sineter dia	US\$ 341.04	+Info
473537	Zir Care Soft, small lense, 12x2mm	US\$ 108.78	+Info
473538	Zir Care Soft, lense, 20x2mm	US\$ 129.36	+Info
473539	Zir Care Soft, small lense, 20x4mm	US\$ 130.83	+Info

Zirconia & Emax Grinders - BESQUAL

Zirconia & Emax Grinders - BESQUAL



Zirconia Contouring Kit
Z-Kit | #221-000

Z-K Pre-Sintered

Z-K Pre-Sintered

Code	Name	Price	+
469842	Zirconia Contouring Kit	US\$ 39.00	+Info
469843	Pre-Sint Zirc MD Cone 5.6x13mm	US\$ 31.20	+Info
469844	Zirconia Medium Disc 15x0.3mm	US\$ 5.44	+Info
469845	Dual Grit Fine/Medium Egg 7x19mm	US\$ 13.20	+Info
469846	Dual Grit Fine/Medium Egg 10x23mm	US\$ 13.20	+Info
469847	Medium Flame 7x13mm	US\$ 5.44	+Info


Zirconia Diamond

Zirconia Diamond

Code	Name	Price	+
469856	Zircon Diamond Stone KIT 7pk	US\$ 57.33	+Info
469857	Zircon Diam Stone Concave #9 5pk	US\$ 41.78	+Info
469858	Zircon Diam Stone Knife-Edge #10	US\$ 42.63	+Info
469859	Zircon Diam Stone Wheel #11 5pk	US\$ 42.63	+Info
469871	Zircon Diam Stone Taper #20 5pk	US\$ 42.63	+Info
469872	Zircon Diam Stone Barrel #39 5pk	US\$ 36.24	+Info
469873	Zircon Diam Stone Inv Cone #42 5pk	US\$ 42.63	+Info
469874	Zircon Diam Stone Bullet #13	US\$ 42.63	+Info

Zirconia & Emax Polishers - AMANN GIRRBACH

Zirconia & Emax Polishers - AMANN GIRRBACH

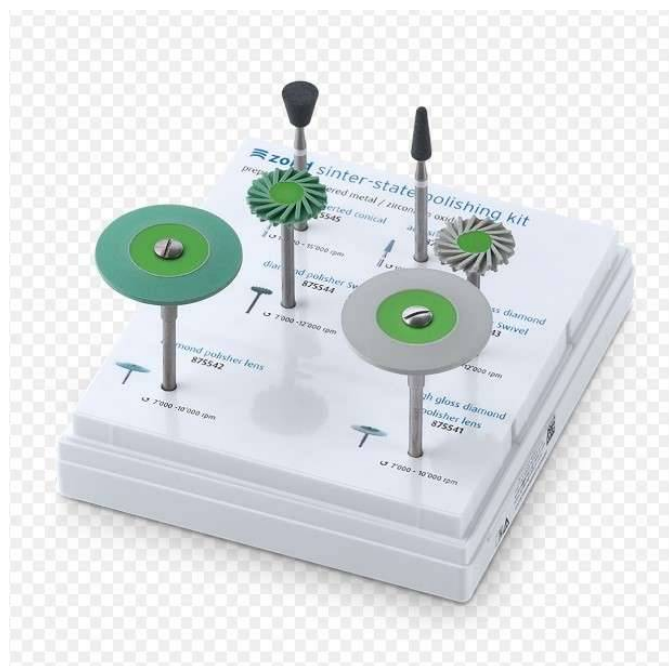

Code	Name	Price	+
468122	Ceramill Zolid FX Starter Kit	US\$ 3,055.64	+Info
468867	universal polisher, 5 pcs.	US\$ 22.77	+Info
468871	polishing disc	US\$ 42.86	+Info
468872	polishing shell	US\$ 40.18	+Info
468873	Zolid preparation dent kit	US\$ 96.44	+Info
468881	high gloss diamond polisher	US\$ 61.61	+Info
468882	diamond polisher	US\$ 61.61	+Info
468883	high gloss diamond polisher	US\$ 41.52	+Info
468884	diamond polisher	US\$ 41.52	+Info
468887	smoothing and pre-polishing	US\$ 41.52	+Info
468888	high gloss polishing	US\$ 41.52	+Info
468890	gloss polishing	US\$ 25.45	+Info
468891	high gloss polishing	US\$ 25.45	+Info
468865	Zolid green-state finishing Kit	US\$ 175.47	+Info
468880	Zolid sinter-state polishing kit	US\$ 207.60	+Info
468893	Zolid polishing dent kit update	US\$ 120.55	+Info

CAD/CAM

3D Accessories

3D Accessories - FORMLABS

3D Accessories - FORMLABS

Code	Name	Price	+
417366	FORMLABS BUILD PLATFORM	US\$ 142.93	+Info

3D Accessories - HARVEST DENTAL

3D Accessories - HARVEST DENTAL

Code	Name	Price	+
475555	InovaPrint Wash 3D Print Resin CI	US\$ 41.47	+Info
475556	InovaPrint Wash 3D Resin Cleaner	US\$ 134.95	+Info

CAD/CAM

3D RESINS

3D NATURALIZATION

3D NATURALIZATION

Code	Name	Price	+
476642	Rodin Palette Naturalization Kit, 1	US\$ 282.69	+Info
476643	Rodin Palette Refill, Blue Light Cu	US\$ 55.89	+Info
476644	Rodin Palette Refill, Yellow Light	US\$ 55.89	+Info
476645	Rodin Palette Refill, Red Light Cur	US\$ 55.89	+Info
476646	Rodin Palette Refill, White Light C	US\$ 55.89	+Info
476647	Rodin Palette Refill, Cream Light C	US\$ 55.89	+Info
476648	Rodin Palette Refill, Orange Light	US\$ 55.89	+Info
476649	Rodin Palette Refill, Violet Light	US\$ 55.89	+Info
476650	Rodin Palette Refill, Dark Brown Li	US\$ 55.89	+Info
476651	Rodin Palette Refill, A-Base Light	US\$ 55.89	+Info
476652	Rodin Palette Refill, B-Base Light	US\$ 55.89	+Info
476653	Rodin Palette Refill, C-Base Light	US\$ 55.89	+Info
476654	Rodin Palette Refill, D-Base Light	US\$ 55.89	+Info
476655	Rodin Palette Refill, Pink Light Cu	US\$ 55.89	+Info
476656	Rodin Palette Refill, Magenta Light	US\$ 55.89	+Info
476657	Rodin Palette, Mixing Well	US\$ 6.74	+Info

C&B MFH Hybrid

C&B MFH Hybrid



RODIN by PAC-DENT

RODIN by PAC-DENT

SCULPTURE

SCULPTURE

1.2 kg

1.2 kg

Code	Name	Price	+
476583	Rodin Sculpture, Shade B4, 1.2kg/Bo	US\$ 647.19	+Info
476591	Rodin Sculpture, Shade D4, 1.2kg/Bo	US\$ 647.19	+Info
476590	Rodin Sculpture, Shade D3, 1.2kg/Bo	US\$ 647.19	+Info
476589	Rodin Sculpture, Shade D2, 1.2kg/Bo	US\$ 647.19	+Info
476588	Rodin Sculpture, Shade C4, 1.2kg/Bo	US\$ 647.19	+Info
476587	Rodin Sculpture, Shade C3, 1.2kg/Bo	US\$ 647.19	+Info
476586	Rodin Sculpture, Shade C1, 1.2kg/Bo	US\$ 647.19	+Info
476585	Rodin Sculpture, Shade A4, 1.2kg/Bo	US\$ 647.19	+Info
476584	Rodin Sculpture, Shade A3.5, 1.2kg/	US\$ 647.19	+Info
476574	Rodin Sculpture, Shade 0M3 (Bleach)	US\$ 647.19	+Info
476582	Rodin Sculpture, Shade B3, 1.2kg/Bo	US\$ 647.19	+Info
476581	Rodin Sculpture, Shade B2, 1.2kg/Bo	US\$ 647.19	+Info
476580	Rodin Sculpture, Shade 0M1, 1.2kg/B	US\$ 647.19	+Info
476579	Rodin Sculpture, Shade C2, 1.2kg/Bo	US\$ 647.19	+Info
476578	Rodin Sculpture, Shade A3, 1.2kg/Bo	US\$ 647.19	+Info
476577	Rodin Sculpture, Shade A2, 1.2kg/Bo	US\$ 647.19	+Info
476576	Rodin Sculpture, Shade A1, 1.2kg/Bo	US\$ 647.19	+Info
476575	Rodin Sculpture, Shade B1, 1.2kg/Bo	US\$ 647.19	+Info

300gr

300gr

Code	Name	Price	+
476619	Rodin Sculpture, Shade B4, 300g/Bot	US\$ 242.19	+Info
476627	Rodin Sculpture, Shade D4, 300g/Bot	US\$ 242.19	+Info
476626	Rodin Sculpture, Shade D3, 300g/Bot	US\$ 242.19	+Info
476625	Rodin Sculpture, Shade D2, 300g/Bot	US\$ 242.19	+Info
476624	Rodin Sculpture, Shade C4, 300g/Bot	US\$ 242.19	+Info
476623	Rodin Sculpture, Shade C3, 300g/Bot	US\$ 242.19	+Info
476622	Rodin Sculpture, Shade C1, 300g/Bot	US\$ 242.19	+Info
476621	Rodin Sculpture, Shade A4, 300g/Bot	US\$ 242.19	+Info
476620	Rodin Sculpture, Shade A3.5, 300g/B	US\$ 242.19	+Info
476610	Rodin Sculpture, Shade 0M3 (Bleach)	US\$ 242.19	+Info
476618	Rodin Sculpture, Shade B3, 300g/Bot	US\$ 242.19	+Info
476617	Rodin Sculpture, Shade B2, 300g/Bot	US\$ 242.19	+Info
476616	Rodin Sculpture, Shade 0M1, 300g/Bo	US\$ 242.19	+Info
476615	Rodin Sculpture, Shade C2, 300g/Bot	US\$ 242.19	+Info
476614	Rodin Sculpture, Shade A3, 300g/Bot	US\$ 242.19	+Info
476613	Rodin Sculpture, Shade A2, 300g/Bot	US\$ 242.19	+Info
476612	Rodin Sculpture, Shade A1, 300g/Bot	US\$ 242.19	+Info

476611 Rodin Sculpture, Shade B1, 300g/Bot

US\$ 242.19 +Info

600gr
600gr

Code	Name	Price	+
476601	Rodin Sculpture, Shade B4, 600g/Bot	US\$ 323.19	+Info
476609	Rodin Sculpture, Shade D4, 600g/Bot	US\$ 323.19	+Info
476608	Rodin Sculpture, Shade D3, 600g/Bot	US\$ 323.19	+Info
476607	Rodin Sculpture, Shade D2, 600g/Bot	US\$ 323.19	+Info
476606	Rodin Sculpture, Shade C4, 600g/Bot	US\$ 323.19	+Info
476605	Rodin Sculpture, Shade C3, 600g/Bot	US\$ 323.19	+Info
476604	Rodin Sculpture, Shade C1, 600g/Bot	US\$ 323.19	+Info
476603	Rodin Sculpture, Shade A4, 600g/Bot	US\$ 323.19	+Info
476602	Rodin Sculpture, Shade A3.5, 600g/B	US\$ 323.19	+Info
476592	Rodin Sculpture, Shade 0M3 (Bleach)	US\$ 323.19	+Info
476600	Rodin Sculpture, Shade B3, 600g/Bot	US\$ 323.19	+Info
476599	Rodin Sculpture, Shade B2, 600g/Bot	US\$ 323.19	+Info
476598	Rodin Sculpture, Shade 0M1, 600g/Bo	US\$ 323.19	+Info
476597	Rodin Sculpture, Shade C2, 600g/Bot	US\$ 323.19	+Info
476596	Rodin Sculpture, Shade A3, 600g/Bot	US\$ 323.19	+Info
476595	Rodin Sculpture, Shade A2, 600g/Bot	US\$ 323.19	+Info
476594	Rodin Sculpture, Shade A1, 600g/Bot	US\$ 323.19	+Info
476593	Rodin Sculpture, Shade B1, 600g/Bot	US\$ 323.19	+Info

SCULPTURE 2.0

SCULPTURE 2.0

1.2 kg
1.2 kg

Code	Name	Price	+
476529	Rodin Sculpture 2.0, Shade B4, 1.2k	US\$ 809.19	+Info
476537	Rodin Sculpture 2.0, Shade D4, 1.2k	US\$ 809.19	+Info
476536	Rodin Sculpture 2.0, Shade D3, 1.2k	US\$ 809.19	+Info
476535	Rodin Sculpture 2.0, Shade D2, 1.2k	US\$ 809.19	+Info
476534	Rodin Sculpture 2.0, Shade C4, 1.2k	US\$ 809.19	+Info
476533	Rodin Sculpture 2.0, Shade C3, 1.2k	US\$ 809.19	+Info
476532	Rodin Sculpture 2.0, Shade C1, 1.2k	US\$ 809.19	+Info
476531	Rodin Sculpture 2.0, Shade A4, 1.2k	US\$ 809.19	+Info
476530	Rodin Sculpture 2.0, Shade A3.5, 1.	US\$ 809.19	+Info
476520	Rodin Sculpture 2.0, Shade 0M3, 1.2	US\$ 809.19	+Info
476528	Rodin Sculpture 2.0, Shade B3, 1.2k	US\$ 809.19	+Info
476527	Rodin Sculpture 2.0, Shade B2, 1.2k	US\$ 809.19	+Info
476526	Rodin Sculpture 2.0, Shade 0M1, 1.2	US\$ 809.19	+Info
476525	Rodin Sculpture 2.0, Shade C2, 1.2k	US\$ 809.19	+Info
476524	Rodin Sculpture 2.0, Shade A3, 1.2k	US\$ 809.19	+Info
476523	Rodin Sculpture 2.0, Shade A2, 1.2k	US\$ 809.19	+Info
476522	Rodin Sculpture 2.0, Shade A1, 1.2k	US\$ 809.19	+Info
476521	Rodin Sculpture 2.0, Shade B1, 1.2k	US\$ 809.19	+Info

300gr
300gr

Code	Name	Price	+
476565	Rodin Sculpture 2.0, Shade B4, 300g	US\$ 323.19	+Info
476573	Rodin Sculpture 2.0, Shade D4, 300g	US\$ 323.19	+Info
476572	Rodin Sculpture 2.0, Shade D3, 300g	US\$ 323.19	+Info
476571	Rodin Sculpture 2.0, Shade D2, 300g	US\$ 323.19	+Info
476570	Rodin Sculpture 2.0, Shade C4, 300g	US\$ 323.19	+Info
476569	Rodin Sculpture 2.0, Shade C3, 300g	US\$ 323.19	+Info
476568	Rodin Sculpture 2.0, Shade C1, 300g	US\$ 323.19	+Info
476567	Rodin Sculpture 2.0, Shade A4, 300g	US\$ 323.19	+Info
476566	Rodin Sculpture 2.0, Shade A3.5, 30	US\$ 323.19	+Info
476556	Rodin Sculpture 2.0, Shade 0M3, 300	US\$ 323.19	+Info
476564	Rodin Sculpture 2.0, Shade B3, 300g	US\$ 323.19	+Info
476563	Rodin Sculpture 2.0, Shade B2, 300g	US\$ 323.19	+Info
476562	Rodin Sculpture 2.0, Shade 0M1, 300	US\$ 323.19	+Info
476561	Rodin Sculpture 2.0, Shade C2, 300g	US\$ 323.19	+Info
476560	Rodin Sculpture 2.0, Shade A3, 300g	US\$ 323.19	+Info
476559	Rodin Sculpture 2.0, Shade A2, 300g	US\$ 323.19	+Info
476558	Rodin Sculpture 2.0, Shade A1, 300g	US\$ 323.19	+Info
476557	Rodin Sculpture 2.0, Shade B1, 300g	US\$ 323.19	+Info

600gr
600gr

Code	Name	Price	+
476547	Rodin Sculpture 2.0, Shade B4, 600g	US\$ 404.19	+Info
476555	Rodin Sculpture 2.0, Shade D4, 600g	US\$ 404.19	+Info
476554	Rodin Sculpture 2.0, Shade D3, 600g	US\$ 404.19	+Info
476553	Rodin Sculpture 2.0, Shade D2, 600g	US\$ 404.19	+Info
476552	Rodin Sculpture 2.0, Shade C4, 600g	US\$ 404.19	+Info
476551	Rodin Sculpture 2.0, Shade C3, 600g	US\$ 404.19	+Info
476550	Rodin Sculpture 2.0, Shade C1, 600g	US\$ 404.19	+Info
476549	Rodin Sculpture 2.0, Shade A4, 600g	US\$ 404.19	+Info
476548	Rodin Sculpture 2.0, Shade A3.5, 60	US\$ 404.19	+Info
476538	Rodin Sculpture 2.0, Shade 0M3, 600	US\$ 404.19	+Info
476546	Rodin Sculpture 2.0, Shade B3, 600g	US\$ 404.19	+Info
476545	Rodin Sculpture 2.0, Shade B2, 600g	US\$ 404.19	+Info
476544	Rodin Sculpture 2.0, Shade 0M1, 600	US\$ 404.19	+Info
476543	Rodin Sculpture 2.0, Shade C2, 600g	US\$ 404.19	+Info
476542	Rodin Sculpture 2.0, Shade A3, 600g	US\$ 404.19	+Info
476541	Rodin Sculpture 2.0, Shade A2, 600g	US\$ 404.19	+Info
476540	Rodin Sculpture 2.0, Shade A1, 600g	US\$ 404.19	+Info
476539	Rodin Sculpture 2.0, Shade B1, 600g	US\$ 404.19	+Info

TITAN
TITAN

1.2 kg

1.2 kg

Code	Name	Price	+
476500	Rodin Titan, Shade 0M3, 1.2kg/Bottl	US\$ 809.19	+Info
476501	Rodin Titan, Shade B1, 1.2kg/Bottle	US\$ 809.19	+Info
476502	Rodin Titan, Shade A1, 1.2kg/Bottle	US\$ 809.19	+Info
476503	Rodin Titan, Shade A2, 1.2kg/Bottle	US\$ 809.19	+Info
476504	Rodin Titan, Shade A3, 1.2kg/Bottle	US\$ 809.19	+Info
476505	Rodin Titan, Shade C2, 1.2kg/Bottle	US\$ 809.19	+Info

300gr

300gr

Code	Name	Price	+
476513	Rodin Titan Shade 0M1 300g/Bottle	US\$ 323.19	+Info
476514	Rodin Titan, Shade 0M3, 300g/Bottle	US\$ 323.19	+Info
476515	Rodin Titan, Shade B1, 300g/Bottle	US\$ 323.19	+Info
476516	Rodin Titan, Shade A1, 300g/Bottle	US\$ 323.19	+Info
476517	Rodin Titan, Shade A2, 300g/Bottle	US\$ 323.19	+Info
476518	Rodin Titan, Shade A3, 300g/Bottle	US\$ 323.19	+Info
476519	Rodin Titan, Shade C2, 300g/Bottle	US\$ 323.19	+Info

600gr

600gr

Code	Name	Price	+
476506	Rodin Titan, Shade 0M1, 600g/Bottle	US\$ 404.19	+Info
476507	Rodin Titan, Shade 0M3, 600g/Bottle	US\$ 404.19	+Info
476508	Rodin Titan, Shade B1, 600g/Bottle	US\$ 404.19	+Info
476509	Rodin Titan, Shade A1, 600g/Bottle	US\$ 404.19	+Info
476510	Rodin Titan, Shade A2, 600g/Bottle	US\$ 404.19	+Info
476511	Rodin Titan, Shade A3, 600g/Bottle	US\$ 404.19	+Info
476512	Rodin Titan, Shade C2, 600g/Bottle	US\$ 404.19	+Info

VarseoSMILE | BEGO

BEGO

Crown Plus

Crown Plus

Code	Name	Price	+
464402	VarseoSmile Crown plus, A1, 250g	US\$ 251.16	+Info
464403	VarseoSmile Crown plus, A2, 250g	US\$ 251.16	+Info
464404	VarseoSmile Crown plus, A3, 250g	US\$ 251.16	+Info
464405	VarseoSmile Crown plus, B1, 250g	US\$ 251.16	+Info
464406	VarseoSmile Crown plus, B3, 250g	US\$ 251.16	+Info
464407	VarseoSmile Crown plus, C2, 250g	US\$ 251.16	+Info
464408	VarseoSmile Crown plus, D3, 250g	US\$ 251.16	+Info
475033	VarseoSmile Crown plus, BL, 250g	US\$ 251.16	+Info

TriniQ

TriniQ

Code	Name	Price	+
476872	VarseoSmileTriniQA1,500g	US\$ 315.00	+Info
476873	VarseoSmileTriniQA2,500g	US\$ 315.00	+Info
476874	VarseoSmileTriniQA3,500g	US\$ 315.00	+Info
476875	VarseoSmileTriniQB1,500g	US\$ 315.00	+Info
476876	VarseoSmileTriniQC2,500g	US\$ 315.00	+Info
476877	VarseoSmileTriniQ 0M1,500g	US\$ 315.00	+Info
476878	VarseoSmileTriniQ 0M3,500g	US\$ 315.00	+Info
476879	VarseoSmileTriniQ A3.5,500g	US\$ 315.00	+Info
476880	VarseoSmileTriniQ B2,500g	US\$ 315.00	+Info
476881	VarseoSmileTriniQ D3,500g	US\$ 315.00	+Info

AMANN GIRRBACH

AMANN GIRRBACH

Code	Name	Price	+
468910	NextDent C&B MFH / BL	US\$ 591.19	+Info
468911	NextDent C&B MFH / N1.5	US\$ 591.19	+Info
468912	NextDent C&B MFH / N1	US\$ 591.19	+Info
468913	NextDent C&B MFH / N2.5	US\$ 591.19	+Info
468914	NextDent C&B MFH / N2	US\$ 591.19	+Info
468915	NextDent C&B MFH / N3	US\$ 591.19	+Info

NEXTDENT

NEXTDENT

Code	Name	Price	+
416955	C&B MFH Bleach Color BL1 Lt	US\$ 546.25	+Info
416960	C&B MFH Bleach Color N1 Lt	US\$ 546.25	+Info
416961	C&B MFH Bleach Color N2 Lt	US\$ 546.25	+Info
416962	C&B MFH Bleach Color N1.5 Lt	US\$ 546.25	+Info
416963	C&B MFH Bleach Color N3 Lt	US\$ 546.25	+Info
416964	C&B MFH Bleach Color N2.5 Lt	US\$ 546.25	+Info

HATHO

POLISHING KIT

Thermoplastics
like *never* before



Packing Diamond D®
Denture Acrylic with
Tecnoflask



Keystone Industries: Laboratory



FREE SHIPPING
FOR ALL DOMESTIC ACCOUNTS!



US
DENTAL
DEPOT

+1 954 874 6325
www.usdentaldepot.com



ALLOY PRINTING

Printing Alloy



MEDILOY | BEGO

BEGO

Code	Name	Price	+
470780	Mediloy S-Co 5kg	US\$ 1,858.43	+Info

RITON 3D

RITON 3D

Code	Name	Price	+
474873	Cobalt-chromium alloy powder	US\$ 200.00	+Info
474874	Titanium alloy powder	US\$ 400.00	+Info

GLAZE

GLAZE

Code	Name	Price	+
476628	Rodin All-Purpose Glaze, 15g/Bottle	US\$ 39.69	+Info
476629	Rodin All-Purpose Glaze, 60g/Bottle	US\$ 96.39	+Info

C&B Temporary

C&B Temporary



VarseoSMILE TEMP | BEGO

BEGO

Code	Name	Price	+
464399	VarseoSmile Temp, A2, 250g	US\$ 149.00	+Info
464400	VarseoSmile Temp, A3, 250g	US\$ 149.00	+Info
464401	VarseoSmile Temp, C2, 250g	US\$ 149.00	+Info

ACKURETTA

ACKURETTA

Code	Name	Price	+
474972	CURO Crown 500g - A1	US\$ 155.93	+Info
474973	CURO Crown 500g - A2	US\$ 155.93	+Info
474974	CURO Crown 500g - A3.5	US\$ 155.93	+Info
474975	CURO Crown 500g - B3	US\$ 155.93	+Info

RITON 3D

RITON 3D

Code	Name	Price	+
474889	Temporary C&B	US\$ 340.00	+Info

DENTURE BASE 3D Resin

Denture



RODIN by PAC-DENT

DENTURE Base

Code	Name	Price	+
476630	Rodin Denture Base, Original Pink,	US\$ 323.19	+Info
476631	Rodin Denture Base , Meharry, 1kg/B	US\$ 323.19	+Info
476632	Rodin Denture Base, Light Pink, 1kg	US\$ 323.19	+Info
476633	Rodin Denture Base, Classic Pink, 1	US\$ 323.19	+Info

AMANN GIRRBACH

AMANN GIRRBACH

Code	Name	Price	+
468916	NextDent Denture 3D+ / Dark Pink	US\$ 438.30	+Info
468917	NextDent Denture 3D+ / Light Pink	US\$ 438.30	+Info
468918	NextDent Denture 3D+ /Opaque Pink	US\$ 438.30	+Info
468919	NextDent Denture 3D+ /red Pink	US\$ 521.77	+Info
468920	NextDent Denture 3D+ /Translucent P	US\$ 438.30	+Info

Denture TEETH

Denture Teeth

Code	Name	Price	+
420179	Denture Teeth A2 Resin 1 L	US\$ 565.58	+Info
420180	Denture Teeth A1 Resin 1 L	US\$ 565.58	+Info
464314	Denture Teeth A3 Resin 1 L	US\$ 565.58	+Info

NEXTDENT

NEXTDENT

Code	Name	Price	+
475036	DENTURE 3D+ Translucent Pink	US\$ 421.20	+Info
475037	DENTURE 3D+ Light Pink	US\$ 421.20	+Info
475038	DENTURE 3D+ Classic Pink	US\$ 421.20	+Info

RITON 3D

RITON 3D

Code	Name	Price	+
474887	Denture	US\$ 460.00	+Info

Gingiva MASK 3D Resin

Gingiva MASK



KeyMASK | KEYPRINT

KEYSTONE

Code	Name	Price	+
470563	KeyMask 0.5kg	US\$ 163.10	+Info

ACKURETTA

ACKURETTA

Code	Name	Price	+
474967	CURO Gingiva 500g - Pink	US\$ 155.93	+Info

AMANN GIRRBACH

AMANN GIRRBACH

Code	Name	Price	+
468921	NextDent Gingiva Mask	US\$ 268.84	+Info

NEXTDENT

NEXTDENT

Code	Name	Price	+
420394	NextDent Gingiva Mask	US\$ 290.01	+Info

RITON 3D

RITON 3D

Code	Name	Price	+
------	------	-------	---

474885	Gingiva mask	US\$ 164.00	+Info
--------	--------------	-------------	-----------------------

KITS / SET

KITS / SET



Code	Name	Price	+
470580	Lab resin intro kit, 0.5 kg	US\$ 366.94	+Info

MODELS

Model



KeyMODEL | KEYPRINT

KEYSTONE

Code	Name	Price	+
470570	KeyModel, 1kg	US\$ 157.29	+Info
470571	KeyModel 0.5kg	US\$ 102.83	+Info
470572	KeyModel Ultra Ivory 1kg	US\$ 80.99	+Info
470573	KeyModel Ultra Gray, 1kg	US\$ 80.99	+Info
470574	KeyModel Ultra Ivory 4kg	US\$ 278.25	+Info
470575	KeyModel Ultra Gray, 4kg	US\$ 278.25	+Info

NEXTDENT

NEXTDENT

Code	Name	Price	+
------	------	-------	---

419998	NextDent 2.0 Peach	US\$ 215.60	+Info
419999	NextDent 2.0 Grey	US\$ 215.60	+Info
420148	NextDent 2.0 White	US\$ 215.60	+Info

RODIN by PAC-DENT

PAC-DENT

Code	Name	Price	+
476640	Rodin Model (TAN), 1kg/Bottle	US\$ 161.19	+Info

VarseoModel | BEGO

BEGO

Code	Name	Price	+
419822	VarseoWax Model	US\$ 269.37	+Info

ACKURETTA

ACKURETTA

Code	Name	Price	+
474959	CURO Element Model 5kg - Beige	US\$ 425.25	+Info
474960	CURO Element Model 5kg - White	US\$ 425.25	+Info
474961	CURO Model 5kg - Yellow	US\$ 728.00	+Info
474962	CURO Element Model 500g - Beige	US\$ 51.98	+Info
474963	CURO Element Model 500g - White	US\$ 51.98	+Info
474964	CURO Model 500g - Yellow	US\$ 93.56	+Info
474965	CURO Model 1kg - Yellow	US\$ 165.38	+Info

AMANN GIRRBACH

AMANN GIRRBACH

Code	Name	Price	+
468922	NextDent Model 2.0 / Grey	US\$ 228.07	+Info
468923	NextDent Model 2.0 / Peach	US\$ 228.06	+Info
468924	NextDent Model 2.0 / White	US\$ 228.06	+Info

ANYCUBIC

ANYCUBIC

Code	Name	Price	+
416684	Anycubic UV Resin Gray 1 lt	US\$ 36.25	+Info

FORMLABS

FORMLABS

Code	Name	Price	+
417367	Model Dental for FORMLAB 1lt	US\$ 215.13	+Info
420178	White Resin 1lt	US\$ 215.13	+Info
500191	Resin Gray for FORMLAB 1lt	US\$ 215.13	+Info

HARVEST DENTAL

HARVEST DENTAL

Code	Name	Price	+
------	------	-------	---

475559	InovaPrint Model 3D Print Resin D	US\$ 194.08	+Info
475560	InovaPrint Model 3D Print Resin O	US\$ 194.08	+Info
475561	InovaPrint Model 3D Print Resin B	US\$ 194.08	+Info

Model STUDY

STUDY Model



Code	Name	Price	+
473957	SIMPLEX study model 800g	US\$ 62.78	+Info

RENFERT

RENFERT

Code	Name	Price	+
473958	SIMPLEX working model 800g	US\$ 62.78	+Info
473959	SIMPLEX aligner model 800g	US\$ 122.68	+Info
473960	SIMPLEX multi-use model 800g	US\$ 138.58	+Info

RITON 3D

RITON 3D

Code	Name	Price	+
474881	Dental model V1.0	US\$ 80.00	+Info
474883	Dental model V3.0	US\$ 80.00	+Info

FREE SHIPPING
FOR ALL DOMESTIC ACCOUNTS!



**US
DENTAL
DEPOT**

+1 954 874 6325
www.usdentaldepot.com



ORTHO IBT

ORTHO IBT



RODIN by PAC-DENT

PAC-DENT

Code	Name	Price	+
476641	Rodin Indirect Bonding Tray, 1kg/Bo	US\$ 201.69	+Info

KeyOrtho IBT | KEYPRINT

KEYSTONE

Code	Name	Price	+
470578	KeyOrtho IBT 1kg	US\$ 358.19	+Info
470579	KeyOrtho IBT .5kg	US\$ 230.09	+Info

ACKURETTA

ACKURETTA

Code	Name	Price	+
474970	CURO IBT 500g - Transparent	US\$ 155.93	+Info

AMANN GIRRBACH

AMANN GIRRBACH

Code	Name	Price	+
468925	NextDent Ortho IBT	US\$ 445.95	+Info

NEXTDENT

NEXTDENT

Code	Name	Price	+
420294	Ortho IBT Flexible material 1kg	US\$ 499.80	+Info

ORTHO Model

ORTHO Model



[KeyOrthoModel](#) | [KEYPRINT](#)

KEYSTONE

Code	Name	Price	+
470564	KeyOrthoModel 1kg	US\$ 157.29	+Info
470565	KeyOrthoModel, 0.5kg Bottle	US\$ 102.83	+Info

ORTHO Splint

ORTHO Splint



[RODIN by PAC-DENT](#)

PAC-DENT

Code	Name	Price	+
476634	Rodin FlexiSplint, Rigid 1kg/Bottle	US\$ 282.69	+Info
476635	Rodin FlexiSplint, Soft, 1kg/Bottle	US\$ 282.69	+Info

[KeySPLINT](#) | [KEYPRINT](#)

KEYSTONE

Code	Name	Price	+
470566	KeySplint Soft 1kg	US\$ 445.55	+Info
470567	KeySplint Soft 0.5kg	US\$ 256.27	+Info
470568	KeySplint Hard 1kg	US\$ 445.55	+Info
470569	Key Splint Hard 0.5 kg	US\$ 256.27	+Info

ACKURETTA

ACKURETTA

Code	Name	Price	+
474971	CURO ProSplint Soft 500g - Transpar	US\$ 155.93	+Info

AMANN GIRRBACH

AMANN GIRRBACH

Code	Name	Price	+
468926	NextDent Ortho Rigid	US\$ 387.33	+Info

NEXTDENT

NEXTDENT

Code	Name	Price	+
420151	NextDent Ortho Rigid Blue 1kg	US\$ 381.98	+Info
420263	Ortho Clear 3D print material 1kg	US\$ 556.22	+Info

RITON 3D

RITON 3D

Code	Name	Price	+
474890	Splint resin	US\$ 340.00	+Info

Surgical GUIDE

Surgical GUIDE



RODIN by PAC-DENT

PAC-DENT

Code	Name	Price	+
476636	Rodin Surgical Guide, 1kg/Bottle	US\$ 201.69	+Info

KeyGUIDE | KEYPRINT

KEYSTONE

Code	Name	Price	+
470561	Key Guide 1kg	US\$ 232.96	+Info
470562	KeyGuide 0.5kg	US\$ 151.48	+Info

ACKURETTA

ACKURETTA

Code	Name	Price	+
474968	CURO Guide 500g - Transparent	US\$ 155.93	+Info

AMANN GIRRBACH

AMANN GIRRBACH

Code	Name	Price	+
468927	NextDent SG / Orange	US\$ 342.74	+Info

FORMLABS

FORMLABS

Code	Name	Price	+
500190	Dental SG Resin for FORMLAB 1lt	US\$ 431.70	+Info

NEXTDENT

NEXTDENT

Code	Name	Price	+
420393	NextDent SG	US\$ 399.84	+Info

RITON 3D

RITON 3D

Code	Name	Price	+
474884	Surgical guide	US\$ 116.00	+Info

VarseoSURGICAL | BEGO

BEGO

Code	Name	Price	+
419823	VarseoWax Surgical Guide	US\$ 314.42	+Info

TRAYS

Rodin Bite Tray Is a biocompatible Class I dental material suitable for printing custom bite trays and impression trays. It is rigid and is compatible with all types of impression materials.



RODIN by PAC-DENT

PAC-DENT

Code	Name	Price	+
476637	Rodin Bite Tray, 1kg/Bottle	US\$ 201.69	+Info

KeyTRAY | KEYPRINT

KEYSTONE

Code	Name	Price	+
470576	Key Tray 1kg	US\$ 203.84	+Info
470577	KeyTray 0.5kg	US\$ 138.32	+Info

ACKURETTA

ACKURETTA

Code	Name	Price	+
474969	CURO Tray 500g - Orange	US\$ 155.93	+Info

AMANN GIRRBACH

AMANN GIRRBACH

Code	Name	Price	+
468931	NextDent Tray / Blue	US\$ 342.74	+Info
468932	NextDent Tray / Pink	US\$ 342.74	+Info

RITON 3D

RITON 3D

Code	Name	Price	+
474886	Individual impression tray	US\$ 230.00	+Info

TRY IN

Rodin Try-In material Is intended for producing Try-In denture models prior to printing the final denture base and teeth.



AMANN GIRRBACH

AMANN GIRRBACH

Code	Name	Price	+
468928	NextDent Try-In / T10	US\$ 332.30	+Info
468929	NextDent Try-In / T11	US\$ 415.37	+Info
468930	NextDent Try-In / T12	US\$ 415.37	+Info

RODIN by PAC-DENT

PAC-DENT

+ Price	Code 476639	Name Rodin Try-In, 1kg/Bottle	US\$ 201.69	+Info
---------	-------------	-------------------------------	-------------	-------

WAX Casting

WAX Casting



RODIN by PAC-DENT

PAC-DENT

Code	Name	Price	+
476638	Rodin Castable, 1kg/Bottle	US\$ 201.69	+Info

KeyCAST | KEYPRINT

KEYSTONE

Code	Name	Price	+
470559	KeyCast 1kg	US\$ 232.96	+Info
470560	KeyCast, .5kg	US\$ 151.48	+Info

ACKURETTA

ACKURETTA

Code	Name	Price	+
474966	CURO ProCast 500g - Red	US\$ 155.93	+Info

AMANN GIRRBACH

AMANN GIRRBACH

Code	Name	Price	+
468909	NextDent Cast / Purple	US\$ 268.85	+Info

ANYCUBIC

ANYCUBIC

Code	Name	Price	+
416686	Transl Green Castable Resin 1lt	US\$ 58.29	+Info

FORMLABS

FORMLABS

Code	Name	Price	+
419996	Resin Castable for FORMLAB-2 1lt	US\$ 431.70	+Info
500189	Resin Castable for FORMLAB-3 1lt	US\$ 431.70	+Info

HARVEST DENTAL

HARVEST DENTAL

Code	Name	Price	+
475558	InovaPrint Press 3D Print Resin O	US\$ 287.03	+Info

NEXTDENT

NEXTDENT

Code	Name	Price	+
420150	NextDent Cast - Purple	US\$ 290.01	+Info

RITON 3D

RITON 3D

Code	Name	Price	+
474888	Casting wax resin	US\$ 170.00	+Info

VarseoWAX | BEGO

BEGO

Code	Name	Price	+
419821	VarseoWax CAD/Cast	US\$ 269.37	+Info

HATHO POLISHING KIT

Thermoplastics
like *never* before



Packing Diamond D®
Denture Acrylic with
Tecnoflask



Keystone Industries: Laboratory



HATHO POLISHING KIT

Thermoplastics
like *never* before



Packing Diamond D®
Denture Acrylic with
Tecnoflask



Keystone Industries: Laboratory



US
DENTAL
DEPOT

FREE SHIPPING
FOR ALL DOMESTIC ACCOUNTS!



+1 954 874 6325

www.usdentaldepot.com



SPLINTS

SPLINTS Nightguards



ANYCUBIC

ANYCUBIC

Code	Name	Price	+
416685	Anycubic UV Resin Clear 1lt	US\$ 79.68	+Info

FORMLABS

FORMLABS

Code	Name	Price	+
419990	Dental LT Clear 1 L	US\$ 571.89	+Info
421112	Clear Resin 1lt	US\$ 215.13	+Info

CAD/CAM

ACCESSORIES for CAD/CAM MILLS



Accessories for Milling units - AMANN GIRRBACH

Accessories for Milling units - AMANN GIRRBACH



3D Accessories - AMANN GIRRBACH

3D Accessories - AMANN GIRRBACH

Code	Name	Price	+
124256	AG UNIVERSAL BONDING KIT	US\$ 78.71	+Info
420392	Blank Holder 98mm (5X)	US\$ 2,143.15	+Info
420988	Chuck 3mm Motion/ Motion2	US\$ 381.07	+Info

Other Models

Other Models

Code	Name	Price	+
470856	Ceramill Cooling 1lt	US\$ 112.11	+Info

Accessories for Milling units - CEREC

Accessories for Milling units - CEREC

Code	Name	Price	+
123891	CEREC MCXL Filter Replac each	US\$ 72.45	+Info
419964	MILLING Chamber Door Cpl Replacem	US\$ 795.75	+Info
420868	CEREC COMPACT Motor Replac each	US\$ 51.20	+Info

Accessories for Milling units - ECCO

Accessories for Milling units - ECCO

Code	Name	Price	+
113793	D-Form 87mm Block Holder	US\$ 12.80	+Info

Accessories for Milling units - ROLAND

Accessories for Milling units - ROLAND



SPINDLES

Code	Name	Price	+
420295	ROLAND DWX-4 Spindle Replacement Kit #ZS-4D - works on 52D 52DCI 51 DWX-4 (NOT ...	US\$ 1,436.85	+Info

DWX 52s & 53s

DWX dry 52s & 53s

Code	Name	Price	+
470861	Sensor Pin for DWX dry	US\$ 37.25	+Info

DWX-42W

DWX-42W

Code	Name	Price	+
464805	ABUTMENT KIT AK-1 FOR 42W	US\$ 1,239.94	+Info
464806	NT TRADING FIXTURE FOR 42W	US\$ 472.38	+Info
464807	GEOMEDI FIXTURE FOR 42W	US\$ 472.38	+Info
464808	MEDENTIKA FIXTURE FOR 42W	US\$ 472.38	+Info

466305	DWX-42W NEW Coolant Tank Filter 2Pk	US\$ 155.94	+Info
469625	Roland Collect for 42W	US\$ 582.35	+Info
469809	6-Pin (Multi-Pin) Adapter - 42W	US\$ 448.19	+Info
473667	DWX-42W Spindle Replacement Kit	US\$ 0.00	+Info

Other Models

Other Models

Code	Name	Price	+
421107	Tool Holder for Roland Mills	US\$ 77.68	+Info
421267	Acc-Clear Tubing	US\$ 16.54	+Info
468942	NEW C-Clamp Adaptor 52D	US\$ 448.19	+Info
474934	Disc Adapter STICKERS	US\$ 25.93	+Info

Roland

Roland

Code	Name	Price	+
464536	Roland Collect for 52dci	US\$ 463.53	+Info

ROLAND ACCESSORIES

Code	Name	Price	+
420349	ROLAND Grease 10ml can each	US\$ 23.53	+Info
420308	ROLAND Coolant Additive for DWX Wet Mills #ZAW-1000D - 1000 ml	US\$ 57.70	+Info
420309	ROLAND Chelating Agent for DWX Wet Mills #ZCH-250D - 250 ml	US\$ 119.95	+Info
420303	ROLAND Cleaning Tool for DWX-52DC #ZCT-4D	US\$ 78.35	+Info
420298	ROLAND TIPS TIME reduction KIT #CA-DK1-FX	US\$ 159.33	+Info
420296	ROLAND DWX Adapter for Model Time reduction Kit CA-52DC	US\$ 179.95	+Info
420304	ROLAND Disc Adapter 6-Pack for Roland 52DCI #ZMA-52DC-6 (disc holder for milli...	US\$ 269.95	+Info
420305	ROLAND Disc Adapter Rack #ZR-52D	US\$ 299.95	+Info
420307	ROLAND Pin-Type Material Adapter for Roland #ZV-52D	US\$ 487.45	+Info
420297	ROLAND TIME reduction KIT #CA-DK1	US\$ 1,249.95	+Info
420299	ROLAND DWX Time reduction Kit For Models #CA-MK1	US\$ 1,249.95	+Info

Accessories for Milling units - USDD

Accessories for Milling units - USDD

Code	Name	Price	+
420351	DESIGN STATION	US\$ 1,215.98	+Info

Accessories for Milling units - YETI

Accessories for Milling units - YETI

Code	Name	Price	+
474390	LENOVO Laptop , 15"	US\$ 1,232.04	+Info


CAD/CAM

ACCESSORIES for SINTERING

Sintering Beads



DIGITECH
DIGITECH

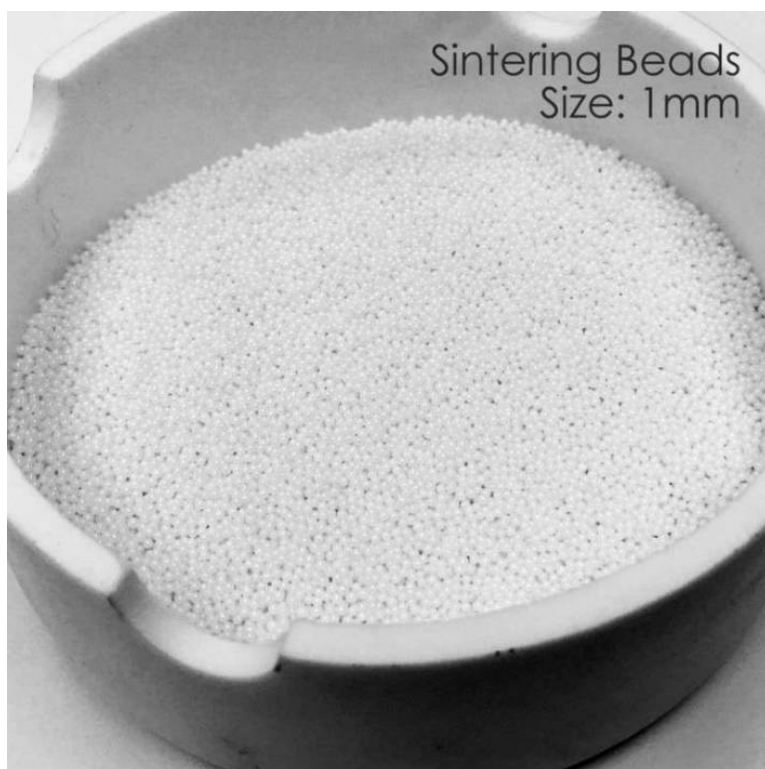

Code	Name	Price	+
115833	Un-Sintered Alum Beads 1mm 1oz	US\$ 5.22	+Info
116114	Un-Sintered Alum Beads 3mm 1oz	US\$ 11.90	+Info
116114-500	Sintering Zirconia Beads 3mm 500gr	US\$ 75.00	+Info
421078	Un-Sintered Alum Beads 1mm 200gr	US\$ 31.92	+Info
421079	Un-Sintered Alum Beads 2mm 200gr	US\$ 31.92	+Info

WHIP MIX

WHIP MIX

Code	Name	Price	+
120530	Sintering Beads 8oz Jar	US\$ 93.44	+Info
120572	Tempil Pellet Kit Pk 5 1300f	US\$ 38.42	+Info

Zirconia Sintering Trays



AMANN GIRRBACH

AMANN GIRRBACH

Other Models

Other Models

Code	Name	Price	+
467257	pliers for sinterbowl Therm 3	US\$ 91.74	+Info

Stackable Bowl

Stackable Bowl

Code	Name	Price	+
467245	Stackable sintering bowl	US\$ 356.52	+Info
467246	pliers for sintering shell, stackab	US\$ 86.64	+Info
467255	Stackable sintering bowl (ZOLID)	US\$ 318.53	+Info
467258	Therm 3 sinterbowl stackable	US\$ 356.52	+Info
467260	Therm 3 Zolid sinterbowl stackable	US\$ 349.17	+Info

DIGITECH

DIGITECH

Code	Name	Price	+
112243	Sagger Firing Platform w/lid 100mm	US\$ 51.07	+Info
116059	Closed Firing Crucible w/lid	US\$ 36.00	+Info
116601	Open Sagger Firing Platform	US\$ 62.50	+Info

WHIP MIX

WHIP MIX

Code	Name	Price	+
120533	Sintering Tray W/lid	US\$ 356.58	+Info
120534	Sintering Tray Spatula	US\$ 100.40	+Info

CAD/CAM

SCANNING SPRAY



SCAN-IT - TAUB

SCAN-IT-TAUB

A masking agent, SCAN-IT creates a proper surface color and reflectivity, allowing more accurate detail capture of ANY scanning surface. Adds a negligible micron coating thickness to the scanning surface. Use any die stone of your choice with accurate results. Easily washes or steams off the scanned surface and leaves no residue.



Code	Name	Price	+
122497	TAUB - Scan-it Thinner - 1/2 Oz - #43-4000	US\$ 9.79	+Info
122496	TAUB - Scan-it Liquid - 1/2oz (12 Cc) - #43-2000	US\$ 39.32	+Info

Scanning Spray - DENSYS

Scanning Spray - DENSYS

Code	Name	Price	+
401872	Speed Scan Spray 60ml	US\$ 84.00	+Info

Scanning Spray - KEYSTONE

Scanning Spray - KEYSTONE

Code	Name	Price	+
114803	CAD/CAM Scan Spray 88gr	US\$ 67.34	+Info

Scanning Spray - RENFERT

Scanning Spray - RENFERT

Code	Name	Price	+
120687	Renfert Scanspray 1x200ml 6.76	US\$ 53.20	+Info
473955	Scanspray extra fine 200ml	US\$ 79.81	+Info

Scanning Spray - YETI

Scanning Spray - YETI

Code	Name	Price	+
417095	Digiscan Spray 300 ml	US\$ 40.32	+Info

CAD/CAM

SOFTWARE



CAM SOFTWARE - HYPERDENT

CAM SOFTWARE - HYPERDENT

Features:

hyperDENT is the world leading software system in the dental CAM market and is valued among experts for its modular product structure. Our customers are provided with CAM solutions tailored to their specific application and implementation needs. Regardless of company size and whether for dental or practice labs, FOLLOW-ME! offers suitable hyperDENT product bundles for a perfect interaction between CAM and machines.

hyperDENT Computer Monitors

hyperDENT - Optimal Machine Workload.
hyperDENT can be utilized on a myriad of milling machines. The integrated milling cycles from the industrial segment are highly efficient and guarantee maximum process stability and top indication quality. High-grade security mechanisms enable a reliable production process for generating complex and stable tool paths. For example, the patented peeling strategies for glass ceramic, which offer perfect surface quality along with maximized tool life. When combined with the Template Generator Module, the milling strategy parameters can be adjusted and optimized for a specific machine. This allows our customers to get the most out of their machines, saving them time and money.

hyperDENT Overview


- * Open, highly automated and flexible system
- * Complete production process coverage with a single CAM software
- * Simple, intuitive operation
- * Maximum precision and time efficiency throughout the entire process
- * Material-independent
- * Flexible addition of software modules
- * Maximum freedom in creating templates
- * Additive manufacturing and milling all-in-one
- * Automatic material offset for post milling during the hybrid process
- * Automatic part identification using ID tags
- * Optimum surface finish quality due to proven tool path calculation
- * Consistent software enhancement and development
- * Worldwide premium support service

Computer Specs Recommendations

Processor (min.): Intel Core i5 2.8GHz

Processor (recomm.): Intel Core i7 3.2GHz

RAM memory (recomm.): 8 GB

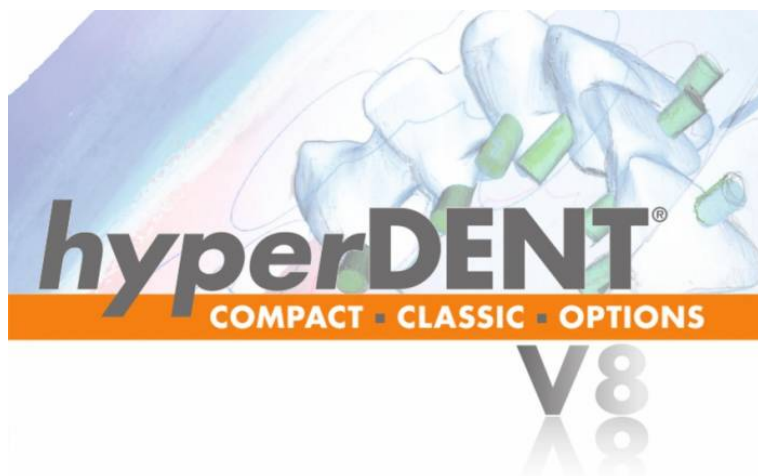

Graphic card: NVIDIA (A perfect function with other graphics cards cannot be guaranteed.)

VRAM Graphic card (min.): 2 GB

VRAM Graphic card (min.): 4 GB

HDD: 60 GB

For the calculations we recommend the best i7 you can afford, but there is little benefit for i9 currently.

32GB of RAM improves multi-tasking like 2-3+instances of hyperDENT plus other necessary programs like the CAD.

You need an NVIDIA card, technically there is a slight improvement for Quadro cards on CAD/Simulations, but a good Geforce works just fine as well.

PACKAGES

Code	Name	Price	+
420226	HYPERDENT - Compact 4X + PP + 1 year	US\$ 0.00	+Info
420225	HYPERDENT - Compact 5X + Roland PPs + MMU + 1 year # X5-1 - # X51	US\$ 0.00	+Info

OPTIONS

Code	Name	Price	+
420230	HYPERDENT - Additional PP for 4-Axis Compact	US\$ 682.50	+Info
420229	HYPERDENT - Additional PP for 5-Axis Compact	US\$ 1,365.00	+Info

CAM Software - UP3D

CAM Software - UP3D

Code	Name	Price	+
473694	Dental CAM Software UPCAM	US\$ 0.02	+Info

CAM Software - USDD

CAM Software - USDD

Code	Name	Price	+
420261	Complete Install & Basic Training	US\$ 1,248.00	+Info
420262	Extended 1 Day service	US\$ 525.01	+Info
420343	Complete Install & NOTraining	US\$ 472.50	+Info
420353	Travel Expenses	US\$ 160.00	+Info

Design Software for CAD - ROLAND

D-CLINIC Design Soft for CAD - ROLAND



Code	Name	Price	+
420366	DESIGN SOFTWARE d-CLINIC	US\$ 3,144.80	+Info

Design Software for CAD - UP3D

Design Soft for CAD - UP3D



Code	Name	Price	+
416557	UPCAD SOFTWARE	US\$ 1,099.98	+Info

CAD/CAM

ZIRCONIA COLORING



FULL CONTOUR COLORING ZIRCONIA LIQUIDS - DIGITECH by UPCERA

Full Contour Coloring Liquid Ecco Specially Designed For Super Translucent Zirconia
FDA CERTIFIED

DIGITECH Full Contour Liquid is made for the coloration of milled zirconia structures.

Coloring is applied before sintering, not for shading ceramic structures.
Alcohol based.

Available in 16 Vita shades:

- * A1 - A2 - A3.5 - A4
- * B1 - B2 - B3 - B4
- * C1 - C2 - C3 - C4
- * D2 - D3 - D4
- * O1- O2 for darker areas like fissures and occlusal valleys.
- * G1-G2-G3 for translucent like areas, like Cuspids and incisal borders.

DIGITECH Full Contour Liquid has no effect on the physical properties of DIGITECH Zirconia blocks such as:

- * Sintering density.
- * Flexural strength.
- * Fracture toughness and hardness.

Excess liquid will affect the stability of restoration and will corrode the heating element if **NOT PROPERLY DRIED WELL**

Note:

- 1 Make sure the restoration is fully dried before coloring.
- 2 Use palette or zirconia crucible to place the restoration (Sold Separately)
- 3 Use non-metal tweezers to operate, do not use metal tweezers. Gloves are recommended during peration.
- 4 The used color liquid should be kept separate from original bottle to avoid



contamination

5 Do not dilute with any water to avoid uneven coloring; use the diluting supplied in the package.

6 Keep the bottle sealed and store between the temperatures of 20-30C.

Drying:

1 - Infrared Drying: Keep a distance of 80-100mm between lamp and restoration. In case of brushing method, 30 minutes for crown and bridges. In case of dipping method, 30 minutes for crown and 30 to 60 minutes for bridge.

2 - Porcelain Furnace Drying: Set the temp range 200-250C with chamber half closed, holding time is the same above.

3 Sintering Oven Drying: Set the temperature at 200C with holding times same as above.

Sintering:

1 - 1480C is recommended, holding time of 2 hours.

But if a better translucency is expected, the temperature could be raised to 1530C to 1600C, but consider that a slight color fade could be possible, in case of that, adjust it with external staining, or use a deeper shade in coloring process.

2 - Place the restoration in the crucible with inside-up (internal part of the crowns facing up) to avoid the sintering beads getting into the restoration.

Make sure the restoration gets maximum supporting points from the sintering bead.

In case of long bridge remain attached to sintered with a supportive base to lower the odds of distortion.

DENTINE COLORS (Vita Shades)

Code	Name	Price	+
114427	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia A1 - 50 ml bottle - ...	US\$ 24.00	+Info
114428	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia A2 - 50ml bottle - T...	US\$ 24.00	+Info
114429	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia A3 - 50ml bottle - T...	US\$ 24.00	+Info
114430	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia A3.5- 50ml bottle - ...	US\$ 24.00	+Info
114431	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia A4 - 50ml bottle- To...	US\$ 23.52	+Info
114432	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia B1 - 50ml bottle - T...	US\$ 24.00	+Info
114433	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia B2 - 50ml bottle - T...	US\$ 23.52	+Info
114434	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia B3- 50ml bottle - To...	US\$ 23.52	+Info
114435	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia B4- 50ml bottle- To ...	US\$ 24.00	+Info
114436	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia C1- 50ml bottle - To...	US\$ 24.00	+Info
114437	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia C2- 50ml bottle - To...	US\$ 23.52	+Info
114438	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia C3- 50ml bottle - To...	US\$ 24.00	+Info
114439	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia C4 - 50ml bottle - T...	US\$ 24.00	+Info
114440	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia D2 - 50ml Bottle - T...	US\$ 23.52	+Info
114441	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia D3 - 50ml bottle - T...	US\$ 23.52	+Info
114443	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia D4 - 50ml bottle - T...	US\$ 24.00	+Info

GUM / PINK coloring

Code	Name	Price	+
115085	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia Gum PINK P2 - 30ml b...	US\$ 18.40	+Info
115086	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia Gum Dark PINK P3 - 3...	US\$ 18.40	+Info
469788	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia Gum PINK P1 - 30ml b...	US\$ 18.40	+Info

ENAMELS



Code	Name	Price	+
115083	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia Incisal Grey G1 - 30...	US\$ 18.40	+Info
115084	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia Incisal Grey G2 - 30...	US\$ 18.40	+Info

Master Kit / Set

Full Contour

Code	Name	Price	+
470296	Full Cont Liq Gum P1 30ml	US\$ 512.80	+Info

OCCLUSAL (Ochres)



Code	Name	Price	+
115081	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia Ochre O1 - 30ml bott...	US\$ 18.03	+Info
115082	DIGITECH - Full Contour Coloring Liquid for Dental Zirconia Ochre O2 - 30ml bott...	US\$ 16.56	+Info

Coloring for Zirconia - AMANN GIRRBACH

Coloring for Zirconia - AMANN GIRRBACH

Code	Name	Price	+
468206	Ceramill Liquid CL1	US\$ 109.58	+Info
468207	Ceramill Liquid CL2	US\$ 109.58	+Info
468208	Ceramill Liquid CL3	US\$ 109.58	+Info

468209	Ceramill Liquid CL4	US\$ 109.58	+Info
--------	---------------------	-------------	-----------------------

COLOR LIQUID PRETAU - ZIRKON ZAHN



Code	Name	Price	+
473764	Color Liquid Prettau Aquarell B3 50	US\$ 145.00	+Info
473772	Color Liquid Prettau Aquarell Incis	US\$ 145.00	+Info
473771	Color Liquid Prettau Aquarell Incis	US\$ 145.00	+Info
473770	Color Liquid Prettau AquarellIncisa	US\$ 145.00	+Info
473769	Colour Liquid Prettau AquarellTissu	US\$ 145.00	+Info
473768	Colour Liquid Prettau AquarellTissu	US\$ 145.00	+Info
473767	Colour Liquid Prettau AquarellTissu	US\$ 145.00	+Info
473766	Color Liquid Prettau Aquarell C3 50	US\$ 145.00	+Info
473765	Color Liquid Prettau Aquarell C1 50	US\$ 145.00	+Info
473755	Colour Liquid WaterbasedA1 - 100ML	US\$ 79.55	+Info
473763	Color Liquid Prettau Aquarell B2 50	US\$ 145.00	+Info
473762	Color Liquid Prettau Aquarell B1 50	US\$ 145.00	+Info
473761	Color Liquid Prettau Aquarell A3.5	US\$ 145.00	+Info
473760	Color Liquid Prettau Aquarell A3 50	US\$ 145.00	+Info
473759	Color Liquid Prettau Aquarell A2 50	US\$ 145.00	+Info
473758	Color Liquid Prettau Aquarell A1 50	US\$ 145.00	+Info
473756	Color Liquid B1 100ml WaterBased	US\$ 79.55	+Info

Coloring for Zirconia - JSP

Coloring for Zirconia - JSP

Code	Name	Price	+
500082	COLORING LIQUID CUP - 10pk	US\$ 11.04	+Info

REFILLABLE STAIN BRUSH - SHARK

REFILLABLE STAIN BRUSH - SHARK

Features:

Rechargeable brushes for Zirconia Size S M L

Package: 1 unit



Code	Name	Price	+
123019	SHARK - Refillable Stain Brush Pen for Coloring Zirconia - Medium each	US\$ 10.75	+Info

123020	SHARK - Refillable Stain Brush Pen for Coloring Zirconia - Large each	US\$ 11.20	+Info
--------	---	------------	-----------------------

DIP & BRUSH K2 CAD/CAM ZIRKON LIQUID - YETI

YETI - DIP & BRUSH K2 CAD/CAM ZIRKON LIQUID

With K2 Dip & Brush it is possible to individualize parts of frameworks (Brush) or to colorize whole frameworks by dipping (Dip) before sinterisation. Water based colours with 100 ml volume are available in 7 dentin and 2 enamel colours.

Dipping technique:

Dentin colorization of frameworks by dipping them is applicable mainly in layering technique. It allows the reduction of zircon whiteness. Thus, less light is reflected and it leads to outstanding aesthetic results particularly in the thin areas.

Painting technique:

Painting technique (Brush) is suitable mainly for monolithic restorations which can be individualized with different dentin and enamel colours before sinterization. It makes the finalization with stains and glaze powder easier after sinterization.

Quantity: 100ml



Code	Name	Price	+
417031	YETI Dip & Brush K2 CAD CAM Zirkon Liquid A1. 100ml. #397-5001	US\$ 66.26	+Info
417032	YETI Dip & Brush K2 CAD CAM Zirkon Liquid A2. 100ml #397-5002	US\$ 73.82	+Info
417033	YETI Dip & Brush K2 CAD CAM Zirkon Liquid A3. 100ml. #397-5003	US\$ 78.53	+Info
417034	YETI Dip & Brush K2 CAD CAM Zirkon Liquid A3.5 100ml. #397-5004	US\$ 66.26	+Info
417035	YETI Dip & Brush K2 CAD CAM Zirkon Liquid B2. 100ml. #397-5007	US\$ 66.26	+Info
417036	YETI Dip & Brush K2 CAD CAM Zirkon Liquid C2. 100ml. #397-5011	US\$ 66.26	+Info
417037	YETI Dip & Brush K2 CAD CAM Zirkon Liquid D2. 100ml. #397-5014	US\$ 78.53	+Info
417039	YETI Dip & Brush K2 CAD CAM Zirkon Liquid GR3. 100ml. #397-5019	US\$ 79.51	+Info

Catalog Index / Indice

MILLING DISCS & BLOCKS	1
CAD/CAM EQUIPMENT	90
CAD/CAM ABRASIVES	104
3D Accessories	126
3D RESINS	127
ACCESSORIES for CAD/CAM MILLS	147
ACCESSORIES for SINTERING	149
SCANNING SPRAY	152
SOFTWARE	152
ZIRCONIA COLORING	155